

Dexcom G5[™] Mobile CGM Enabled



User Guide



t:slim X2 Insulin Pump with Dexcom G5 Mobile CGM User Guide

Congratulations on the purchase of your new t:slim $X2^{\text{\tiny M}}$ System. Your decision to use insulin pump therapy and continuous glucose monitoring is a sign of your commitment to your diabetes care. We recognize and respect the importance of your decision.

This User Guide is designed to assist you with the features and functions of the t:slim X2 System. It provides important warnings and cautions on proper operation as well as technical information to ensure your safety. It also provides step-by-step instructions on how to properly program, manage and care for your t:slim X2 System.

Changes in equipment, software, or procedures occur periodically; information describing these changes will be included in future editions of this User Guide. Please contact your local Tandem Diabetes Care representative to obtain a replacement copy of the User Guide that is the correct version for your pump.

No part of this publication may be reproduced, stored in a retrieval system, nor transmitted in any form or by any means, electronic or mechanical, without the prior written permission of Tandem Diabetes Care.

Tandem Diabetes Care, Inc. 11075 Roselle Street San Diego, CA 92121 USA +1 858 255 6269 www.tandemdiabetes.com

If you have any questions, contact your local Tandem Diabetes Care representative.

www.tandemdiabetes.com/contact

© 2018 Tandem Diabetes Care, Inc. All rights reserved. Tandem Diabetes Care, touch simplicity, and t:slim are registered trademarks in the US, and t:slim X2 is a trademark, registered trademark, and/or the subject of pending trademark applications in other territories worldwide.

Software Version Pendleton

DEXCOM, DEXCOM G5, DEXCOM G5 MOBILE CGM, and SEVEN/SEVEN PLUS are trademarks of Dexcom, Inc. HUMALOG is a trademark of Eli Lilly and Company. NOVOLOG and NOVORAPID are trademarks of Novo Nordisk A/S. COMFORT is a trademark of Unomedical A/S. The *Bluetooth*® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Tandem Diabetes Care is under license.

All other trademarks and copyrights are the property of their respective owners.

EC REP

MDSS GmbH Schiffgraben 41 30175 Hannover Germany



Section 1 – Before You Begin

Chapto	er 1 – Introduction	
1.1	System Description	18
1.2	About this User Guide	18
1.3	Important User Information	19
1.4	Important Pediatric User Information	19
1.5	Conventions of this Guide	21
Chapte	er 2 – Important Safety Information	
2.1	Indications for Use	24
2.2	Contraindications	24
2.3	System Warnings	25
2.4	System Precautions	29
2.5	Potential Benefits From Using the System	33
2.6	Possible Risks From Using the System	
2.7	Working with your Healthcare Provider	36
2.8	Emergency Kit	36
2.9	Verification of Proper Functionality	
2.10	Wireless Co-existence and Data Security	37

Chapter 3 – Getting to Know Your t:slim X2 System

3.1	What your t:slim X2 System Package Includes
3.2	System Terminology
3.3	Explanation of System Symbols 4-
3.4	Explanation of System Icons
3.5	Explanation of System Colors
3.6	Screen Lock
3.7	CGM Screen Lock
3.8	Home Screen
3.9	CGM Home Screen 5
3.10	Status Screen 5
3.11	Bolus Screen 6
3.12	Options Screen
3.13	My Pump Screen 6
3.14	My CGM Screen
3.15	Number Keypad Screen
3.16	Letter Keypad Screen

Section 2 – Key Pump Features

Chapter	r 4 – Getting Started
4.1 4.2 4.3 4.4 4.5 4.6 4.7 4.8	Charging the t:slim X2 Pump .74 Using the Touch Screen .75 Turning on the t:slim X2 Pump Screen .76 Unlocking the t:slim X2 Pump Screen .76 Selecting Your Language .76 Editing Time .77 Editing Date .77 Turning the Pump Off .77
	r 5 – Infusion Site Care and Loading Cartridge
5.1 5.2 5.3 5.4 5.5	Infusion Site Selection and Care
5.6	Setting Site Reminder

Chapter 6 – Personal Profiles 6.1 6.2 Programming a New Personal Profile94 6.3 6.4 6.5 6.6 6.7 6.8 Chapter 7 – Bolus 7.1 7.2 7.3 7.4 7.5 7.6 Chapter 8 – Stop/Resume Insulin 8.1 8.2

Section 3 – Additional Pump Operations

Chapte	r 9 – Temporary Basal Rate	
9.1 9.2	Setting a Temp Rate	10 10
Chapte	r 10 – Quick Bolus	
10.1 10.2	Setting Up Quick Bolus	14 14
Chapte	r 11 – t:slim X2 Pump Settings	
11.1 11.2 11.3 11.4	Display Settings1Mobile Device Connection1t:slim X2 Sound Volume1Turn Feature Lock ON or OFF1	18 19
Chapte	r 12 – t:slim X2 Pump Info and History	
12.1 12.2	t:slim X2 Pump Info	22 22

Chapte	r 13 – t:slim X2 Pump Reminders
13.1 13.2 13.3 13.4 13.5	Low BG Reminder124High BG Reminder125After Bolus BG Reminder125Missed Meal Bolus Reminder126Site Reminder126
Chapte	r 14 – User Settable Alerts & Reminders
14.1	Low Insulin Alert
14.2	Auto-Off Alarm
_	ion 4 – Pump Safety Alerts & Reminders
Sect	

Table of Contents

15.7	Incomplete Fill Tubing Alert	139
15.8	Incomplete Fill Cannula Alert	140
15.9	Incomplete Setting Alert	141
15.10	Basal Rate Required Alert	142
15.11	Max Hourly Bolus Alert	143
15.12	Max Bolus Alert 1	144
15.13	Max Bolus Alert 2	145
15.14	Max Basal Alert 1	146
15.15	Max Basal Alert 2	147
15.16	Min Basal Alert 1	148
15.17	Min Basal Alert 2	149
15.18	Connection Error Alert	150
15.19	Pairing Code Timeout	151
15.20	Power Source Alert	152
15.21	Data Error Alert	153
Chapte	r 16 – t:slim X2 Pump Alarms	
16.1	Resume Pump Alarm	157
16.2	Low Power Alarm 1	158
16.3	Low Power Alarm 2	
16.4	Empty Cartridge Alarm	
16.5	Cartridge Error Alarm	

16.6	Cartridge Removal Alarm
16.7	Temperature Alarm
16.8	Occlusion Alarm 1
16.9	Occlusion Alarm 2
16.10	Screen On/Quick Bolus Button Alarm
16.11	Altitude Alarm
16.12	Reset Alarm
Chapte	er 17 – t:slim X2 Pump Malfunction
17.1	Malfunction
Sect	tion 5 - Getting Started with Your CGM System
Chapte	er 18 – Dexcom G5 Mobile CGM Overview
18.1	CGM System Overview
18.2	Sensor Overview
18.3	Receiver (t:slim X2 Pump) Overview
18.4	Transmitter Overview

Table of Contents

Unapte	er 19 – CGM Seπings	
19.1 19.2 19.3 19.4 19.5	About Bluetooth Disconnecting from the Dexcom Receiver Entering Your Transmitter ID Setting CGM Volume CGM Info	180 180 181
Chapte	er 20 – Setting CGM Alerts	
20.1 20.2 20.3 20.4 20.5 20.6	Setting Your High Glucose Alert and Repeat Feature Setting Your Low Glucose Alert and Repeat Feature Rate Alerts Setting Your Rise Alert Setting Your Fall Alert Setting Your Out of Range Alert	187 187 188 188
Chapte	er 21 – Inserting the Sensor and Placing the Transmitter	
21.1 21.2 21.3 21.4	Before You Start	195 195 195
21.5	Inserting the CGM Sensor	196

21.6	Attaching the CGM Transmitter	198
Chapte	r 22 – Starting a CGM Sensor Session	
22.1 22.2	Starting Sensor	
Chapte	r 23 – Calibrating Your CGM System	
23.1 23.2 23.3 23.4 23.5	Calibration Overview 2 Startup Calibration 2 Calibration Blood Glucose Value and Correction Bolus 2 12 Hour Calibration Update 2 Other Reasons You May Need to Calibrate 2	207 208 208
Chapte	r 24 – Viewing CGM Data on Your t:slim X2 Pump	
24.1 24.2 24.3 24.4 24.5	Overview	213 214 217

Table of Contents

Chapter	25 – Ending Your CGM Sensor Session	
25.2	Automatic Sensor Shut-Off	222
Secti	on 6 - CGM Alerts, Errors and Troubleshooting	
Chapter	26 – CGM Alerts and Errors	
26.1	Startup Calibration Alert	229
26.2	Second Startup Calibration Alert	230
26.3	12 Hour Calibration Alert	
26.4	Incomplete Calibration	232
	Calibration Timeout	
26.6	Wait 15 Minute Calibration Error Alert	234
	Calibration Required Alert	
	CGM High Alert	
	CGM Low Alert	
26.10	CGM Fixed Low Alert	238
26.11	CGM Rise Alert	239
	CGM Rapid Rise Alert	
	CGM Fall Alert	

26.14	CGM Rapid Fall Alert	242
26.15	Unknown Sensor Reading	
26.16	Out of Range Alert	
26.17	Low Transmitter Battery Alert	245
26.18	Transmitter Error	246
26.19	Failed Sensor Error	247
26.20	CGM System Error	248
Chapte	er 27 – CGM Troubleshooting	
27.1	CGM Pairing Troubleshooting	252
27.2	Calibration Troubleshooting	252
27.3	Sensor Insertion Troubleshooting	252
27.4	Unknown Sensor Reading Troubleshooting	
27.5	Out of Range/No Antenna Troubleshooting	254
27.6	Failed Sensor Troubleshooting	
27.7	Sensor Inaccuracies	255
Sect	tion 7 - Living With and Caring for Your System	
Chapte	er 28 – Lifestyle Issues and Travel	
28.1	Overview	260

Chapter 29 – Taking Care of Your t:slim X2 System		
29.1	Overview	
Section 8 - Technical Specifications		
Chapte	r 30 – Technical Specifications	
30.1 30.2 30.3 30.4 30.5 30.6 30.7 30.8 30.9 30.10 30.11 30.12	Overview268t:slim X2 System Specifications269t:slim X2 Pump Options and Settings278t:slim X2 Pump Performance Characteristics281Dexcom G5 Mobile CGM Transmitter Performance Characteristics283Dexcom G5 Mobile CGM Sensor Performance Characteristics284Electromagnetic Compatibility326Electromagnetic Emissions328Electromagnetic Immunity329Distances Between the t:slim X2 Pump and RF Equipment333Quality of Wireless Service335FCC Notice Concerning Interference335	
Index.		

THIS PAGE IS INTENTIONALLY LEFT BLANK

Section 1

Before You Begin

Chapter 1

Introduction

1.1 System Description

The t:slim X2™ Insulin Pump is made up of the t:slim X2 Insulin Pump and the t:slim 3mL (300 units) cartridge. The t:slim X2 Insulin Pump delivers insulin in two ways: continuous, or basal insulin delivery, and bolus insulin delivery to cover carbohydrates eaten (food bolus) and to lower high blood glucose (correction bolus). The disposable cartridge is filled with up to 300 units of U-100 insulin and attached to the pump. The cartridge is replaced every few days.

The t:slim X2 Insulin Pump and the Dexcom CGM Mobile Transmitter (sold separately) paired together may be referenced as the "System."

The Dexcom G5 Mobile Sensor is a disposable device that is inserted under the skin to continuously monitor glucose levels for up to 7 days. The Dexcom G5 Mobile Transmitter connects to the sensor pod and wirelessly sends readings to the pump display, in place of the Dexcom Receiver, every 5 minutes. The display shows sensor glucose readings, trend graph, direction and

rate of change arrows. The sensor is discarded after a session of up to 7 days. The transmitter is reusable and is replaced about every 3 months.

The sensor measures glucose in the fluid under the skin—not in blood, and sensor readings are not identical to readings from a blood glucose meter. You still need a blood glucose meter to calibrate your CGM on a regular basis to help ensure the accuracy of sensor glucose readings.

The t:slim X2 Insulin Pump can be used for basal and bolus insulin delivery with or without Dexcom G5 Mobile CGM. If the Dexcom G5 Mobile Sensor and Transmitter are not used, sensor glucose readings will not be sent to the pump display and you will not receive any sensor glucose alerts.

The Dexcom G5 Mobile Sensors and Transmitter are sold and shipped separately by Dexcom.

■ NOTE: Device Connections

The Dexcom G5 Mobile CGM only allows pairing with one medical device at a time (either the t:slim X2 Pump or the Dexcom receiver), but you

can still use the Dexcom Mobile App and your t:slim X2 Pump simultaneously using the same transmitter ID.

1.2 About this User Guide

This User Guide covers important information on how to operate your t:slim X2 System. It provides step-by-step instructions to help you properly program, manage and care for the System. It also provides important warnings and cautions on proper operation and technical information to ensure your safety.

In this Guide, the t:slim X2 Pump may be referred to as "your pump" or "your t:slim X2 Pump." The Dexcom G5 Mobile Transmitter may be referred to as "your transmitter." The Dexcom G5 Mobile Sensor may be referred to as "your sensor." Together, the Dexcom G5 Mobile Transmitter and Dexcom G5 Mobile Sensor may be referred to as "your CGM."

The User Guide is organized into sections. Section 1 provides important information you need to know before you start using the System. Sections 2-

4 cover instructions for using the t:slim X2 Insulin Pump. Sections 5–6 cover instructions for using Dexcom G5 Mobile CGM with your t:slim X2 Insulin Pump.

Pump screens used in this Guide to demonstrate how to use features are examples only. They should not be considered suggestions for your individual needs.

1.3 Important User Information

Review all instructions in this User Guide before using the System.

If you are not able to use the System according to the instructions in this User Guide, you may be putting your health and safety at risk.

If you are currently using the t:slim X2 Pump without Dexcom G5 Mobile CGM, or if you are currently using Dexcom G5 Mobile CGM with the Dexcom G5 Mobile Receiver, it is still very important that you review all instructions in this User Guide before using the combined System.

Pay special attention to Warnings and Precautions in this User Guide. Warnings and Precautions are identified with a **A** symbol.

If you still have questions after reading this User Guide, contact your local Tandem Diabetes Care representative.

1.4 Important Pediatric User Information

The following recommendations are meant to help younger users and their caregivers program, manage and care for the System.

Younger children may inadvertently press or tap the pump, leading to unintentional delivery of insulin.

It is the responsibility of the healthcare provider and caregiver to determine if the user is appropriate for treatment with this device.

We recommend reviewing the Quick Bolus and Feature Lock capabilities of the Tandem pump and determining how they best fit with your care plan. These features are detailed in Chapters 10 Quick Bolus and 11 t:slim X2 Pump Settings of this User Guide.

Inadvertent dislodgement of the infusion site may occur more frequently with children so consider securing the infusion site and tubing.

A WARNING

DO NOT allow small children (either pump users or non-users) to ingest small parts, such as the rubber USB port cover and cartridge components. Small parts could pose a choking hazard. If ingested or swallowed, these small component pieces may cause internal injury or infection.

A WARNING

The System includes parts (such as the USB cable and infusion set tubing) that could pose a strangulation or asphyxiation hazard. Always use the appropriate length of infusion set tubing and arrange cables and tubing to minimize the risk of strangulation. **ENSURE** that these parts are stored in a secure place when not in use.

A WARNING

For patients who do not self-manage their disease, the Feature Lock function should **ALWAYS**

Chapter 1 – Introduction

be **ON** when the pump is not being used by a caregiver. The Feature Lock function is intended to prevent inadvertent button presses that may lead to insulin delivery or changes in the pump settings. These changes can potentially lead to hypoglycemic or hyperglycemic events.

A WARNING

For patients whose insulin administration is managed by a caregiver, ALWAYS turn off the Quick Bolus feature to avoid inadvertent bolus delivery. If the Feature Lock is turned on, the Quick Bolus feature is automatically disabled. Inadvertent button presses or tampering with the insulin pump could result in over delivery or under delivery of insulin. This can cause very low or very high blood glucose which could result in serious injury or death.

1.5 Conventions of this Guide

The following are conventions used in this Guide (such as terms, icons, text formatting, and other conventions) along with their explanations.

Convention	Explanation
Bolded Text	Text that is in bold and in a different font than the rest of the sentence or step indicates an onscreen or physical button name.
Touchscreen	The front glass screen of your pump, which displays all programming, operating, and alarm/alert information.
Тар	Quickly and lightly touch the screen with your finger.
Press	Use your finger to depress a physical button (the Screen On/Quick Bolus Button is the only physical/hardware button on your t:slim X2 Pump).
Hold	Keep pressing a button or touching an icon or menu until its function is complete.
Menu	A list of options on your touch screen that allow you to perform specific tasks.
Icon	An image on your touch screen that indicates an option or item of information, or a symbol on the back of your t:slim X2 Pump or its packaging.

THIS PAGE IS INTENTIONALLY LEFT BLANK

Chapter 2

Important Safety Information

2.1 Indications for Use

Indications for Use

The t:slim X2 System consists of the t:slim X2 Insulin Pump paired with the Dexcom G5 Mobile CGM (Continuous Glucose Monitor).

The t:slim X2 Insulin Pump is intended for the subcutaneous delivery of insulin, at set and variable rates, for the management of diabetes mellitus in persons requiring insulin. The t:slim X2 Insulin Pump can be used solely for continuous insulin delivery and as part of the t:slim X2 System to receive and display continuous glucose measurements from the Dexcom G5 Mobile Sensor and Transmitter.

When used with the Dexcom G5 Mobile CGM, the t:slim X2 System is designed to replace fingerstick blood glucose testing for diabetes treatment decisions. The t:slim X2 System aids in the detection of episodes of hyperglycemia and hypoglycemia, facilitating both acute and long-term therapy adjustments, which may minimize these

excursions. Interpretation of the t:slim X2 System results should be based on the trends and patterns seen with several sequential readings over time.

The t:slim X2 System is indicated for use in individuals 6 years of age and greater.

The t:slim X2 System is intended for single patient use and requires a prescription.

The device is indicated for use with U-100 insulin.

2.2 Contraindications

Contraindications

The t:slim X2 System is not intended for anyone unable or unwilling to:

- Test blood glucose (BG) levels as recommended by your healthcare provider
- Demonstrate adequate carbohydrate-counting skills (preferred, not required)

- Maintain sufficient diabetes selfcare skills
- See your healthcare provider(s) regularly

You must also have adequate vision and/or hearing in order to recognize your System alerts.

The t:slim X2 Pump, Dexcom G5 Mobile Transmitter, and Dexcom G5 Mobile Sensor must be removed before Magnetic Resonance Imaging (MRI), Computed Tomography (CT) scan, or diathermy treatment. Exposure to MRI, CT, or diathermy treatment can damage the System.

Taking medications with acetaminophen/paracetamol while wearing the sensor may falsely raise your sensor glucose readings. The level of inaccuracy depends on the amount of acetaminophen/paracetamol active in your body and may be different for each person.

2.3 System Warnings

WARNINGS: t:slim X2 Insulin Pump

DO NOT start to use your t:slim X2 System before reading the User Guide. Failure to follow the instructions in the User Guide can result in over delivery or under delivery of insulin. This can cause very low or very high blood glucose which could result in serious injury or death. If you have questions or need further clarification on your System use, ask your healthcare provider or contact your local Tandem Diabetes Care representative.

DO NOT start to use your t:slim X2 System before you have been appropriately trained on its use by a certified t:slim X2 System trainer. Consult with your healthcare provider for your individual training needs for the entire t:slim X2 System. Failure to complete the necessary training on the System could result in serious injury or death.

DO NOT use any other insulin with your System other than U-100. Use of insulin

with lesser or greater concentration can result in under delivery or over delivery of insulin. This can cause very high or a very low blood glucose.

DO NOT put any other drugs or medications inside your System cartridge. The System is designed only for Continuous Subcutaneous Insulin Infusion (CSII) using Humalog or NovoLog/NovoRapid insulin. Use of other drugs or medications can damage the pump and result in injury if infused.

DO NOT start to use your System before consulting with your healthcare provider to determine which features are most appropriate for you. Only your healthcare provider can determine and help vou adjust vour Basal Rate(s). Carb Ratio(s), Correction Factor(s), Target BG, and duration of insulin action. In addition, only your healthcare provider can determine your CGM settings and how you should use your sensor trend information to help you manage vour diabetes. Incorrect settings can result in over delivery or under delivery of insulin. This can cause very low or very high blood glucose.

DO be prepared to inject insulin with an alternative method if delivery is interrupted for any reason. Your t:slim X2 System is designed to deliver insulin reliably, but because it uses only rapidacting insulin, you will not have longacting insulin in your body. Failure to have an alternative method of insulin delivery can lead to very high blood glucose or Diabetic Ketoacidosis (DKA).

DO use only cartridges and insulin infusion sets with matching connectors and follow the instructions for use. Failure to do so may result in over delivery or under delivery of insulin and may cause very low or very high blood glucose.

DO NOT place your infusion set on any scars, lumps, moles, stretch marks or tattoos. Placing your infusion set in these areas can cause swelling, irritation or infection. This can affect insulin absorption and cause high or low blood glucose.

DO NOT ignore infusion set cannula fractures. Infusion set cannulas may fracture on rare occasions. If an infusion set cannula breaks and no portion of it is visible above the skin, do not attempt

Chapter 2 – Important Safety Information

to remove it. Seek professional medical help if you have symptoms of infection or inflammation—redness, swelling or pain—at the insertion site. If you experience a broken infusion set cannula, please report this to your local Tandem Diabetes Care representative.

NEVER fill your tubing while your infusion set is connected to your body. Always ensure that the infusion set is disconnected from your body before filling the tubing. Failure to disconnect your infusion set from your body before filling the tubing can result in over delivery of insulin. This can cause serious injury or death from very low blood glucose.

DO NOT reuse cartridges or use cartridges other than those manufactured by Tandem Diabetes Care, Inc. Use of cartridges not manufactured by Tandem Diabetes Care or reuse of cartridges may result in over delivery or under delivery of insulin. This can cause very low or very high blood glucose.

ALWAYS twist the tubing connector between the cartridge tubing and the infusion set tubing an extra quarter of a turn to ensure a secure connection. A loose connection can cause insulin to leak, resulting in under delivery of insulin. This can cause high blood glucose.

DO NOT disconnect the tubing connector between the cartridge tubing and the infusion set tubing. If the connection comes loose, disconnect the infusion set from your body before tightening. Failure to disconnect before tightening can result in over delivery of insulin. This can cause low blood glucose.

DO NOT remove or add insulin from a filled cartridge after loading onto the pump. This will result in an inaccurate display of the insulin level on the Home Screen and you could run out of insulin before the pump detects an empty cartridge. This can cause very high blood glucose, or Diabetic Ketoacidosis (DKA).

DO NOT deliver a bolus until you have reviewed the calculated bolus amount on the pump display. If you dose an insulin amount that is too high or too low, this could cause very high or very low blood glucose. You can always

adjust the insulin units up or down before you decide to deliver your bolus.

DO NOT allow small children (either pump users or non-users) to ingest small parts, such as the rubber USB port cover and cartridge components. Small parts could pose a choking hazard. If ingested or swallowed, these small component pieces may cause internal injury or infection.

The System includes parts (such as the USB cable and infusion set tubing) that could pose a strangulation or asphyxiation hazard. ALWAYS use the appropriate length of infusion set tubing and arrange cables and tubing to minimize the risk of strangulation. ENSURE that these parts are stored in a secure place when not in use.

For patients who do not self-manage their disease, the Feature Lock function should ALWAYS be ON when the pump is not being used by a caregiver. The Feature Lock function is intended to prevent inadvertent button presses that may lead to insulin delivery or changes in the pump settings. These

changes can potentially lead to hypoglycemic or hyperglycemic events.

For patients whose insulin administration is managed by a caregiver, AI WAYS turn off the Quick Bolus feature to avoid inadvertent bolus delivery.

If the Feature Lock is turned on, the Quick Bolus feature is automatically disabled. Inadvertent button presses or tampering with the insulin pump could result in over delivery or under delivery of insulin. This can cause very low or very high blood glucose which could result in serious injury or death.



WARNINGS: Radiology and Medical Procedures and your t:slim X2 System

ALWAYS notify the provider/technician about your diabetes and your t:slim X2 System. If you need to discontinue use of the System for medical procedures, follow your healthcare provider's instructions to replace missed insulin when you reconnect to the pump. Check your blood glucose before disconnecting from the pump and again when you reconnect and treat high

blood alucose levels as recommended by your healthcare provider.

DO NOT expose your pump, transmitter, or sensor to:

- X-ray
- Computed Tomography (CT) scan
- Magnetic Resonance Imaging (MRI)
- Positron Emission Tomography (PET) scan
- Other exposure to radiation

The t:slim X2 System is magnetic resonance (MR) Unsafe. You must take off vour pump, transmitter, and sensor and leave them outside the procedure room if you are going to have any of the above procedures.

In addition to the above, DO NOT expose your pump, transmitter, or sensor to:

Pacemaker/Automatic Implantable Cardioverter Defibrillator (AICD) placement or reprogramming

- Cardiac Catheterization
- Nuclear Stress Test

You must take off your pump, transmitter, and sensor and leave them outside the procedure room if you are going to have any of the above medical procedures.

There are other procedures where you should proceed with caution:

- Laser Surgery Your System can usually be worn during the procedure. However, some lasers can create interference and cause the System to alarm.
- General Anesthesia Depending on the equipment being used, you may or may not need to remove your System. Be sure to ask your healthcare provider.

There is no need to disconnect for electrocardiograms (EKGs) or colonoscopies. If you have questions, contact vour local Tandem Diabetes Care representative.

Chapter 2 – Important Safety Information

WARNINGS: Using Dexcom G5 Mobile with your t:slim X2 **Insulin Pump**

DO NOT ignore symptoms of high and low glucose. If your sensor glucose alerts and readings do not match your symptoms, measure your blood glucose with a blood glucose meter even if your sensor is not reading in the high or low range.

CALIBRATE your CGM at least once every 12 hours. Calibrating less often than every 12 hours might cause sensor glucose readings to be inaccurate and glucose alerts to become unreliable. This could result in you missing severe hypoglycemia (low blood glucose) or hyperalycemia (high blood glucose) events.

DO NOT ignore sensor fractures. Sensors may fracture on rare occasions. If a sensor breaks and no portion of it is visible above the skin, do not attempt to remove it. Seek professional medical help if you have symptoms of infection or inflammation—redness, swelling, or pain—at the insertion site. If you experience a broken sensor, please report

this to your local Tandem Diabetes Care representative.

DO NOT use Dexcom G5 Mobile CGM in pregnant women or persons on dialysis. The System is not approved for use in pregnant women or persons on dialysis and has not been evaluated in these populations. Sensor glucose readings may be inaccurate in these populations and could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

DO NOT use Dexcom G5 Mobile CGM in critically ill patients. It is not known how different conditions or medications common to the critically ill population may affect the performance of the System. Sensor glucose readings may be inaccurate in critically ill patients, and solely relying on the sensor glucose alerts and readings for treatment decisions could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

DO NOT insert the sensor in sites other than the abdomen (belly) or upper but-

tocks (for ages 6-17 only). Other sites have not been studied and are not approved. Use in other sites might cause sensor glucose readings to be inaccurate and could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood alucose) events.

DO NOT expect CGM alerts until after the 2-hour startup. You will NOT get any sensor glucose readings or alerts until after the 2-hour startup ends AND you complete the startup calibration. During this time you might miss severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

DO NOT use your transmitter if it is damaged/cracked. This could create an electrical safety hazard or malfunction, which might cause electrical shocks.

STORE the Dexcom G5 Mobile CGM sensor at temperatures between 2°C (36°F) to 25°C (77°F) for the length of the sensor's shelf life. You may store the sensor in the refrigerator if it is within this temperature range. The sensor should not be stored in a freezer. Storing the sensor improperly might cause the sensor glucose readings to be inaccurate, and you might miss severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

DO NOT allow young children to hold the sensor, transmitter or transmitter kit box without adult supervision. The sensor and transmitter include small parts that may pose a choking hazard. Keep the transmitter kit box away from young children; it contains a magnet that should not be swallowed.

System Precautions



PRECAUTIONS: t:slim X2 Insulin Pump

DO NOT open or attempt to repair your System. The System is a sealed device that should be opened and repaired only by Tandem Diabetes Care. Modification could result in a safety hazard. If your System seal is broken, the System is no longer watertight and the warranty is voided.

CHANGE your infusion set every 48 to 72 hours as recommended by your healthcare provider. Wash your hands with anti-bacterial soap before handling the infusion set and thoroughly clean the insertion site on your body to avoid infection. Contact your healthcare provider if you have symptoms of infection at your insulin infusion site.

Al WAYS remove all air bubbles from the System before beginning insulin delivery. Ensure there are no air bubbles when drawing insulin into the filling syringe, hold the pump with the white fill port pointed up when filling the tubing, and ensure that there are no air bubbles in the tubing when filling. Air in the system takes space where insulin should be and can affect insulin delivery.

CHECK your infusion site daily for proper placement and leaks, REPLACE your infusion set if you notice leaks around the site. Improperly placed sites or leaks around the infusion site can result in under delivery of insulin.

CHECK your infusion set tubing daily for any leaks, air bubbles, or kinks, Air in the tubing, leaks in the tubing, or kinked tubing may restrict or stop insulin delivery and result in under delivery of insulin.

CHECK the tubing connection between your cartridge tubing and infusion set tubing daily to ensure it is tight and secure. Leaks around the tubing connection can result in under delivery of insulin.

DO NOT change your infusion set before bedtime or if you will not be able to test your blood glucose 1-2 hours after the new infusion set is placed. It is important to confirm that the infusion set is inserted correctly and delivering insulin. It is also important to respond quickly to any problems with the insertion to ensure continued insulin delivery.

ALWAYS check that your cartridge has enough insulin to last through the night before going to bed. If you are sleeping, you could fail to hear the Empty Cartridge Alarm and miss part of your basal insulin delivery.

CHECK your System's personal settings regularly to ensure they are cor-

Chapter 2 – Important Safety Information

rect. Incorrect settings can result in over delivery or under delivery of insulin. Consult with your healthcare provider as needed.

ALWAYS make sure that the correct time and date are set on your System. When using AM/PM, always check that the AM/PM setting is accurate. Not having the correct time and date setting may affect safe insulin delivery.

CONFIRM that the screen display turns on, you can hear audible beeps, feel the pump vibrate, and see the green LED light blinking around the edge of the Screen On Button when you connect a power source to the USB port. These features are used to notify you about alerts, alarms, and other conditions that require your attention. If these features are not working, discontinue use of the System and contact your local Tandem Diabetes Care Customer representative.

CHECK your System regularly for potential alarm conditions that may display. It is important to be aware of conditions that may affect insulin delivery and require your attention so you can respond as soon as possible.

DO NOT use the vibrate feature for alerts and alarms during sleep unless otherwise directed by your healthcare provider. Having the volume for alerts and alarms set to high will help ensure that you don't miss an alert or alarm.

ALWAYS look at the screen to confirm correct programming of the bolus amount when you first use the Quick Bolus feature. Looking at the screen will ensure that you are correctly using the beep/vibration commands to program the intended bolus amount.

DO NOT use your pump if you think it might be damaged due to dropping it or hitting it against a hard surface. Check that the System is working properly by plugging a power source into the USB port and confirming that the display turns on, you hear audible beeps, feel the pump vibrate, and see the green LED light blinking around the edge of the Screen On Button. If you are unsure about potential damage, discontinue use of the System and con-

tact your local Tandem Diabetes Care representative.

AVOID exposure of your System to temperatures below 5°C (40°F) or above 37°C (99°F). Insulin can freeze at low temperatures or degrade at high temperatures. Insulin that has been exposed to conditions outside of the manufacturer's recommended ranges can affect the safety and performance of the System.

AVOID submersing your pump in fluid beyond a depth of 1 meter or for more than 30 minutes (IPX7 rating). If your pump has been exposed to fluid beyond these limits, check for any signs of fluid ingress. If there are signs of fluid entry, discontinue use of the pump and contact your local Tandem Diabetes Care representative.

AVOID areas where there may be flammable anesthetics or explosive gases. The System is not suitable for use in these areas and there is a risk of explosion. Remove your pump if you need to enter these areas.

MAKE SURE to not move further than the length of the USB cable when you are connected to the pump and to a charging source. Moving further than the length of the USB cable may cause the cannula to be pulled out of the infusion site. For this reason it is recommended not to charge the pump while sleeping.

DISCONNECT your infusion set from your body while on high-speed/high gravity amusement park thrill rides. Rapid changes in altitude or gravity can affect insulin delivery and cause injury.

DISCONNECT your infusion set from your body before flying in an aircraft without cabin pressurization or in planes used for aerobatics or combat simulation (pressurized or not). Rapid changes in altitude or gravity can affect insulin delivery and cause injury.

CONSULT your healthcare provider about lifestyle changes such as weight gain or loss, and starting or stopping exercise. Your insulin needs may change in response to lifestyle changes. Your basal rate(s) and other settings may need adjustment.

CHECK your blood glucose using a blood glucose meter following a gradual elevation change of up to 300 meters, such as when snow skiing or driving on a mountain road. Delivery accuracy can vary up to 15% until 3 units of total insulin have been delivered or elevation has changed by more than 300 meters. Changes in delivery accuracy can affect insulin delivery and cause injury.

ALWAYS check with your healthcare provider for specific guidelines if you want or need to disconnect from the pump for any reason. Depending on the length of time and reason you are disconnecting, you may need to replace missed basal and/or bolus insulin. Check your blood glucose before disconnecting from the pump and again when you reconnect, and treat high blood glucose (BG) levels as recommended by your healthcare provider.

ENSURE that your personal insulin delivery settings are programmed into the pump before you use the System if you receive a warranty replacement. Failure to enter your insulin delivery settings could result in over delivery or under delivery of insulin. This can cause

very low or very high blood glucose which could result in serious injury or death. Consult your healthcare provider as needed.

Interference with your System's electronics by cell phones can occur if worn in close proximity. It is recommended that your pump and cell phone be worn at least 16 centimeters apart.

ALWAYS dispose of used System components such as cartridges, syringes, needles, infusion sets, and CGM sensors following the instructions from your healthcare provider. Wash your hands thoroughly after handling used System components.



PRECAUTIONS: Using Dexcom G5 Mobile CGM with your t:slim X2 Insulin Pump

DO NOT open the sensor package until you have washed your hands with soap and water, and let them dry. You may contaminate the insertion site and suffer an infection if you have dirty hands while inserting the sensor.

Chapter 2 – Important Safety Information

DO NOT insert the sensor until you have cleaned the skin with a topical antimicrobial solution, such as isopropyl alcohol, and allowed the skin to dry. Inserting into unclean skin might lead to infection. Do not insert the sensor until the cleaned area is dry so the sensor adhesive will stick better.

AVOID using the same spot repeatedly for sensor insertion. Rotate your sensor placement sites, and do not use the same site for two sensor sessions in a row. Using the same site might cause scarring or skin irritation.

AVOID inserting the sensor in areas that are likely to be bumped, pushed or compressed, or areas of skin with scarring, tattoos, or irritation as these are not ideal sites to measure glucose. Insertion in those areas might affect accuracy and could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

AVOID injecting insulin or placing an infusion set within 8 centimeters of the sensor. The insulin might affect sensor accuracy and could result in you miss-

ing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

DO NOT use the sensor if its sterile package has been damaged or opened. Using an unsterile sensor might cause infection.

To calibrate the CGM System, DO enter the exact blood glucose value that your blood glucose meter displays within 5 minutes of a carefully performed blood glucose measurement. Do not enter sensor glucose readings for calibration. Entering incorrect blood glucose values, blood glucose values obtained more than 5 minutes before entry, or sensor glucose readings might affect sensor accuracy and could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

DO NOT calibrate if your blood glucose is changing at a significant rate, typically more than 0.11 mmol/L per minute. Do not calibrate when your receiver screen is showing the rising or falling single arrow or double arrow, which indicates that your blood glucose is rapidly rising

or falling. Calibrating during significant rise or fall of blood glucose may affect sensor accuracy and could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

The System accuracy may be affected when your glucose is changing at a significant rate (e.g., 0.11 to 0.17 mmol/L or more than 0.17 mmol/L each minute), such as during exercise or after a meal.

AVOID separating the transmitter and pump by more than 6 meters. The transmission range from the transmitter to the pump is up to 6 meters without obstruction. Wireless communication does not work well through water so the range is much less if you are in a pool, bathtub, or on a water bed, etc. Types of obstruction differ and have not been tested. If your transmitter and pump are farther than 6 meters apart or are separated by an obstruction, they might not communicate or the communication distance may be shorter and result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

DO NOT use alternative blood glucose site testing (blood from your palm or forearm, etc.) for calibration. Alternative site blood glucose values may be different than those taken from a fingerstick blood glucose value and may not represent the timeliest blood glucose value. Use a blood glucose value taken only from a fingerstick for calibration. Alternative site blood glucose values might affect sensor accuracy and result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

ENSURE that your transmitter ID is programmed into the pump before you use the System if you receive a warranty replacement pump. The pump cannot communicate with the transmitter unless the transmitter ID is entered. If the pump and transmitter are not communicating, you will not receive sensor glucose readings and you might miss severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

DO NOT discard your transmitter. It is reusable. The same transmitter is used

for each session until you have reached the end of the transmitter battery life.

The Dexcom G5 Mobile Sensor, Transmitter, and Receiver are not compatible with the SEVEN/SEVEN PLUS Transmitter and Receiver. Different generations will not connect with each other and will not work.

2.5 Potential Benefits From Using the System

- The t:slim X2 System provides an automated way to deliver basal and bolus insulin. Delivery can be finetuned based on up to 6 customizable Personal Profiles, each with up to 16 time-based settings for Basal Rate, Carb Ratio, Correction Factor, and Target BG. In addition, the Temp Rate feature allows you to program a temporary basal rate change for up to 72 hours.
- The t:slim X2 System gives you the option of delivering a bolus all at once, or delivering a percentage over an extended period of time

without navigating to different menus. You can also program a bolus more discreetly using the Quick Bolus feature, which can be used without looking at the pump, and can be programmed in increments of either units of insulin or grams of carbohydrate.

- From the bolus screen, the "calculator within a calculator" feature allows you to enter multiple carbohydrate values and add them together. The System's bolus calculator will recommend a bolus based on the entire amount of carbohydrates entered, which can help eliminate guesswork.
- The System keeps track of the amount of active insulin from food and correction boluses (IOB). When programming additional food or correction boluses, the pump will subtract the amount of IOB from the recommended bolus if your BG is below the target set in your active Personal Profile. This can help prevent insulin stacking, which can lead to hypoglycemia.

Chapter 2 – Important Safety Information

- You can program a number of reminders that will prompt you to retest your BG after a low or high BG is entered, as well as a "Missed Meal Bolus Reminder" which will alert you if a bolus isn't entered during a specified period of time. If activated, these can help reduce the likelihood that you will forget to check your blood glucose or bolus for meals.
- You have the ability to view a variety of data right on your screen, including the time and amount of your last bolus, your total insulin delivery by day, as well as broken into basal, food bolus, and correction bolus.
- When paired with Dexcom G5
 Mobile Transmitter and Sensor,
 your t:slim X2 System can receive
 CGM readings every 5 minutes,
 which are displayed as a trend
 graph on the Home Screen. You
 can also program your pump to
 alert you when your CGM readings
 are above or below a given level, or
 are rising or falling quickly. Unlike
 the readings from a standard blood

- glucose meter, CGM readings allow you to view trends in real time, as well as capture information when you would otherwise be unable to check your blood sugar, such as while you are asleep. This information can be useful for you and your healthcare provider when considering changes to your therapy. In addition, the programmable alerts can help you to spot potential low or high blood glucose sooner than you would using a only a blood glucose meter.
- CGM use has been shown to increase time in your target glucose range, without increasing time spent above or below your target range. Real-time CGM can help improve diabetes control (lower A1C values, reducing alveemic variability and time spent in low and high blood glucose ranges)1, 2, 3 which can help reduce diabetes related complications.4,5 These benefits can be seen especially with using real-time CGM at least 6 days per week² and can be sustained over time.⁶ In some cases, patients perceived an increase in their quality

- of life and peace of mind when using real-time CGM as well as reporting a high satisfaction with CGM.⁷
- ¹ Garg S, Zisser H, Schwartz S, Bailey T, Kaplan R, Ellis S, Jovanovic L. Improvement in glycemic excursions with a transcutaneous, real-time continuous glucose sensor: a randomized controlled trial. Diabetes Care. 2006; 29:44-50.
- ² JDRF CGM Study Group. Continuous Glucose Monitoring and Intensive Treatment of Type 1 Diabetes. NEJM 2008;359:1464-76.
- ³ Battelino. Effect of continuous glucose monitoring of hypoglycemia in type 1 diabetes. Diabetes Care 2011; 34(4): 795-800.
- ⁴The Diabetes Control and Complications Research Group. The effect of intensive treatment of diabetes on the development and progression of long-term complications of insulin-dependent diabetes mellitus. N Eng J Med. 1993; 329:997-1036.

⁵ Ohkubo Y, Kishikawa H, Araki E, et al. Intensive insulin therapy prevents progression of diabetic microvascular complications in Japanese patients with non-insulin dependent diabetes mellitus: a randomized prospective 6-year study. Diabetes Res Clin Pract. 1995; 28:103-117.

⁶ JDRF CGM Study Group. Sustained Benefit of Continuous Glucose Monitoring on A1c, Glucose Profiles, and Hypoglycemia in Adults With Type 1 Diabetes, Diabetes Care 2009; 32: 2047-2049.

⁷JDRF CGM Study Group. Quality-of-Life Measures in Children and Adults With Type 1 Diabetes. Diabetes Care 2010; 33; 2175-2177.

2.6 Possible Risks From Using the System

As with any medical device, there are risks associated with using the t:slim X2 System. Many of the risks are common to insulin therapy in general, but there are additional risks associated with continuous insulin infusion and continuous insulin infusion and continuous.

ous glucose monitoring. Reading your User Guide and following the Instructions for Use are critical for the safe operation of your system. Consult your healthcare provider about how these risks may impact you.

Risks associated with using the pump functionality of the System

Inserting and wearing an infusion set might cause infection, bleeding, pain or skin irritations (redness, swelling, bruising, itching, scarring or skin discoloration).

There is a remote chance that an infusion set cannula fragment could remain under your skin if the cannula breaks while you are wearing it. If you think a cannula has broken under your skin, contact your healthcare provider and contact your local Tandem Diabetes Care representative.

Other risks associated with infusion sets include occlusions and air bubbles in the tubing, which can affect insulin delivery.

A PRECAUTION

At very low basal rates checking your blood glucose frequently may give you an early indication of an occlusion.

Risks that could result from pump failure include the following:

- possible hypoglycemia (low blood glucose) from over-delivery of insulin due to a hardware defect
- hyperglycemia (high blood glucose) and ketosis possibly leading to Diabetic Ketoacidosis (DKA) due to pump failure resulting in cessation of insulin delivery due to either a hardware defect or software anomaly.

Risks associated with using the CGM functionality of the System

Inserting the sensor and wearing the adhesive patch might cause infection, bleeding, pain and skin irritations (readiness, swelling, bruising, itching, scarring, or skin discoloration).

There is a remote chance that a sensor fragment could remain under your skin

if the sensor breaks while you are wearing it. If you think a sensor has broken under your skin, contact your healthcare provider and call your local Tandem Diabetes Care representative.

Other risks associated with CGM use include the following:

- You will not get sensor glucose alerts when the alert function is turned off, your transmitter and pump are out of range, or when your pump is not showing sensor glucose readings. You might not notice alerts if you are unable to hear them or feel the vibration.
- There are a number of risks as a result of the fact that the Dexcom G5 Mobile CGM takes readings from fluid below the skin (interstitial fluid) instead of blood. There are differences in how glucose is measured in the blood compared to how it is measured in interstitial fluid, and glucose is absorbed into the interstitial fluid slower than it is absorbed into the blood, which can cause CGM readings to lag behind

readings from a blood glucose meter.

2.7 Working with your Healthcare Provider

Any clinical language presented in this User Guide is based on the assumption that you have been educated by your healthcare provider on certain terms and how they apply to you in your diabetes management. Your healthcare provider can help you establish diabetes management guidelines that best fit your lifestyle and needs.

Monitor your blood glucose (BG) with the guidance of your healthcare provider. According to the American Association of Diabetes Educators' white paper "Insulin Pump Therapy: Guidelines for Successful Outcomes," patients should routinely check their BG levels at least 4 times daily (optimally 6 to 8 times daily) in order to detect hyperglycemia (high blood glucose) and hypoglycemia (low blood glucose) early. Undetected hyperglycemia or hypoglycemia can result without proper monitoring.

Consult your healthcare provider before using the System to determine which features are most appropriate for you. Only your healthcare provider can determine and help you adjust your basal rate(s), insulin-to-carbohydrate ratio(s), correction factor(s), blood glucose (BG) target, and duration of insulin action. In addition, only your healthcare provider can determine your CGM settings and how you should use your sensor trend information to help you manage your diabetes.

2.8 Emergency Kit

Make sure that you always have an insulin syringe and vial of insulin or insulin pen with you as a backup for emergency situations. You should also always have an appropriate emergency kit with you. Talk with your healthcare provider regarding what items this kit should include.

Supplies to carry every day:

 Blood glucose testing supplies: meter, strips, control solution, lancets, meter batteries

- Fast-acting carbohydrate to treat low blood glucose
- Extra snack for longer coverage than fast-acting carbohydrate
- Glucagon emergency kit
- Rapid-acting insulin vial and syringes or insulin pen
- Infusion sets (minimum of 2)
- Insulin pump cartridges (minimum of 2)
- Infusion site preparation products (antiseptic wipes, skin adhesive)
- Diabetes identification card or jewelry

2.9 Verification of Proper Functionality

A power supply (AC adapter with micro-USB connector) is provided as part of the System. Before using your System, ensure that the following occur when you connect a power supply into the USB port of your t:slim X2 Pump:

- You hear an audible alert
- Your see the green light illuminate from the edge around the Screen On/Quick Bolus Button
- You feel a vibratory alert
- You see a charge symbol (lightning bolt) on the battery level indicator

In addition, before using the System, ensure the following:

- Press the Screen On/Quick Bolus Button to turn the screen on so that you can see the display
- When the display screen is on, the touch screen responds to your finger tap

A PRECAUTION

CONFIRM that the screen display turns on, you can hear audible beeps, feel the pump vibrate, and see the green LED light blinking around the edge of the Screen On Button when you connect a power source to the USB port. These features

are used to notify you about alerts, alarms, and other conditions that require your attention. If these features are not working, discontinue use of the System and contact your local Tandem Diabetes Care representative.

2.10 Wireless Co-existence and Data Security

The t:slim X2 System is designed to work safely and effectively in the presence of wireless devices typically found at home, work, retail stores, and places of leisure where daily activities occur. See Section 30.11 Quality of Wireless Service for more information.

The t:slim X2 System is designed to send and accept Bluetooth® Low Energy (BLE) communication. Communication is not established until you enter the appropriate credentials into your pump.

The t:slim X2 System and system components ensure data security via proprietary means and ensure data integrity

Chapter 2 – Important Safety Information

using error checking processes, such as cyclic redundancy checks.

Chapter 3

Getting to Know Your t:slim X2 System

3.1 What your t:slim X2 System Package Includes

Your t:slim X2 System should include the following items:

- 1. t:slim X2 Insulin Pump
- 2. Pump Case
- 3. t:slim X2 System User Guide
- 4. USB Cable
- 5. Wall Power USB Adapter
- 6. AC Power Adapter Plugs
- 7. Cartridge Removal Tool

If you use a Dexcom G5 Mobile CGM, the Sensors and Transmitter are sold and shipped separately by Dexcom.

If any of these items are missing, contact your local Tandem Diabetes Care representative.

Your t:slim X2 Pump is shipped from Tandem Diabetes Care with a clear

screen protector. Do not remove the screen protector.

Your t:slim X2 Pump comes from Tandem Diabetes Care with a protective cover in the place where the cartridge is normally inserted. This cover must be removed and replaced with a cartridge prior to initiating insulin delivery.

Your t:slim X2 Pump also includes consumable components that may require replacement during the life of your pump, including:

- Pump case(s)/clip(s)
- Screen protector
- USB rubber door
- USB cable

Supply Reordering

To order cartridges, infusion sets, supplies, accessories, screen protectors, sensors and transmitter, please contact your local Tandem Diabetes Care representative.

3.2 System Terminology

Pump Terminology

Basal

Basal is a slow continuous delivery of insulin, which keeps BG levels stable between meals and during sleep. It is measured in units per hour (units/hr).

BG

BG is the abbreviation for blood glucose, which is the level of glucose in the blood, measured in mmol/L.

BG Target

BG target is a specific blood glucose goal, an exact number, not a range. When a BG is entered in the t:slim X2 Pump, the calculated insulin bolus will be adjusted up or down as needed to attain this target.

Bolus

A bolus is a quick dose of insulin that is usually delivered to cover food eaten or a high BG. With the t:slim X2 Pump it can be delivered as a Standard, a Correction, an Extended, or a Quick Bolus.

Cannula

The cannula is the part of the infusion set that is inserted under the skin through which insulin is delivered.

Carb

Carbs or Carbohydrates are sugars and starches that the body breaks down into glucose and uses as an energy source, measured in grams.

Carb Ratio

The carb ratio is the number of grams of carbohydrate that 1 unit of insulin will cover. Also known as insulin-to-carbohydrate ratio.

Correction Bolus

A correction bolus is given to correct high BG.

Correction Factor

A correction factor is the amount of BG that is lowered by 1 unit of insulin. Also known as the Insulin Sensitivity Factor (ISF).

Extended Bolus

An extended bolus is a bolus that is delivered over a period of time. It is commonly used to cover food that takes longer to digest. When administering an extended bolus with your t:slim X2 Pump, enter the DELIVER NOW portion to dose a percentage of insulin immediately and the remaining percentage over a period time.

Grams

Grams are a unit of measurement for carbohydrates.

Insulin Duration

Insulin duration is the amount of time that insulin is active and available in the body after a bolus has been delivered. This also relates to the calculation for Insulin on Board (IOB).

Insulin On Board (IOB)

IOB is the insulin that is still active (has the ability to continue to lower the BG) in the body after a bolus has been delivered.

Load

Load refers to the process of removing, filling, and replacing a new cartridge and infusion set.

Personal Profile

A personal profile is a personalized

group of settings that defines the delivery of basal and bolus insulin within specific time segments throughout a 24 hour period.

Quick Bolus

Quick bolus (using the Quick Bolus Button) is a way to deliver a bolus by following beep/vibration commands without navigating through or viewing the t:slim X2 Pump screen.

Temp Rate

Temp rate is an abbreviation for a temporary basal rate. It is used to increase or decrease the current basal rate for a short period of time to accommodate special situations. 100% is the same basal rate as programmed. 120% means 20% more and 80% means 20% less than the programmed basal rate.

Units

Units are the measurement for insulin.

USB Cable

USB is the abbreviation for Universal Serial Bus. The USB cable connects into the t:slim X2 Pump's micro USB port.

CGM Terminology

Alternate Site BG Testing

Alternate site BG testing is when you take a blood glucose value on your meter using a blood sample from an area on your body other than your fingertip. Do not use alternate site testing to calibrate your sensor.

Applicator

The applicator is a disposable piece that comes attached to the sensor pod and inserts the sensor under the skin. There is a needle inside the applicator that is removed after you insert the sensor.

Calibration

Calibration is when you enter blood glucose values from a blood glucose meter into the System. Calibrations are needed for your System to show continuous glucose readings and trend information.

CGM

Continuous glucose monitoring.

Glucose Data Gaps

Glucose data gaps occur when your System is unable to provide a sensor glucose reading.

Glucose Trends

Glucose trends let you see the pattern of your glucose levels. The trend graph shows where your glucose levels have been during the time shown on the screen and where your glucose levels are now.

HypoRepeat

HypoRepeat is an optional alert setting that keeps repeating the fixed low alert every 5 seconds until your sensor glucose value rises above 3.1 mmol/L or you confirm it. This alert can be helpful if you want extra awareness for severe lows.

mmol/L

Millimoles per liter. The standard unit of measure for sensor glucose readings.

Receiver

When the Dexcom G5 Mobile CGM is used with the t:slim X2 Pump to display CGM readings, the insulin pump replaces the receiver. A smartphone

with the Dexcom Mobile App may be used in addition to the t:slim X2 Insulin Pump to receive sensor readings.

Rise and Fall (Rate of Change) Alerts Rise and fall alerts occur based on how much and how fast your glucose levels rise or fall.

RF

RF is the abbreviation for radio frequency. RF transmission is used to send glucose information from the transmitter to the pump.

Safety Lock

The safety lock keeps the needle inside the applicator before you are ready to insert the sensor. It also helps you snap the transmitter out of the sensor pod after your sensor session ends.

Sensor

The sensor is the part that includes an applicator and wire. The applicator inserts the wire under your skin, and the wire measures glucose levels in your tissue fluid.

Sensor Pod

The sensor pod is the small plastic base of the sensor attached to your skin that holds the transmitter in place.

Startup Period

The startup period is the 2-hour period after you tell the System you inserted a new sensor. Sensor glucose readings are not provided during this time.

System Reading

A System reading is a sensor glucose reading shown on your pump. This reading is in mmol/L units and is updated every 5 minutes.

Transmitter

The transmitter is the part that snaps into the sensor pod and wirelessly sends glucose information to your pump.

Transmitter ID

The transmitter ID is a series of numbers and/or letters that you enter into your pump to let it communicate with the transmitter.

Transmitter Latch

The transmitter latch is a small disposable piece that snaps the transmitter into the sensor pod. It is removed after the transmitter is snapped in.

Trend (Rate of Change) Arrows

Trend arrows show how fast your glucose levels are changing. There are 7 different arrows that show when your glucose direction and speed change.

3.3 Explanation of System Symbols

The following are symbols (and their descriptions), which you may find on your t:slim X2 System and/or its packaging. These symbols tell you about the proper and safe use of the system.

Symbol	Meaning
\triangle	Caution; Consult Manual for Important Safety Documentation
	See Instructions for Use
SN	Serial Number of Device
REF	Part Number
IPX7	Watertight Equipment (protected against the effects of temporary immersion in water)
†	Type BF Applied Part (patient isolation, not defibrillator protected)
	Manufacturer
RyOnly	For sale by or on the order of a physician only (U.S.)

Symbol	Meaning
	Use By Date
2	Do Not Re-Use
STERILE R	Sterile by Radiation
PYROGEN	Non Pyrogenic
-20 °C -4 °F	Two-sided Temperature Limits
LOT	Lot Number
((\(\sigma\))	Non-ionizing Radiation
_W	Date of Manufacture

Explanation of System Symbols (Part 1 of 2)

Symbol	Meaning
===	Direct Current (DC) voltage
~	Alternating Current
IP28	Temporary Submersion
IP22	Vertically Falling Drops
%	Two-Sided Humidity Limits
X	European Commission Directive 2006/66/EC
	Electrical Equipment Designed Primarily for Indoor Use
*	Keep Dry

Symbol	Meaning
	Class II Equipment
EC REP	Authorized Representative in the European Community
CE	Marking Certifies That the Device Meets the European Council Directive 93/42/EEC
	Do Not Use if Package is Damaged
SB	Ship By Date
MR	MR Unsafe
\rightarrow	Input

Explanation of System Symbols (Part 2 of 2)

3.4 Explanation of System Icons

The following icons may appear on your t:slim X2 Pump screen:

Symbol	Meaning
80%	The amount of charge left in the pump battery.
Y	CGM sensor session is active, and the transmitter is communicating with the pump.
Y	CGM sensor session is active, but the transmitter is not communicating with the pump.
	CGM calibration is required.
Ţ	A system reminder, alert, error, or alarm is active.
В	Basal insulin is programmed and being delivered.
~	Accept button. Tap to continue to the next screen or to answer yes to a message on the pump screen.
~	Save button. Tap to save settings on the screen.
-	New button. Tap to add a new item.

Symbol	Meaning
235 u	Amount of insulin remaining in the reservoir.
Т	A temporary basal rate is active.
0	A basal rate of 0 u/hr is active.
T	A temporary basal rate of 0 u/hr is active.
	A bolus is being delivered.
<u> </u>	All insulin deliveries are stopped.
×	Cancel button. Tap to cancel the current operation.
×	Decline button. Tap to exit the screen or answer no to a message on the pump screen.
—	Back button. Tap to navigate to the previous screen.

Pump Icon Definitions (Part 1 of 2)

Symbol	Meaning
×	Delete button. Tap to delete characters or digits on a keypad.
	The associated setting is turned ON.
②	CGM sensor session is active, but the transmitter is not communicating with the pump.
	The CGM sensor has failed. See Section 24.4 Additional CGM Status Symbols.
	The CGM sensor session has ended. See Section 24.4 Additional CGM Status Symbols.
Z •	Wait 15 minutes calibration error.
	Startup calibration is required (2 BG values).
	Additional startup calibration is required
	CGM calibration is required.
	Space button. Tap to enter a space on the character keypad.

Symbol	Meaning
=	Total button. Tap to total values on a keypad.
	The associated setting is turned OFF.
mmol/L	Unknown sensor reading.
oexcom)	Transmitter error.
	Sensor startup 0–30 minutes.
	Sensor startup 31–60 minutes.
	Sensor startup 61–90 minutes.
	Sensor startup 91–119 minutes.
**	Bluetooth®.

Pump Icon Definitions (Part 2 of 2)

3.5 Explanation of System Colors



Red LED

1 red blink every 30 seconds indicates a malfunction or alarm condition.



Yellow LED

1 yellow blink every 30 seconds indicates an alert or reminder condition.



Green LED

- 1 green blink every 30 seconds indicates the pump is functioning normally.
- 3 green blinks every 30 seconds indicate the pump is charging.



Orange Highlight

When editing settings, changes are highlighted in orange for review before saving.

THIS PAGE IS INTENTIONALLY LEFT BLANK

3.6 Screen Lock

- 1. Time and Date Display: Displays the current time and date.
- Alert Icon: Indicates a reminder, alert or alarm is active behind the lock screen.
- Battery Level: Displays the level of battery power remaining. When connected for charging, the charging icon (lightning bolt) will display.
- 4. 1-2-3: Unlocks pump screen.
- Insulin On Board (IOB): Amount and time remaining of any active insulin on board.
- 6. Active Bolus Icon: Indicates a bolus is active.
- 7. Status: Displays current system settings and insulin delivery status.
- 8. **Insulin Level**: Displays the current amount of insulin in the cartridge.


9. **Tandem Logo:** Returns to the Home Screen.



3.7 CGM Screen Lock

- 1. Time and Date Display: Displays the current time and date.
- Antenna: Indicates communication status between pump and transmitter.
- Battery Level: Displays the level of battery power remaining. When connected for charging, the charging icon (lightning bolt) will display.
- 4. High Glucose Alert Setting.
- Glucose Target Range.
- Low Glucose Alert Setting.
- Plot of Most Recent Sensor Glucose Readings.
- 8. 1-2-3: Unlocks pump screen.
- 9. Active Bolus Icon: Indicates a bolus is active.


- 10. Status: Displays current system settings and insulin delivery status.
- 11. **Insulin Level**: Displays the current amount of insulin in the cartridge.
- 12. Most Recent 5-Minute Glucose Reading.
- 13. **Trend Arrow:** Indicates direction and rate of change.
- 14. Trend Graph Time (HRS): 1, 3, 6, 12 and 24 hour views available.
- Insulin On Board (IOB): Amount and time remaining of any active insulin on board.

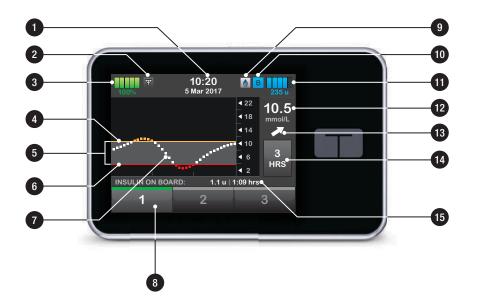

3.8 Home Screen

- Battery Level: Displays the level of battery power remaining. When connected for charging, the charging icon (lightning bolt) will display.
- USB Port: Port to charge your t:slim X2 Pump battery. Close the cover when not in use.
- 3. Bolus: Program and deliver a bolus.
- Options: Stop/Resume insulin delivery, manage Pump and CGM Settings, program a Temp Rate, Load cartridge, and view History.
- Insulin On Board (IOB): Amount and time remaining of any active insulin on board.
- 6. Time and Date Display: Displays the current time and date.
- 7. **Status:** Displays current system settings and insulin delivery status.

- 8. **Insulin Level:** Displays the current amount of insulin in the cartridge.
- 9. **Tandem Logo:** Returns to the Home Screen.
- 10. Cartridge Tubing: Tubing that is attached to the cartridge.
- 11. **Tubing Connector:** Connects the cartridge tubing to the infusion set tubing.
- 12. Screen On/Quick Bolus Button:
 Turns the t:slim X2 Pump screen
 on/off or programs a Quick Bolus (if
 activated).
- LED Indicator: Illuminates when connected to a power supply and indicates proper functionality.

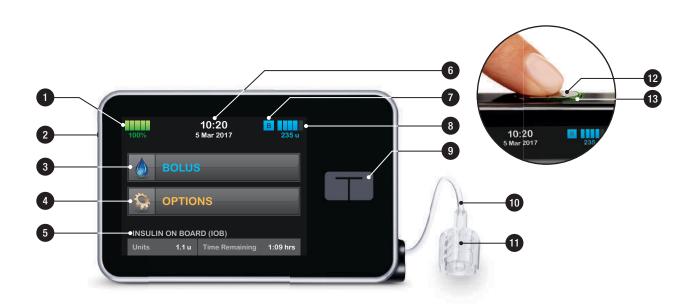

3.9 CGM Home Screen

- 1. Time and Date Display: Displays the current time and date.
- Antenna: Indicates communication status between pump and transmitter.
- Battery Level: Displays the level of battery power remaining. When connected for charging, the charging icon (lightning bolt) will display.
- High Glucose Alert Setting.
- Glucose Target Range.
- 6. Low Glucose Alert Setting.
- Plot of Most Recent Sensor Glucose Readings.
- Options: Stop/Resume insulin delivery, manage Pump and CGM Settings, program a Temp Rate, Load cartridge, and view History.
- 9. Bolus: Program and deliver a bolus.

- 10. Status: Displays current system settings and insulin delivery status.
- 11. **Insulin Level:** Displays the current amount of insulin in the cartridge.
- 12. Most Recent 5-Minute Glucose Reading.
- 13. **Trend Arrow:** Indicates direction and rate of change.
- 14. Trend Graph Time (HRS): 1, 3, 6, 12 and 24 hour views available.
- Insulin On Board (IOB): Amount and time remaining of any active insulin on board.


To view CGM information on the full screen:

From the Home Screen tap anywhere on the CGM trend graph.



Tap the "minimize" icon to return to the Home Screen.



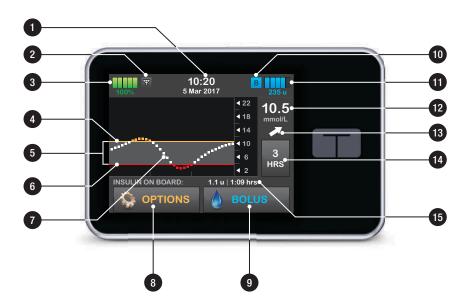

3.10 Status Screen

The status screen can be accessed from the lock screen and the Home Screen. It is for display only; no changes can be made from this screen.

- 1. **Profile:** Displays current active Personal Profile.
- Basal Rate: Displays current basal rate being delivered. (If a Temp Rate is active, it is displayed in units/hr.)
- 3. Last Bolus: Displays the amount, date and time of last bolus.
- 4. Carbohydrates: Indicates whether feature is on or off.
- 5. **Up/Down Arrow:** Indicates there is more information.
- Correction Factor: Displays current correction factor used to calculate a bolus.
- 7. Carb Ratio: Displays current carb ratio used to calculate a bolus.

- 8. Target BG: Displays current BG target used to calculate a bolus.
- Insulin Duration: Displays current insulin duration setting used to calculate insulin on board.
- 10. Last Calibration: Displays date and time of last calibration.
- 11. Time Sensor Started: Displays date and time of last time sensor started.
- 12. Transmitter Battery: Displays transmitter battery status.
- 13. Mobile Connection: Displays whether the mobile connection is turned on or off, whether a mobile device is paired with the pump, and if so whether the device is actively connected to the pump.


The Mobile App may not yet be available in your area.



3.11 Bolus Screen

- 1. Returns to the Home Screen.
- 2. Carbs: Enter grams of carb.
- Units: Displays total units calculated. Tap to enter a bolus request or change (override) a calculated bolus.
- View Calculation: Displays how the insulin dose was calculated using the current settings.
- Add BG: Enter blood glucose level.
 If you have an active CGM session and your CGM value and trend arrow are available, this value is auto-populated.
- 6. Moves to next step.
- 7. Returns to the Home Screen.
- 8. Insulin: Enter units of insulin.
- 9. **Units:** Displays total units calculated. Tap to enter a bolus


- request or change (override) a calculated bolus.
- View Calculation: Displays how the insulin dose was calculated using the current settings.
- 11. Add BG: Enter blood glucose level.
- 12. Moves to next step.

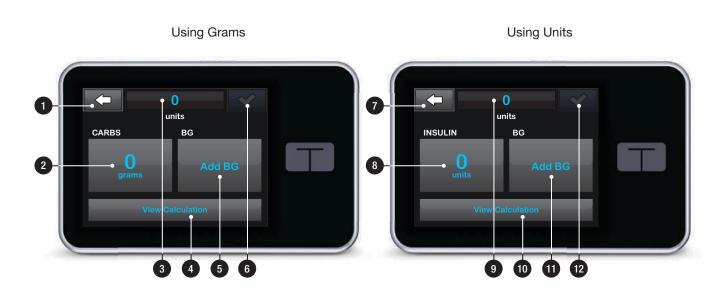

3.12 Options Screen

- 1. Returns to the Home Screen.
- Stop Insulin: Stops insulin delivery.
 If insulin delivery is stopped,
 RESUME INSULIN will be displayed.
- My Pump: Personal Profiles, Alerts & Reminders, Pump Settings, and Pump Info.
- 4. My CGM: Start/Stop Sensor, Calibrate CGM, CGM Alerts, and CGM Settings.
- 5. **Temp Rate:** Programs a temporary basal rate.
- 6. **Up/Down Arrow:** Indicates there is more information.
- 7. Load: Change Cartridge, Fill Tubing, Fill Cannula, and Site Reminder.
- 8. Pump History: Displays historical log of pump events.

CGM History: Displays historical log of CGM events.

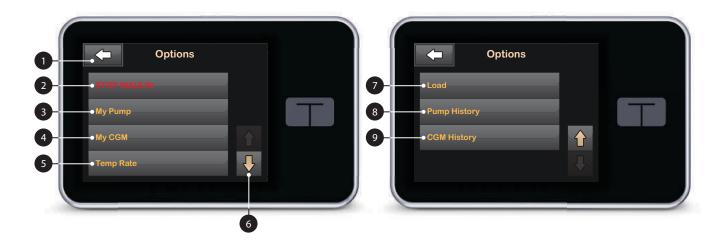

3.13 My Pump Screen

- Personal Profiles: A group of settings that define basal and bolus delivery.
- 2. Alerts & Reminders: Customize Pump Alerts and Reminders.
- 3. Pump Settings: Customize Quick Bolus, Display Settings, Bluetooth Settings, Sound Volume, turn Feature Lock on or off, and set Time and Date.
- Pump Info: Displays t:slim X2 Pump serial number and other technical information.




3.14 My CGM Screen

- Start Sensor: Starts a CGM session. If sensor is active, STOP SENSOR will be displayed.
- 2. Calibrate CGM: Enter a calibration blood glucose value. Only active when sensor session is active.
- 3. **CGM Alerts:** Customize CGM Alerts.
- 4. Transmitter ID: Enter the transmitter ID.
- 5. **CGM Info:** View the CGM information.



3.15 Number Keypad Screen

- 1. Value Entered.
- 2. Returns to previous screen.
- 3. Keypad Numbers.
- 4. Allows numbers to be added on gram screen. If in units, this displays as a decimal point.
- 5. Completes task and saves information entered.
- 6. **Units/Grams:** Value of what is entered.

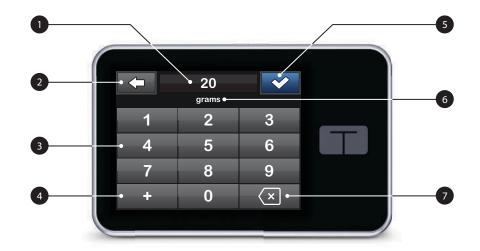

3.16 Letter Keypad Screen

- 1. Name of Profile.
- 2. Returns to previous screen.
- 3. Space: Enters a space.
- 4. 123: Changes keypad mode from letters (ABC) to numbers (123).
- 5. Saves entered information.
- Letters: Tap once for first letter displayed, 2 quick taps for middle letter, and 3 quick taps for third letter.
- 7. Deletes last letter or number entered.



Section 2

Key Pump Features

Chapter 4

Getting Started

4.1 Charging the t:slim X2 Pump

The t:slim X2 Pump is powered by an internal lithium polymer rechargeable battery. A full charge will last up to 7 days with average use, or 5 days when using CGM with average use. Please be aware that the battery life on a single charge can vary considerably depending on individual usage, including insulin delivered, display-on time, and frequency of reminders, alerts and alarms.

Accessories for charging from wall and automobile outlets, as well as from a PC USB port are included with the pump. Use only the accessories provided with the System to charge your t:slim X2 Pump. If you lose any of the accessories, or need a replacement, contact your local Tandem Diabetes Care representative.

The battery level indicator is displayed in the upper left portion of the Home Screen. The charge amount will increase or decrease by 5% at a time (for example, you will see 100%, 95%, 90%, 85%). When the charge amount is less than 5%, it will begin decreasing

1% at a time (for example, you will see 4%, 3%, 2%, 1%).

When you first receive your t:slim X2 Pump, you will need to connect it to a charging source before it can be used. Charge the pump until the battery level indicator on the upper left portion of the Home Screen reads 100% (initial charge can take up to 2.5 hours).

The t:slim X2 Pump continues to operate normally while charging. You do not need to disconnect from the pump while charging.

A PRECAUTION

MAKE SURE to not move further than the length of the USB cable when you are connected to the pump and to a charging source. Moving further than the length of the USB cable may cause the cannula to be pulled out of the infusion site. For this reason it is recommended not to charge the pump while sleeping.

If you choose to disconnect from the pump while charging, check with your healthcare provider for specific guidelines. Depending on the length of time you are disconnected, you may need to replace missed basal and/or bolus

insulin. Check your blood glucose before disconnecting from the pump and again when you reconnect.

To charge the pump from an AC Power Outlet:

- 1. Plug the included USB cable into the AC power adapter.
- 2. Plug the AC power adapter into a grounded AC power outlet.
- 3. Plug the other end of the cable into the micro USB port on the pump.

To charge the pump using the optional Car Power USB Adapter:

- 1. Plug the USB cable into the Car Power USB Adapter.
- Plug the Car Power USB Adapter into a grounded auxiliary power outlet.
- 3. Plug the other end of the cable into the micro USB port on the pump.

A WARNING

When using an optional Car Power USB Adapter, the charger must be connected to an isolated, battery powered 12 Volt system, such as an automobile. Connecting the DC vehicle adapter charger to 12 Volt DC that is generated by a power supply connected to alternating current (AC) mains is prohibited.

To charge the pump using a USB port on a Personal Computer (PC):

Ensure that the PC complies with the IEC 60950-1 (or equivalent) safety standard.

- 1. Plug the included USB cable into your computer.
- 2. Plug the other end of the cable into the micro USB port on the pump.

Depending on your computer, charging time will vary. The pump will display a CONNECTION ERROR ALERT message if it is not charging properly.

When you charge the t:slim X2 Pump, you will notice the following:

• The screen illuminates

- An audible alert
- The LED (edge around the Screen On/Quick Bolus Button) blinks green
- A vibrating alert.
- A charge symbol (lightning bolt) on the battery level indicator appears

A PRECAUTION

CONFIRM that the screen display turns on, you can hear audible beeps, feel the pump vibrate, and see the green LED light blinking around the edge of the Screen On Button when you connect a power source to the USB port. These features are used to notify you about alerts, alarms, and other conditions that require your attention. If these features are not working, discontinue use of the System and contact your local Tandem Diabetes Care representative.

Charging Tips

Tandem Diabetes Care recommends periodically checking the battery level indicator, charging the pump for a short period of time every day (10 to 15 minutes), and also avoiding frequent full discharges.

■ NOTE: Fully Discharged Battery

If the battery is fully discharged, the screen may not power on immediately when connected to a charging source. The LED around the Screen On/Quick Bolus Button will blink green until there is enough charge to power on the touch screen.

4.2 Using the Touch Screen

To turn on your t:slim X2 Pump screen, first press the Screen On/Quick Bolus Button, then use the pad of your finger to quickly and lightly tap on the screen. Do not use your finger nail or other object to interact with the screen. It will not activate the screen or its functions.

Your t:slim X2 Pump is designed to give you quick and easy access to the functions that you will use in your day-to-day diabetes management whether basic or advanced.

The t:slim X2 Pump has several safety features to prevent unintentional interaction with the touch screen. The screen must be unlocked by tapping 1–2–3 in sequence. On all screens, if three non-active areas of the touch screen

Chapter 4 – Getting Started

are tapped before an active area is tapped, the screen will turn off to prevent accidental button presses.

■ NOTE: Touch Screen Tips

When using the t:slim X2 Pump, tap the Tandem Logo to return to the Home Screen or tap
to return to the previous screen.

4.3 Turning on the t:slim X2 Pump Screen

To turn on your t:slim X2 Pump screen, press the Screen On/Quick Bolus Button, located on the top of the pump, once.

The Screen Lock screen will be displayed.

NOTE: Turning off the Pump Screen

Turn off the pump screen by pressing the Screen On/Quick Bolus button before placing the pump back in its case or any pocket/clothing. Always position the pump screen away from the skin when worn under clothing.

The pump continues to function normally when the screen is not on.

4.4 Unlocking the t:slim X2 Pump Screen

The screen lock screen appears anytime you turn on the screen, and after a bolus or temp rate is requested. To unlock the screen:

- 1. Press Screen On/Quick Bolus Button.
- 2. Tap 1.
- Tap 2.
- 4. Tap 3.
- The pump screen is now unlocked. The last screen that was viewed will be displayed.

You must tap 1–2–3 in sequential order to unlock the pump. If you do not tap 1–2–3 in sequential order, the pump will force you to restart the unlock sequence from the beginning.

4.5 Selecting Your Language

The language selection screen displays when you unlock the pump screen for the first time, or when you unlock the pump screen after turning the pump off.

To select your language:

 Tap the circle next to the language you want to display. Tap the down or up arrow to see additional language selections.



2. Tap to save the selection and continue with pump setup.

4.6 Editing Time

After powering up your t:slim X2 Pump for the first time, set the current time and date. Refer back to this section if you need to edit the time for either traveling in a different time zone or adjusting for Daylight Savings Time.

A PRECAUTION

ALWAYS make sure that the correct time and date are set on your System. When using AM/ PM, always check that the AM/PM setting is accurate. Not having the correct time and date setting may affect safe insulin delivery.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- 3. Tap Pump Settings.
- 4. Tap Time and Date.
- Tap Edit Time.
- 6. Tap Time.

- Using the onscreen keypad, enter the hour and minutes. Verify and tap
- 3. Tap **Time of Day** to set AM or PM. Verify the correct time is set and tap

Any edits to Time or Date will not be saved until you tap .

1.7 Editing Date

- 1. From the Time and Date screen tap Edit Date.
- 2. Tap Month.
- Find and tap the current month displayed on the right. Use Up/ Down Arrow to view months not displayed.
- 4. Tap Day. Using the onscreen keypad enter the current day. Verify and tap ...
- 5. Tap Year.

- Using the onscreen keypad enter the current year. Verify and tap
- 7. Verify the correct date is set and tap
- 8. Tap **Tandem Logo** to return to the Home Screen.

.8 Turning the Pump Off

To turn the pump off completely for storage, shipping, or for reset:

- 1. Plug the pump into a power source.
- 2. Press and hold the Screen On/ Quick Bolus Button for 30 seconds.

THIS PAGE IS INTENTIONALLY LEFT BLANK

Chapter 5

Infusion Site Care and Loading Cartridge

Chapter 5 - Infusion Site Care and Loading Cartridge

5.1 Infusion Site Selection and Care

A WARNING

D0 use only cartridges and insulin infusion sets with matching connectors and follow their instructions for use. Failure to do so may result in over delivery or under delivery of insulin and may cause very low or very high blood glucose.

▲ WARNING

DO NOT ignore infusion set cannula fractures. Infusion set cannulas may fracture on rare occasions. If an infusion set cannula breaks and no portion of it is visible above the skin, do not attempt to remove it. Seek professional medical help if you have symptoms of infection or inflammation—redness, swelling or pain—at the insertion site. If you experience a broken infusion set cannula, please report this to your local Tandem Diabetes Care representative.

A WARNING

DO NOT place your infusion set on any scars, lumps, moles, stretch marks or tattoos. Placing your infusion set in these areas can cause swelling, irritation or infection. This can affect insulin absorption and cause high or low blood glucose.

A PRECAUTION

CHECK your infusion site daily for proper placement and leaks. REPLACE your infusion set if you notice leaks around the site. Improperly placed sites or leaks around the infusion site can result in under delivery of insulin.

General Guidelines

Site Selection

- Your infusion set can be worn anywhere on your body that you would normally inject insulin.
 Absorption varies from site to site.
 Discuss options with your healthcare provider.
- The most commonly used sites are the abdomen, upper buttocks, hips, upper arms, and upper legs.
- The abdomen is the most popular site because of access to fatty tissue. If using the abdominal area, AVOID:
 - Areas that would constrict the site such as the belt line, waistline, or where you would normally bend.

- Areas 5 centimeters (2 inches) around your belly button.
- Any scars, moles, stretch marks, or tattoos.
- Areas within 8 centimeters (3 inches) of your CGM sensor site.

Site Rotation

A PRECAUTION

CHANGE your infusion set every 48–72 hours as recommended by your healthcare provider. Wash your hands with anti-bacterial soap before handling the infusion set and thoroughly clean the insertion site on your body to avoid infection. Contact your healthcare provider if you have symptoms of infection at your insulin infusion site.

- The infusion set must be replaced and rotated every 48–72 hours, or more often if needed.
- With experience, you will find areas that not only provide better absorption, but are more comfortable. Keep in mind, using the same areas may cause scarring

or lumps which can affect insulin absorption.


 Consult your healthcare provider to establish a rotation schedule that best fits your needs.

Keep it Clean

- When changing your infusion set, use clean techniques to avoid an infection.
- Wash your hands, use antiseptic wipes or infusion site preparation products, and keep the area clean.
- Site preparation products that have both an antiseptic and an adhesive are encouraged.

5.2 Cartridge Instructions for Use

For complete cartridge labeling, consult the Cartridge Instructions for Use included in the t:slim™ Cartridge box.

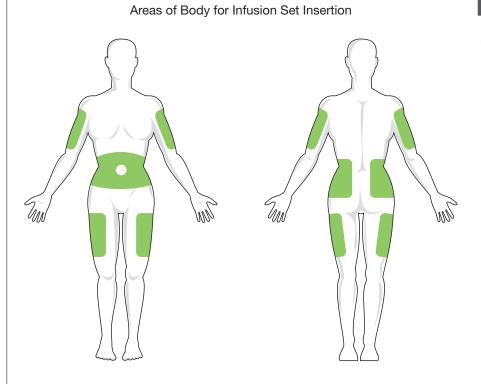

Chapter 5 – Infusion Site Care and Loading Cartridge

5.3 Filling and Loading a t:slim Cartridge

This section describes how to fill the cartridge with insulin and load the cartridge into your t:slim X2 Pump. The single-use disposable cartridge can hold up to 300 units (3.0 mL) of insulin.

A WARNING

DO NOT use any other insulin with your System other than U-100 Humalog or NovoLog/NovoRapid. Only Humalog and NovoLog/NovoRapid have been tested and found to be compatible for use in the System. Use of insulin with lesser or greater concentration can result in under delivery or over delivery of insulin. This can cause very high or a very low blood glucose.

A WARNING

DO NOT reuse cartridges or use cartridges other than those manufactured by Tandem Diabetes Care. Use of cartridges not manufactured by Tandem Diabetes Care or reuse of cartridges may result in over delivery or under delivery of insulin. This can cause very low or very high blood glucose.

Before you begin, make sure you have the following items:

- 1 unopened cartridge
- 3.0 mL syringe and fill needle
- Vial of Humalog or NovoLog/ NovoRapid insulin
- Alcohol prep swab
- 1 new infusion set
- Infusion set Instructions for Use


■ NOTE: Fill Tubing Beep Volume

The pump will beep or vibrate, depending on your pump settings, while the cartridge is filling with insulin. To change the Fill Tubing sound setting, see 11.3 t:slim X2 Sound Volume.

■ NOTE: Removing the Cartridge

DO NOT remove the used cartridge from the pump during the load process until prompted on the t:slim X2 Pump screen.

The illustration identifies the connector and insulin fill port used in the cartridge filling process.

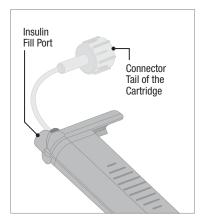

A PRECAUTION

CHANGE your cartridge every 48–72 hours as recommended by your healthcare provider. Wash your hands with anti-bacterial soap before handling the infusion set and thoroughly clean the insertion site on your body to avoid infection. Contact your healthcare provider if you have symptoms of infection at your insulin infusion site.

Instructions for Drawing Insulin from Vial into Syringe


A PRECAUTION

ALWAYS remove all air bubbles from the System before beginning insulin delivery. Ensure there are no air bubbles when drawing insulin into the filling syringe, hold the pump with the white fill port pointed up when filling the tubing, and ensure that there are no air bubbles in the tubing when filling. Air in the system takes space where insulin should be and can affect insulin delivery.

The fill estimate displayed on the pump is the amount of insulin available for delivery. It does not include insulin needed to fill the tubing (up to 30 units) and a small amount of insulin that is not available for delivery. When filling the syringe, add approximately 45 units to the amount of insulin you want available for delivery.


For example, the pump requires a minimum of 50 units available for delivery after fill tubing has been completed. Fill the syringe with approximately 95 units to have enough to fill your tubing and still have 50 units available for delivery.

- Inspect the needle and syringe package for any signs of damage. Discard any damaged product.
- 2. Wash your hands thoroughly.
- 3. Wipe the rubber septum of the insulin vial with an alcohol swab.
- Remove the needle and syringe from their packaging. Securely twist needle onto syringe. Safely remove protective cap from needle by pulling outward.
- 5. Draw air into syringe up to the amount of insulin desired.

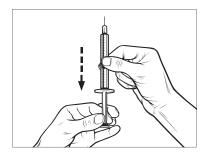

6. With insulin vial upright, insert needle into vial. Inject air from


syringe into vial. Maintain on syringe plunger.

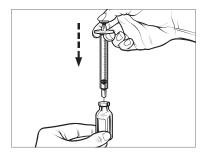


7. With needle still inserted into vial, turn vial and syringe upside down. Release syringe plunger. Insulin will begin to flow from the vial into the syringe.

Chapter 5 – Infusion Site Care and Loading Cartridge

8. Slowly pull back the plunger to the desired amount of insulin.

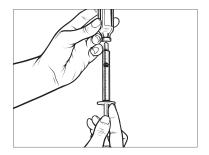


 While the filling needle is still in the vial and upside down, tap the syringe so that any air bubbles rise to the top. Then slowly push the plunger upwards, forcing any air bubbles back into the vial.

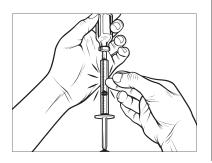


- 10. Check the syringe for air bubbles and do one of the following:
 - If there are air bubbles present, repeat step 9.
 - If no air bubbles are present, remove the filling needle from the vial.

Instructions for Filling the Cartridge


- Inspect the cartridge package for any signs of damage. Discard any damaged product.
- 2. Open the package and remove the cartridge.
- Hold the cartridge upright and gently insert the needle into the white insulin fill port on the cartridge. The needle is not intended to go all the way in, so do not force it.

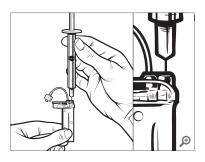

4. Keeping the syringe vertically aligned with the cartridge, and the needle inside the fill port, pull back on the plunger until it is fully retracted. This will remove any residual air from the cartridge. Bubbles will rise toward the plunger.

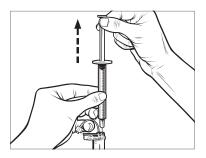


 Make sure the needle is still in the fill port and release the plunger.
 Pressure will pull the plunger to its neutral position but it will NOT push any air back inside the cartridge.

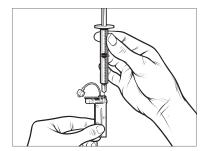


- 6. Withdraw the needle from the fill port.
- 7. Turn the syringe upright and pull down on the plunger. Flick the


barrel to make sure that any air bubbles rise to the top.

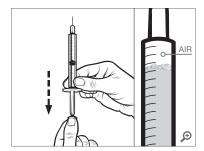

8. Gently press on the plunger to remove air bubbles until insulin fills the needle hub and you see a drop of insulin at the tip of the needle.

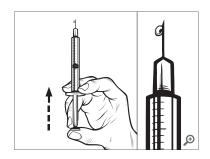

 Re-insert the needle in the fill port and slowly fill the cartridge with insulin. It is normal to feel some back pressure as you slowly press on the plunger.

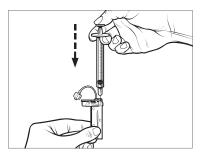

- 10. Maintain on the plunger while you remove the needle from the cartridge. Check the cartridge for leaks. If you detect insulin leaking, discard the cartridge and repeat entire process with a new cartridge.
- 11. Consult your local waste management authorities for proper disposal instructions of used needles, syringes, cartridges, and infusion sets.

Chapter 5 - Infusion Site Care and Loading Cartridge

Instructions on How to Install a Cartridge

If this is the very first time you are loading the cartridge, remove the shipping canister (which is not for human use) from the back of the pump.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap the Down Arrow.
- 3. Tap Load.
- During the load sequence, the Tandem Logo is disabled.

Tapping it will not return to the Home Screen.

- 4. Tap Change Cartridge.
- Screen will display that all insulin deliveries will be stopped. Tap to continue.
 - **NOTE: First Time Use**

This screen will not be displayed if this is the first time loading a new cartridge and you have not started actively pumping.

- 6. Disconnect the infusion set from your body and tap to continue.
- "Preparing for Cartridge" screen is displayed.
- Remove the used cartridge. If needed, place the cartridge removal tool or the edge of a coin in the slot at the bottom of the cartridge and twist to aid in the removal of the cartridge.
- Place bottom of the cartridge at the end of the pump. Make sure cartridge is lined up to both guide tracks.



9. Push on the circular fill port next to the cartridge tubing to slide the

cartridge onto the pump. Tap UNLOCK when completed.



- 10. Tap to continue.
- Detecting Cartridge screen is displayed.
- After completing the cartridge change, the pump will automatically prompt you to fill the tubing.

A WARNING

DO NOT remove or add insulin from a filled cartridge after loading onto the pump. This will result in an inaccurate display of the insulin level on the Home Screen and you could run out of insulin before the pump detects an empty cartridge. This can cause very high blood glucose, or Diabetic Ketoacidosis (DKA).

5.4 Filling Tubing

Filling the Infusion Set Tubing with Insulin

This section describes how to fill the infusion set tubing with insulin after you change the cartridge.

To fill the tubing without changing the cartridge, from the Home Screen tap OPTIONS, tap Down Arrow, tap Load, tap Fill Tubing and then follow the instructions.

- Tap NEW if you installed a new cartridge.
- Tap FILL if you did not install a new cartridge and want to continue with filling the tubing.

A WARNING

NEVER fill your tubing while your infusion set is connected to your body. Always ensure that the infusion set is disconnected from your body before filling the tubing. Failure to disconnect your infusion set from your body before filling the tubing can result in over delivery of insulin.

This can cause serious injury or death from very low blood glucose.

A PRECAUTION


ALWAYS remove all air bubbles from the System before beginning insulin delivery. Ensure there are no air bubbles when drawing insulin into the filling syringe, hold the pump with the white fill port pointed up when filling the tubing, and ensure that there are no air bubbles in the tubing when filling. Air in the system takes space where insulin should be and can affect insulin delivery.

A PRECAUTION

CHECK your infusion set tubing daily for any leaks, air bubbles, or kinks. Air in the tubing, leaks in the tubing, or kinked tubing may restrict or stop insulin delivery and result in under delivery of insulin.

- 1. Verify that the infusion set is disconnected from your body.
- 2. To begin connecting the infusion set tubing to the tubing connector on the cartridge, ensure that the new infusion set package is not damaged, and remove the sterile tubing from the package. If the

- package is damaged or opened, discard of properly and use another tubing set.
- Remove the infusion set tubing cap from the tubing connector. Be careful to keep the tubing connector away from unclean areas.
- Attach the infusion set tubing to the tubing connector on the cartridge tubing. Twist clockwise until finger tight and then twist another quarter of a turn to ensure a secure connection.

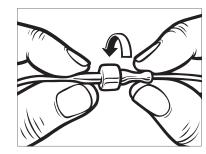


AWARNING

ALWAYS twist the tubing connector between the cartridge tubing and the infu-

Chapter 5 – Infusion Site Care and Loading Cartridge

sion set tubing an extra quarter of a turn to ensure a secure connection. A loose connection can cause insulin to leak, resulting in under delivery of insulin. This can cause high blood glucose.

- 5. Tap 🕶
- Hold the pump vertically to ensure any air in the cartridge will be dispelled first. Tap START. The pump will beep and vibrate regularly while the tubing is being filled.

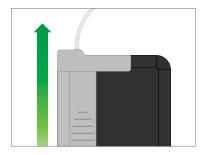

√ "Starting Fill" screen is displayed.

The following are approximate amounts of insulin to fill different tubing lengths:

- 15-20 units for 60 cm tubing
- 20-25 units for 80 cm tubing
- 25–30 units for 110 cm tubing
- 7. Tap STOP after 3 drops of insulin are seen at the end of the infusion set tubing.
- ✓ Stopping Fill screen is displayed.
- Detecting Insulin screen is displayed.
- 8. Verify that drops are seen and tap DONE.
- If you do not see drops, tap FILL.
 The Fill Tubing screen appears,
 repeat steps 3 to 5 until you see 3 drops of insulin at the end of the tubing.

If you did not tap STOP, the "Max fill amount reached!" screen will appear. Do one of the following:

 If you are finished filling the tubing, tap DONE.

- If you want to fill the tubing with more than 30 units, tap FILL to go back to the Fill Tubing screen.
- ✓ Fill Tubing is complete.

■ NOTE: Initial Display of Insulin

After tubing fill is complete, when the pump returns to the Home Screen, an estimate of how much insulin is in the cartridge is displayed in the upper right portion of the screen. You will see one of the following on the screen:

+ 40 u	More than 40 units detected
	in the cartridge
+ 60 u	More than 60 units detected
	in the cartridge
+ 120 u	More than 120 units detected
	in the cartridge
+ 180 u	More than 180 units detected
	in the cartridge
+ 240 u	More than 240 units detected
	in the cartridge

After 10 units are delivered, an actual number of units remaining in the cartridge will be displayed on the Home Screen

The amount of insulin remaining displayed on the Home Screen will decrease 5 units at a time (for example, you will see 140, 135, 130, 125). When less than 40 units remain, it will begin decreasing 1 unit at a time (for example, you will see 40, 39, 38, 37) until there is 1 unit remaining.

5.5 Filling Cannula

Filling the Infusion Set Cannula with Insulin

This section describes how to fill the infusion set cannula with insulin after you fill the tubing.

To fill the cannula without filling the tubing, from the Home Screen, tap OPTIONS, tap Down Arrow, tap Load, tap Fill Cannula and then follow the instructions below.

If you are using a steel needle infusion set, there is no cannula; skip this section.

To Fill the Cannula:

1. Tap Fill Cannula.

- Insert a new infusion set and connect filled tubing to site, then tap
- 3. Tap Edit Fill Amount.
- ✓ The cannula fill amount displayed is based on your last cannula fill amount. Filling stops at this amount.
- 4. Select amount needed for cannula fill.
 - Refer to your infusion set instructions for use for proper cannula fill amount.
 - If the amount needed is not listed, tap Other amount and use the onscreen keypad to enter a value between 0.1 to 1.0 unit.
- 5. Tap START.
- ✓ "Starting Fill" screen is displayed.
- ✓ After fill is complete, "Stopping Fill" screen is displayed.

■ NOTE: Stopping Fill

You can tap STOP at any time during the fill process if you want to stop filling the cannula.

- The screen will return to the Load menu if the Site Reminder is turned off.
- Tap to resume insulin if finished. Or tap Site Reminder to set reminder. If Site Reminder is on, the pump will automatically display the Site Reminder Screen (refer to next section).

5.6 Setting Site Reminder

This section describes how to set the Site Reminder after you fill the cannula.

To set the Site Reminder without filling the cannula, from the Home Screen, tap OPTIONS, tap Down Arrow, tap Load, tap Site Reminder then follow the instructions below.

A PRECAUTION

DO NOT change your infusion set before bedtime or if you will not be able to test your blood

Chapter 5 – Infusion Site Care and Loading Cartridge

glucose 1–2 hours after the new infusion set is placed. It is important to confirm that the infusion set is inserted correctly and delivering insulin. It is also important to respond quickly to any problems with the insertion to ensure continued insulin delivery.

- Tap if correct. Tap Edit
 Reminder if settings need to be changed.
- 2. Tap Remind Me In and select the number of days (1–3).
- ✓ The default for the Site Reminder is set for 3 days
- 3. Tap Remind Me At. Use the onscreen keypad to enter time and tap ...
- Tap Time of Day to change AM or PM. Tap ... Verify Site Reminder is set correctly and tap ...
- ✓ "Setting Saved" screen is displayed.
- ✓ Load screen is displayed.
- Tap ______.

- ✓ A reminder to test BG in 1 to 2 hours will display.
- 6. Tap ____.

■ NOTE: First Time Use

If this is the first time using your t:slim X2 Pump and a Personal Profile has not been defined, a screen will notify you that a profile must be activated to resume insulin.

Tap ok.

 RESUMING INSULIN screen is temporarily displayed.

Chapter 6

Personal Profiles

6.1 Personal Profiles Overview

WARNING

DO NOT start to use your System before consulting with your healthcare provider to determine which features are most appropriate for you. Only your healthcare provider can determine and help you adjust your Basal Rate(s), Carb Ratio(s), Correction Factor(s), Target BG, and duration of insulin action. In addition, only your healthcare provider can determine your CGM settings and how you should use your sensor trend information to help you manage your diabetes. Incorrect settings can result in over delivery or under delivery of insulin. This can cause very low or very high blood glucose.

A Personal Profile is a group of settings that define basal and bolus delivery within specific time segments throughout a 24-hour period. Each profile can be personalized with a name. Within a Personal Profile the following can be set:

 Timed Settings: Basal Rate, Correction Factor, Carb Ratio and Target BG. Bolus Settings: Insulin Duration, Max Bolus and Carbohydrates setting (on/off).

The t:slim X2 Pump uses the settings in your active profile to calculate the delivery of basal insulin, food boluses and correction boluses based on your Target BG. If you only define a basal rate in Timed Settings, your pump will only be able to deliver basal insulin and standard and extended boluses. Your pump will not calculate correction boluses.

Up to 6 different Personal Profiles can be created and up to 16 different time segments can be set in each Personal Profile. Having several Personal Profiles provides more flexibility for your body and lifestyle. For example, you could have "Weekday" and "Weekend" profiles if you have different insulin delivery needs on weekdays and weekends, based on schedule, food intake, activity, and exercise, etc.

6.2 Creating a New Profile

Creating Personal Profiles

You can create up to 6 Personal Profiles; however, only 1 can be active at a time. In the Personal Profiles screen, the active profile is positioned at the top of the list and is marked as ON. When you create a Personal Profile, you can set any or all of the following Timed Settings:

- Basal Rate (your basal rate in units/ hr)
- Correction Factor (amount 1 unit of insulin lowers BG)
- Carb Ratio (grams of carbs covered by 1 unit of insulin)
- Target BG (your ideal BG level, measured in mmol/L)

Although you do not need to define every setting, some pump features require certain settings to be defined and activated. When you are creating a new profile, your pump prompts you to

set up any required settings before you can continue.

The ranges you can set for Timed Settings are:

- Basal (range: 0 and 0.1 to 15 units/ hr)
- Correction Factor (range: 1 unit:0.1 mmol/L to 1 unit:33.3 mmol/L)
- Carb Ratio (range: 1 unit:1 gram to 1 unit:300 grams)

Below a carb ratio of 1:10, increments can be entered in 0.1 g. For example a carb ratio of 1:8.2 can be programmed.

 Target BG (range: 3.9 mmol/L to 13.9 mmol/L)

In addition, you can set any or all of the following Bolus Settings:

 Carbs (on indicates entering grams of Carb; off indicates entering units of insulin)

- Insulin Duration (how long a bolus lowers your BG)
- Max Bolus (the maximum amount for a single bolus)

The default settings and ranges for Bolus Settings are as follows:

- Carbs (default: off if no Carb Ratio is defined)
- Insulin Duration (default: 5 hrs; range: 2 to 8 hrs)
- Max Bolus (default: 10 units; range: 1 to 25 units)

Insulin Duration and Insulin on Board (IOB)

Your t:slim X2 Pump remembers how much insulin you have taken from previous boluses. It does this by relying on the Insulin Duration. The Insulin Duration reflects the amount of time that insulin is actively lowering your BG.


Consult your healthcare provider to accurately set your Insulin Duration.

While the Insulin Duration setting reflects how long insulin from previous boluses lowers your BG, the IOB feature reflects how much insulin is remaining in your body from previous boluses. IOB is always displayed on the Home Screen and is used in bolus delivery calculations when applicable. When a BG is entered during bolus programming, your t:slim X2 Pump will consider any active IOB and calculate an adjusted bolus if necessary.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- 3. Tap Personal Profiles.
- 4. Tap **NEW** to create a new profile.
- Using the onscreen keypad, enter a profile name (up to 16 characters) and tap

Tap once for first letter displayed, 2 quick taps for middle letter; and 3 quick taps for the third letter.

6. Tap Press to Setup to begin setting insulin delivery settings.

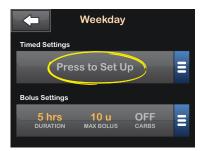

6.3 Programming a New Personal Profile

Once the Personal Profile has been created, the settings must be programmed. The first time segment will start at 00:00 (midnight).

- You must program a basal rate in order to have a Personal Profile that you can activate.
- Be sure to tap after entering or changing a value.

A PRECAUTION

ALWAYS confirm that the decimal point placement is correct when entering your Personal Profile information. Incorrect decimal point placement can prevent you from getting the proper insulin amount that your healthcare provider has prescribed for you.

Timed Settings



- 1. Once the new profile has been created, Tap Basal.
- 2. Using the onscreen keypad, enter your basal rate and tap ...
- 3. Tap Correction Factor.
- 4. Using the onscreen keypad, enter your correction factor (the mmol/L

- that 1 unit of insulin will lower BG) and tap .
- 5. Tap Carb Ratio.
- 6. Using the onscreen keypad, enter your insulin-to-carbohydrate ratio (the grams of carb to be covered by 1 unit of insulin) and tap
- 7. Tap Target BG.
- 8. Using the onscreen keypad, enter your target BG and tap .
- 9. Review entered values and tap
- 10. Confirm Settings.
 - Tap if entered data is correct.
 - Tap to make changes.
- 11. Tap to set the Bolus Settings.

Tap to create additional time segments.



Adding More Time Segments

When adding more time segments, any settings that you entered in the previous time segment are copied and appear in the new segment. This allows you to simply adjust only the specific settings you want, rather than have to enter them all over again.

- 1. On the Add Segment screen, tap Start Time.
- Using the onscreen keypad, enter the time (hour and minutes) that you want the segment to begin, and tap

- 3. On the Add Segment screen, tap Time of Day to select AM or PM.
- Once a time segment is set beyond 12:00 PM, the default will change to PM.
- 4. Tap 🕶
- Repeat steps 1 to 6 from the Creating a New Profile section above for each segment you want to create (up to 16).

To find time segments in the list that are not displayed on the first screen, tap the **Down Arrow**.

Bolus Settings

1. Tap the Bolus Settings Panel.



2. Tap Insulin Duration.



3. Using the onscreen keypad, enter the desired time for the duration of insulin action (2–8 hrs) and tap

Chapter 6 – Personal Profiles

- 4. Tap Max Bolus.
- 5. Using the onscreen keypad, enter the desired amount for maximum bolus (1–25 units) and tap ...

■ NOTE: 25 Unit Max Bolus

If you set the max bolus to 25 units and a bolus larger than 25 units is calculated using your Carb Ratio or Correction Factor, after the bolus is delivered a reminder screen will appear. The option of delivering the remaining amount of the bolus up to an additional 25 units will be given (see Section 15.12 Max Bolus Alert 1).

- Tap the Carbohydrates slider (it should move to the right) to turn on and use the carb ratio when calculating boluses.
- 7. Review entered values and tap

 ...
- 8. Confirm Settings.
 - Tap if entered data is correct.
 - Tap x to make changes.

Tap Tandem Logo to return to the Home Screen.

Adding More Personal Profiles

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- 3. Tap Personal Profiles.
- 4. Tap + .
- Name the new profile and repeat steps for Timed Settings and Bolus Settings.

■ NOTE: Carbohydrates Options

If the first profile you created is programmed using a carb ratio, any new profile will also have the Carbohydrates option turned **ON**, but a ratio will still need to be defined.

6.4 Editing or Reviewing an Existing Profile

1. From the Home Screen, tap OPTIONS.

- 2. Tap My Pump.
- 3. Tap Personal Profiles.
- 4. Tap the name of the Personal Profile to edit or review.
- 5. Tap Edit.

■ NOTE: Review Settings

To review settings but bypass editing the settings, skip the remaining steps in this section. You can tap to navigate to the Personal Profiles list or tap the Tandem Logo to return to the Home Screen.

- 6. Tap Timed Settings Panel.
- 7. Tap the desired time segment to edit.
- 8. Tap Basal, Correction Factor,
 Carb Ratio or Target BG to make
 changes as needed and use the
 onscreen keypad to enter changes.
 Tap ...
- View recent changes and tapConfirm Settings.

- Tap if entered data is correct.
- Tap x to make changes.
- Edit other time segments within the Timed Settings by tapping on them and using the same steps described above.
- 11. Tap after editing all of the time segments.
- 12. Tap the Bolus Settings Panel to change DURATION, MAX BOLUS or CARBS as needed. Use the onscreen keypad to enter desired changes. Tap ...
- 13. Confirm Settings.
 - Tap if entered data is correct.
 - Tap x and make changes.
- 14. Tap **Tandem Logo** to return to the Home Screen.

■ NOTE: Adding a Time Segment

To Add a time segment, tap + and enter the desired start time.

■ NOTE: Deleting a Time Segment

To Delete a time segment, tap on the X to the left of the time segment and tap \checkmark to confirm.

6.5 Duplicating an Existing Profile

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- 3. Tap Personal Profiles.
- 4. Tap the name of the Personal Profile to duplicate.
- 5. Tap Duplicate.
- 6. Confirm profile to duplicate by tapping .
- 7. Using the onscreen keypad, enter the name (up to 16 characters) for the new profile and tap ...

- ✓ Profile Duplicated screen is displayed.
- ✓ A new Personal Profile will be created with the same settings as the profile copied.
- Tap the Timed Settings or Bolus Settings Panel to make changes to the new profile.

6.6 Activating an Existing Profile

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- 3. Tap Personal Profiles.
- 4. Tap the name of the Personal Profile to be activated.
 - The Activate and Delete options are disabled for the active profile because the profile is already activated. You cannot delete a profile until you have activated another profile.

Chapter 6 – Personal Profiles

- If you have only 1 profile defined, you do not need to activate it (that profile is automatically activated).
- 5. Tap Activate.
- A screen to confirm the activation request is displayed.
- 6. Tap ____.
- Profile Activated screen is displayed.

6.7 Renaming an Existing Profile

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- 3. Tap Personal Profiles.
- 4. Tap the name of the Personal Profile to be renamed.
- 5. Tap Down Arrow, and then Rename.

- Using the onscreen keypad, rename the profile name (up to 16 characters) and tap
- 7. Tap **Tandem Logo** to return to the Home Screen.

6.8 Deleting an Existing Profile

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- 3. Tap Personal Profiles.
- 4. Tap the name of the Personal Profile to be deleted.

■ NOTE: Active Personal Profile

The active Personal Profile cannot be deleted.

- 5. Tap Delete.
- 6. Tap 🕶
- Profile Deleted screen is displayed.

7. Tap the **Tandem Logo** to return to the Home Screen.

Chapter 7

Bolus

7.1 Bolus Overview

A WARNING

DO NOT deliver a bolus until you have reviewed the calculated bolus amount on the pump display. If you dose an insulin amount that is too high or too low, this could cause very high or very low blood glucose. You can always adjust the insulin units up or down before you decide to deliver your bolus.

A bolus is a quick dose of insulin that is usually delivered to cover food eaten or to correct a high BG.

The minimum bolus size is 0.05 units. The maximum bolus size is 25 units. If you attempt to deliver a bolus that is larger than the amount of insulin in the cartridge, a message screen appears indicating that there is not enough insulin to deliver the bolus.

Your t:slim X2 Pump offers you the ability to deliver different boluses to cover carbohydrate intake (food bolus) and bring your BG back to target (correction bolus). Food and correction boluses can also be programmed together.

If Carbohydrates is turned on in your active personal profile, you will enter grams of carbohydrate and the bolus will be calculated using your Carb Ratio.

If Carbohydrates is turned off in your active personal profile, you will enter units of insulin to request the bolus.

A PRECAUTION

CHECK your System's personal settings regularly to ensure they are correct. Incorrect settings can result in over delivery or under delivery of insulin. Consult your healthcare provider as needed.

7.2 BG Entry and Correction Bolus Calculation

Once the pump knows your BG value, either from the CGM or from manual entry, it will determine whether to recommend that a correction bolus be added to any other bolus requested on the Bolus Screen.

When your blood glucose is:

- Above Target BG: the insulin for the food bolus and the correction bolus will be added together. If IOB is present, it will only be used in the calculation of the correction portion of the bolus.
- Between 3.9 mmol/L and Target BG: You will be given an option to reduce the food bolus to correct for the low blood glucose. In addition, if IOB is present, it will also be used to reduce the bolus calculation.
- Below 3.9 mmol/L: The food bolus will be reduced to automatically correct for the low blood glucose. In addition, if IOB is present, it will also be used to reduce the bolus calculation.

Always treat hypoglycemia (low blood glucose) with fast-acting carbohydrates according to the instructions of you healthcare provider and then re-test your blood glucose to ensure that the treatment was successful.

BG Value Auto-Population with CGM

If you have a CGM session active, and if there is both a CGM value and a CGM trend arrow available on the CGM Home Screen, the correction bolus confirmation screen is displayed if appropriate as soon as you tap the Bolus button from the CGM Home Screen.

BG Value Manual Entry

If you are not using a CGM, or if your CGM value or trend arrow are not available on the Home Screen, the correction bolus confirmation screen is displayed if appropriate after you manually enter your BG value on the Bolus Screen.

- 1. From the Home Screen tap Bolus.
- 2. Tap Add BG.
- Using the onscreen keypad, enter your BG value and tap . The BG value is saved in pump History whether or not a bolus is delivered.

BG Above Target

If your BG is above your Target BG, the pump presents you with the option for the pump to calculate and add a correction bolus to any other bolus you request.



- To accept the correction bolus press . A correction bolus is calculated and will be added to any food bolus you request on the Bolus Screen.
- To decline the correction bolus, press X. No correction bolus will be added to any food bolus you request on the Bolus Screen.

BG Below Target

If your BG is below your Target BG, the pump presents you with the option for the pump to calculate and subtract a correction bolus from any other bolus you request.



- To accept the correction bolus press . A correction bolus is calculated and will be subtracted from any food bolus you request on the Bolus Screen.
- To decline the correction bolus, press . No correction bolus will be subtracted from any food bolus you request on the Bolus Screen.

BG Within Target

If your BG is within your target range, no correction bolus screen is displayed.

You can override the calculated bolus by tapping on the calculated units value and entering the units of insulin you want delivered. The bolus override is always active.



7.3 Food Bolus Using Units

If bolusing using a carb ratio, skip to Section 7.4 Food Bolus Using Grams.

- 1. From the Home Screen, tap BOLUS.
- 2. Tap 0 units.

3. Using the onscreen keypad enter units of insulin to be delivered, then tap ...

APRECAUTION

ALWAYS confirm that the decimal point placement is correct when entering bolus information. Incorrect decimal point placement can prevent you from getting the proper amount of insulin that your health-care provider has prescribed for you.

- 4. Tap to confirm the units of insulin to be delivered.
- 5. Confirm Request.
 - Tap if entered data is correct.
 - Tap x to go back to make changes or view calculations.
- 6. Tap ____.
- The bolus initiated screen is temporarily displayed.

7.4 Food Bolus Using Grams

- 1. From the Home Screen, tap BOLUS.
- 2. Tap 0 grams.
- 3. Using the onscreen keypad enter grams of carb and tap ...
 - To add multiple carb values enter first value, then tap + , enter second value, tap + . Continue until done.
 - To clear the value entered and start over, tap the button.
- Check that the grams of carb are entered in the correct location on the screen.
- 5. Tap to confirm the units of insulin to be delivered.

You can always tap View Calculation to display the Delivery Calculation screen.

Confirm Request.

- Tap if entered data is correct.
- Tap X to go back to make changes or view calculations.
- 7. Tap DELIVER.
- The bolus initiated screen is temporarily displayed.

7.5 Extended Bolus

The Extended Bolus feature allows you to deliver part of the bolus now and part of the bolus slowly over a period of up to 8 hours. This can be helpful for high fat meals such as pizza or if you have gastroparesis (delayed stomach emptying).

When using extended bolus, any correction bolus amount will always be given in the DELIVER NOW portion. Talk with your healthcare provider to determine if this feature is appropriate for you as well as for recommendations on the split between now and later and the duration for the later portion.

Only 1 extended bolus can be active at any given time. However, if the DELIVER LATER portion of an extended bolus is active, you can request another standard bolus.

- 1. From the Home Screen, tap BOLUS.
- 2. Tap 0 grams (or 0 units).
- 3. Using the onscreen keypad enter grams of carb (or units of insulin). Tap ...
- 4. If desired, tap Add BG and using the onscreen keypad enter BG value. Tap
- 5. Tap to confirm the units of insulin to be delivered.

You can always tap View Calculation to display the Delivery Calculation screen.

- 6. Confirm Request.
 - Tap if entered data is correct.

- Tap to go back to make changes or view calculations.
- Tap 50% under DELIVER NOW to adjust the percentage of the food bolus that is to be delivered immediately.

The percentage value for DELIVER LATER is automatically calculated by the pump. The default is 50% NOW and 50% LATER. The default for DURATION is 2 hours.

 Use the onscreen keypad to enter the percentage of the bolus to DELIVER NOW and tap

For the DELIVER NOW portion, the minimum amount is .05 units. If the DELIVER NOW portion is less than .05 units, you will be notified and the DELIVER NOW portion will be set to .05 units.

10. Tap 2 hrs under DURATION.

Chapter 7 – Bolus

The DELIVER LATER portion of the Extended Bolus also has minimum and maximum rates. If you program a DELIVER LATER rate outside of these limits, you will be notified and the duration of the DELIVER LATER portion will be adjusted.

- 11. Use the onscreen keypad to adjust the length of time the bolus is to be delivered, then tap ...
- 12. Tap 🕶.

You can always tap View Units to display the breakdown of units to be delivered NOW versus LATER.

- 13. Confirm Request.
 - Tap if entered data is correct.
 - Tap to go back to make changes or view calculations.
- 14. Tap 🕶
- 15. The bolus initiated screen is temporarily displayed.

7.6 Canceling or Stopping a Bolus

Canceling a Bolus If delivery HAS NOT STARTED:

- 1. Tap 1–2–3 to access the Home Screen.
- 2. Tap X (stop icon) to cancel the bolus.



- The bolus button will remain inactive while the bolus is being canceled.
- Once canceled, the bolus button will become active again on the Home Screen.

Stopping a Bolus if delivery of the BOLUS HAS STARTED:

- 1. Tap 1–2–3 to access the Home Screen.
- 2. Tap X (stop icon) to stop delivery.
- Tap
- √ The BOLUS STOPPED screen is displayed and the units delivered are calculated.
- 4. Units requested and delivered are shown. Tap ok .

Chapter 8

Stop/Resume Insulin

8.1 Stopping Insulin Delivery

You can stop all insulin delivery at any time. When you stop all insulin delivery, any active bolus and any active temp rate are immediately stopped. No insulin delivery can take place while your pump is stopped.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap STOP INSULIN.
- 3. Tap STOP.
- ✓ The All Deliveries Stopped screen appears before returning to the Home Screen showing the status ALL DELIVERIES STOPPED. A red exclamation mark icon also appears to the right of the time and date.

8.2 Resuming Insulin Delivery

If pump screen is not on, press Screen On/Quick Bolus Button once to turn on your t:slim X2 Pump screen.

- 1. Tap 1–2–3 to access the Home Screen.
- 2. Tap RESUME.
- √ The RESUMING INSULIN screen is temporarily displayed.
- OR -
- 1. From the Home Screen, tap OPTIONS.
- 2. Tap RESUME INSULIN.
- 3. Tap RESUME.
- ✓ The RESUMING INSULIN screen is temporarily displayed.

THIS PAGE IS INTENTIONALLY LEFT BLANK

Section 3

Additional Pump Operations

Chapter 9

Temporary Basal Rate

9.1 Setting a Temp Rate

A Temp Rate is used to increase or decrease (by percentage) the current basal rate for a period of time. This feature can be helpful for situations such as exercise or illness.

When you enter the Temp Rate screen, the default values are 100% (current basal rate) and a Duration of 0:15 min. The Temp Rate can be set from a minimum of 0% of current basal rate to a maximum of 250% of current basal rate in increments of 1%.

Duration can be set from a minimum of 15 minutes to a maximum of 72 hours in increments of 1 minute.

If you program a Temp Rate greater than 0% but less than the minimum allowable basal rate of 0.1 units/hr, you will be notified that the selected rate is too low and that it will be set to the minimum allowable rate for delivery.

If you program a Temp Rate more than the maximum allowable basal rate of 15 units/hr, you will be notified that the selected rate is too high and that it will be set to the maximum allowable rate for delivery.

NOTE: Temp Rate and Insulin Suspension

If insulin is manually suspended for any reason during a Temp Rate, the time remaining for the temp rate will continue to tick down during the suspension period.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap Temp Rate.
- Using the onscreen keypad enter desired percentage. The current rate is 100%. An increase is greater than 100% and decrease is less than 100%.
- 4. Tap ...
- Tap Duration. Using the onscreen keypad enter desired length of time for Temp Rate. Tap

You can always tap View Units to see the actual units to be delivered.

- 6. Verify settings and tap .
- √ The "TEMP RATE STARTED" screen is temporarily displayed.
- The Screen Lock screen will be displayed with the icon indicating a Temp Rate is active.
 - An orange "T" means a Temp Rate is active.
 - A red "T" means a Temp Rate of 0 is active.

9.2 Stopping a Temp Rate

To stop an active temp rate:

- 1. From the Home Screen, tap OPTIONS.
- On the Options screen, tap X (stop icon) on the right side of Temp Rate.
- 3. On the confirmation screen, tap

✓ The TEMP RATE STOPPED screen appears before returning to the Options screen. THIS PAGE IS INTENTIONALLY LEFT BLANK

Chapter 10

Quick Bolus

10.1 Setting Up Quick Bolus

Setting up the Quick Bolus function enables you to deliver a bolus by simply pressing a button. It is a way to deliver a bolus by following beep/vibration commands without navigating through or viewing the pump screen.

The default for the Quick Bolus function is off. Quick Bolus can be set to either units of insulin or grams of carbohydrate. The increment options are 0.5, 1.0, 2.0, and 5.0 units; or 2, 5, 10 and 15 grams of carb.

The quick bolus delivery setting (grams of carbohydrate or units of insulin) is independent of the active Personal Profile bolus setting.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- Tap Pump Settings.
- Tap Quick Bolus.

- 5. Tap Increment Type.
- 6. Tap units of insulin or grams of carbohydrate to select.
- 7. Tap Increment Amount.
- 8. Select the preferred increment amount.

■ NOTE: Incrementing Amount

The increment amount is added with each press of the Quick Bolus Button when delivering a quick bolus.

- 9. Review entered values and tap
- 10. Confirm Settings.
 - Tap if entered data is correct.
 - Tap x to go back to make changes.
- 11. Tap **Tandem Logo** to return to the Home Screen.

10.2 Delivering Quick Bolus

If the Quick Bolus function is turned On, you can deliver a bolus without having to look at the t:slim X2 Pump's screen. Simply use the Quick Bolus Button to deliver your bolus. Quick boluses are delivered as standard boluses (there is no BG entry or extended bolus).

A PRECAUTION

ALWAYS look at the screen to confirm correct programming of the bolus amount when you first use the Quick Bolus feature. Looking at the screen will ensure that you are correctly using the beep/vibration commands to program the intended bolus amount.

- Press and hold the Quick Bolus Button. The Quick Bolus screen will appear. Listen for 2 beeps (if Sound Volume is set to beep) or feel for vibrations (if Sound Volume is set to vibrate).
- Press the Quick Bolus Button for each increment until desired amount is reached. The pump will beep/ vibrate for each button press.

- The pump will beep/vibrate once for each increment pressed to confirm desired amount.
- After the pump beeps/vibrates, press and hold the Quick Bolus Button to deliver the bolus.

■ NOTE: Safety Features

If you want to cancel the bolus and return to the Home Screen, tap x on the QUICK BOLUS screen.

If more than 10 seconds have passed with no input, the bolus is canceled and never delivered.

You cannot exceed the Max Bolus setting defined in your active Personal Profile when using the Quick Bolus feature. Once you reach the Max Bolus amount, a different tone will sound to notify you (if Quick Bolus is set to vibrate, the pump will stop vibrating in response to additional button presses to notify you). Look at the screen to confirm the bolus amount.

You cannot exceed 20 button presses when using the Quick Bolus feature. Once you reach 20 button presses, a different tone will sound to notify you (if Quick Bolus is set

to vibrate, the pump will stop vibrating in response to additional button presses to notify you). Look at the screen to confirm the bolus amount.

If you hear a different tone at any point during programming or the pump stops vibrating in response to button presses, look at the screen to confirm the bolus amount. If the Quick Bolus screen does not display the correct bolus amount, use the touch screen to enter bolus information.

✓ The bolus initiated screen is temporarily displayed.

THIS PAGE IS INTENTIONALLY LEFT BLANK

Chapter 11

t:slim X2 Pump Settings

Chapter 11 – t:slim X2 Pump Settings

This chapter contains explanations and instructions for Pump Settings that are not covered in other chapters.

11.1 Display Settings

You can set the Screen Timeout to the length of time you want the screen to stay on before it automatically turns off. The default for the Screen Timeout is 30 seconds. The options are 15, 30, 60, and 120 seconds.

You can always turn the screen off before it automatically times out by pressing the Screen On/Quick Bolus button.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- 3. Tap Pump Settings.
- 4. Tap Display Settings.
- Tap Screen Timeout.
- 6. Select preferred time and tap

7. Tap **Tandem Logo** to return to the Home Screen.

11.2 Mobile Device Connection

You can connect one mobile device to the pump to display pump information on your mobile device.

The mobile app may not yet be available in your area.

Pair a Mobile Device

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- Tap Pump Settings.
- 4. Tap Bluetooth Settings.
- Tap the Mobile Connection slider (it should move to the right) to turn on the mobile connection. The Pair Device button is displayed.



- 6. Tap Pair Device.
- 7. A mobile app reminder screen is displayed.
- 8. Tap to generate your device pairing code.
- Follow the instructions to enter the pairing code into the app on your mobile device.
- 10. Tap ox to close the PAIRING CODE screen. If the pump has successfully paired with a mobile device, the DEVICE PAIRED screen is displayed.

11. Tap ok to close the DEVICE PAIRED screen.

11.3 t:slim X2 Sound Volume

Sound Volume is preset to high. A change to the Sound Volume can be made in Pump Settings.

Sound Volume can be personalized for the Button taps, Quick Bolus, Bolus, Reminders, Alerts, Alarms, and Fill Tubing. Options for Sound Volume include high, medium, low, and vibrate.

A PRECAUTION

DO NOT use the vibrate feature for alerts and alarms during sleep unless otherwise directed by your healthcare provider. Having the volume for alerts and alarms set to high will help ensure that you don't miss an alert or alarm.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- Tap Pump Settings.
- 4. Tap Sound Volume.

- 5. Tap desired option. Use **Up/Down Arrow** to view additional options.
- 6. Select preferred volume.
- 7. Continue to make changes for all Sound Volume options by repeating steps 5 and 6.
- 8. Tap when all changes are complete.

Tap **Tandem Logo** to return to the Home Screen.

11.4 Turn Feature Lock ON or OFF

Feature Lock is preset to off. With the Feature Lock turned on, you cannot deliver a bolus, change any pump settings or access any Personal Profiles.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- 3. Tap Pump Settings.

- 4. Tap **Feature Lock** to turn On or Off and tap ...
- A screen to verify that you want to activate the Feature Lock will be displayed.
- 5. Tap to confirm.
- 6. Tap 🛶.
- 7. Tap **Tandem Logo** to return to the Home Screen.

■ NOTE: Sound Volume

When Feature Lock is turned on, Sound Volume will be set to high until the Feature Lock is turned off.

THIS PAGE IS INTENTIONALLY LEFT BLANK

Chapter 12

t:slim X2 Pump Info and History

12.1 t:slim X2 Pump Info

Your t:slim X2 Pump allows access to information about your pump. In the Pump Info screen you have access to items such as your pump Serial Number and software/hardware versions.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- 3. Tap Pump Info.
- 4. Scroll through the Pump Info using the Up/Down Arrows.
- 5. Tap **Tandem Logo** to return to the Home Screen.

12.2 t:slim X2 Pump History

Pump History displays a historical log of pump events. At least 90 days of data can be viewed in History. When the maximum number of events is reached, the oldest events are removed from the history log and replaced with the most recent events. The following can be viewed in History:

Delivery Summary, Total Daily Dose, Bolus, Basal, Load, BG, Alerts & Reminders, and Complete.

Delivery Summary breaks down total insulin delivery by basal and bolus types into units and percentages. It can be viewed by the selected time period of: Today, 7 Day, 14 Day and 30 Day Average.

Total Daily Dose breaks down basal and bolus delivery into units and percentages for each individual day. You can scroll through each individual day to see your total insulin delivery.

The Bolus, Basal, Load, BG, Alerts & Reminders, and Complete are categorized by date. The event details in each report are listed by time.

The letter "D" (D: Alert) before an Alert or Alarm indicates the time it was declared. The letter "C" (C: Alert) indicates the time it was cleared.

Bolus History shows the bolus request, the bolus start time, and the bolus completion time.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap the Down Arrow.
- 3. Tap Pump History.
- 4. Tap desired option.
- 5. Tap **Tandem Logo** to return to the Home Screen.

Chapter 13

t:slim X2 Pump Reminders

Chapter 13 – t:slim X2 Pump Reminders

Your t:slim X2 Pump lets you know important information about the System with Reminders, Alerts, and Alarms. Reminders are displayed to notify you of an option that you have set (for example, a reminder to check you BG after a bolus). Alerts display automatically to notify you about safety conditions that you need to know (for example, an alert that your insulin level is low). Alarms display automatically to let you know of an actual or potential stopping of insulin delivery (for example, an alarm that the insulin cartridge is empty). Pay special attention to Alarms.

If multiple Reminders, Alerts, and Alarms happen at the same time, Alarms will be displayed first, Alerts will be displayed second, and Reminders will be displayed third. Each must be confirmed separately until all have been confirmed.

Information in this section will help you learn how to respond to Reminders.

Reminders notify you with a single sequence of 3 notes or a single vibration depending on the volume/vibrate setting in Sound Volume. They repeat

every 10 minutes until acknowledged. Reminders do not escalate.

13.1 Low BG Reminder

The Low BG Reminder prompts you to re-test your BG after a low BG value is entered. When turning this reminder on, you need to set a low BG value that triggers the reminder, as well as how much time should pass before the reminder occurs.

The default for this reminder is preset to off. If on, Remind Me Below 3.9 mmol/L, and Remind Me After 15 min, but you can set these values from 3.9 mmol/L to 6.7 mmol/L and 10 to 20 min.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- Tap Alert & Reminders.
- 4. Tap Pump Reminders.
- 5. Tap Low BG.

- 6. Low BG is set to on; to turn off, tap LOW BG.
 - a. Tap Remind Me Below and using the onscreen keypad, enter a Low BG value (from 3.9 mmol/L to 6.7 mmol/L) that you want to trigger the reminder, then tap
 - b. Tap Remind Me After and using the onscreen keypad, enter the time (from 10 to 20 min), then tap
 - c. Tap when all changes are complete.
 - d. Tap **Tandem Logo** to return to the Home Screen.

To Respond to the Low BG Reminder

To clear the reminder, tap or and then check BG using your blood glucose meter.

13.2 High BG Reminder

The High BG Reminder prompts you to re-test your BG after a high BG value is entered. When you turn this reminder on, you need to set a high BG value that triggers the reminder, as well as how much time should pass before the reminder occurs.

The default for this reminder is preset to off. If on, Remind Me Above 11.1 mmol/L, and Remind Me After 120 min, but you can set these values from 8.3 to 16.7 mmol/L and 1 to 3 hrs.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- 3. Tap Alerts & Reminders.
- 4. Tap Pump Reminders.
- 5. Tap High BG.
- 6. High BG is set to on; to turn off, tap High BG.

- a. Tap Remind Me Above and using the onscreen keypad, enter a High BG value (from 8.3 to 16.7 mmol/L) that you want to trigger the reminder, then tap
- b. Tap Remind Me After and using the onscreen keypad, enter the time (from 1 to 3 hours), then tap
- c. Tap when all changes are complete.
- 7. Tap **Tandem Logo** to return to the Home Screen.

To Respond to the High BG Reminder

To clear the reminder tap or and then check BG using your blood glucose meter.

13.3 After Bolus BG Reminder

The After Bolus BG Reminder prompts you to test your BG at a selected time after bolus delivery. When turning this

reminder on, you need to set how much time should pass before the reminder occurs. The default is 1 hour and 30 minutes. It can be set from 1 to 3 hours.

- 1. From the Home Screen, tap OPTIONS
- 2. Tap My Pump.
- 3. Tap Alerts & Reminders.
- 4. Tap Pump Reminders.
- 5. Tap After Bolus BG.
- 6. After Bolus BG is set to on; to turn off, tap After Bolus BG.
- 7. Tap Remind Me After and using the onscreen keypad, enter the time (from 1 to 3 hours) that you want to trigger the reminder, then tap ...
- 8. Tap when all changes are complete.
- 9. Tap **Tandem Logo** to return to the Home Screen.

Chapter 13 – t:slim X2 Pump Reminders

To Respond to the After Bolus BG Reminder

To clear the reminder tap or and then check BG using your blood glucose meter.

13.4 Missed Meal Bolus Reminder

The Missed Meal Bolus Reminder lets you know if a bolus was not delivered during a specified time period. Four separate reminders are available. When programming this reminder you need to select the Days, the Start Time, and End Time for each reminder.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- 3. Tap Alerts & Reminders.
- Tap Pump Reminders.
- Tap Missed Meal Bolus.
- 6. On the Missed Meal Bolus screen, tap which reminder you want to set

(Reminder 1 to 4) and do the following:

- a. Tap Reminder 1 (or 2, 3, 4).
- b. Reminder 1 is set to on; to turnoff, tap Reminder 1.
- c. Tap Selected Days and tap the day(s) you want the reminder to be on, then tap
- d. Tap Start Time, tap Time and using the onscreen keypad enter the start time, then tap
- e. Tap **Time of Day** to select AM or PM, then tap
- f. Tap End Time, tap Time and using the onscreen keypad enter the end time, then tap
- g. Tap Time of Day to select AM or PM, then tap .
- h. Tap when all changes are complete.

7. Tap the **Tandem Logo** to return to the Home Screen.

To Respond to the Missed Meal Bolus Reminder

To clear the reminder tap and deliver a bolus if necessary.

13.5 Site Reminder

The Site Reminder prompts you to change your infusion set. The default for this reminder is preset to off. If on, the reminder can be set for 1 to 3 days and at a time of day selected by you.

For detailed information on the Site Reminder feature, see Section 5.6 Setting Site Reminder.

To Respond to the Site Reminder

To clear the reminder tap and change your infusion set.

Chapter 14

User Settable Alerts & Reminders

14.1 Low Insulin Alert

Your t:slim X2 Pump keeps track of how much insulin remains in the cartridge and alerts you when it is low. The default for this alert is preset to 20 units. You can set this alert setting anywhere between 10 and 40 units. When the insulin amount reaches the set value, the Low Insulin Alert beeps/ vibrates and appears on the screen. After the alert is cleared, the low insulin indicator (a single red bar on the insulin level display on the Home Screen appears).


- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- 3. Tap Alerts & Reminders.
- Tap Pump Alerts.
- 5. Tap Low Insulin.
- Using the onscreen keypad, enter the number of units (from 10 to 40 units) that you want the Low Insulin

Alert value to be set to, and tap

7. Tap when all changes are complete.

To Respond to the Low Insulin Alert

To clear the alert, tap ok.

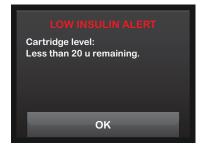

14.2 Auto-Off Alarm

Your t:slim X2 Pump can stop insulin delivery and alert you (or whoever is with you) if there has been no interaction with the pump within a specified period of time. The default for this alarm is preset to 12 hours. You can set it anywhere between 5 and 24 hours, or

off. This alarm notifies you that there has been no interaction with the pump in the specified number of hours and the pump will shut down after 30 seconds.

When the number of hours since you have pressed the Screen On/Quick Bolus Button and tapped any interactive screen option or delivered a Quick Bolus passes the set value, the Auto-Off Alarm beeps and appears on the screen, and insulin delivery stops.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My Pump.
- 3. Tap Alerts & Reminders.
- 4. Tap Pump Alerts.
- 5. Tap Auto-Off. A confirmation screen will appear.
 - Tap to continue.
 - Tap to go back.

- 6. Verify Auto-Off is set to on, then tap Time.
- 7. Using the onscreen keypad, enter the number of hours (from 5 to 24 hrs) that you want the Auto-Off Alarm to be triggered, and tap
- 8. Tap , then tap when all changes are complete.
- 9. Tap **Tandem Logo** to return to the Home Screen.

To Respond to Auto-Off Warning Tap DO NOT SHUT DOWN.




✓ The warning clears and the pump returns to normal operation.

If you do not clear the warning within the 30-second countdown period, the AUTO-OFF ALARM occurs, accompanied by an audible alarm. This alarm notifies you that your pump has stopped delivering insulin.

Auto-Off Alarm Screen

Тар ок .

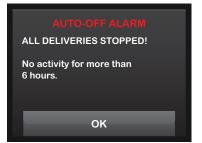

The Home Screen appears, indicating a status of "All Deliveries Stopped."

You must resume delivery to continue therapy. See Section 8.2 Resuming Insulin Delivery.

Section 4

Pump Safety Alerts & Reminders

Chapter 15

t:slim X2 Pump Alerts

Chapter 15 – t:slim X2 Pump Alerts

Your t:slim X2 Pump lets you know important information about the System with Reminders, Alerts, and Alarms. Reminders are displayed to notify you of an option that you have set (for example, a reminder to check you BG after a bolus). Alerts display automatically to notify you about safety conditions that you need to know (for example, an alert that your insulin level is low). Alarms display automatically to let you know of an actual or potential stopping of insulin delivery (for example, an alarm that the insulin cartridge is empty). Pay special attention to Alarms.

If multiple Reminders, Alerts, and Alarms happen at the same time, Alarms will be displayed first, Alerts will be displayed second, and Reminders will be displayed third. Each must be confirmed separately until all have been confirmed.

Information in this section will help you learn how to respond to Alerts.

Alerts notify you with 2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume. They repeat regularly

until acknowledged. Alerts do not escalate.

■ NOTE: CGM Alerts

There is an additional list of Alerts and Errors related to CGM use in Chapter 26 CGM Alerts and Errors.

15.1 Low Insulin Alert

What will I see on the screen?



What does it mean?

5 units or less of insulin remain in the cartridge.

How will the System notify me?

2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes, every 5 minutes until acknowledged.

How should I respond?

Tap over . Change your cartridge as soon as possible to avoid the EMPTY CARTRIDGE ALARM and running out of insulin.

15.2 Low Power Alert 1

What will I see on the screen?



What does it mean?

Less than 25% of battery power remains.

How will the System notify me?

2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes, every 5 minutes until acknowledged.

How should I respond?

Tap ox . Charge your pump as soon as possible to avoid the second LOW POWER ALERT.

■ NOTE: Low Battery Display

Once the LOW POWER ALERT occurs, the low-power indicator (a single red bar on the battery level display on the Home Screen) appears.

15.3 Low Power Alert 2

What will I see on the screen?

LOW POWER ALERT Recharge pump or all deliveries will stop. OK

What does it mean?

Less than 5% of battery power remains. Insulin delivery will continue for 30 minutes and then the pump will power off and insulin delivery will stop.

How will the System notify me?

2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes, every 5 minutes until acknowledged.

How should I respond?


Tap ox. Charge your pump immediately to avoid the LOW POWER ALARM and system power off.

■ NOTE: Low Battery Display

Once the LOW POWER ALERT occurs, the low-power indicator (a single red bar on the battery level display on the Home Screen) appears.

15.4 Incomplete Bolus Alert

What will I see on the screen?

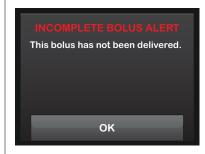

What does it mean?

You started a bolus request but did not complete the request within 90 seconds.

How will the System notify me?

2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes, every 5 minutes until acknowledged.

How should I respond?

Tap ox. The Bolus screen will appear. Continue with your bolus request. Tap if you do not want to continue your bolus request.

15.5 Incomplete Temp Rate Alert

What will I see on the screen?

INCOMPLETE TEMP RATE This temp rate has not been started. OK

What does it mean?

You started to set up a temp rate but did not complete the request within 90 seconds.

How will the System notify me?

2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?


Yes, every 5 minutes until acknowledged.

How should I respond?

Tap . The Temp Rate screen will appear. Continue setting up your temp rate. Tap if you do not want to continue setting up your temp rate.

15.6 Incomplete Cartridge Change Alerts

What will I see on the screen?

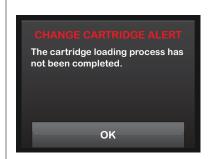

What does it mean?

You selected Change Cartridge from the Load menu but did not complete the process within 3 minutes.

How will the System notify me?

2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes, every 5 minutes until acknowledged.

How should I respond?

Tap ox. Complete the cartridge change process.

15.7 Incomplete Fill Tubing Alert

What will I see on the Screen?

FILL TUBING ALERT The fill tubing process has not been completed. OK

What does it mean?

You selected Fill Tubing from the Load menu but did not complete the process within 3 minutes.

How will the System notify me?

2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes, every 5 minutes until acknowledged.

How should I respond?

Tap OK. Complete the fill tubing process.

15.8 Incomplete Fill Cannula Alert

What will I see on the Screen?



What does it mean?

You selected Fill Cannula from the Load menu but did not complete the process within 3 minutes.

How will the System notify me?

2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes, every 5 minutes until acknowledged.

How should I respond?

Tap OK . Complete the cannula fill process.

15.9 Incomplete Setting Alert

What will I see on the screen?

INCOMPLETE SETTING A setting was being modified, but has not been saved.

What does it mean?

You started to set up a new Personal Profile but did not save or complete the programming within 5 minutes.

How will the System notify me?

2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?


Yes, every 5 minutes until acknowledged.

How should I respond?

Tap OK. Complete programming the Personal Profile.

15.10 Basal Rate Required Alert

What will I see on the screen?

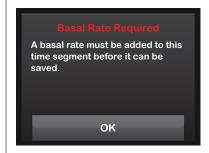

What does it mean?

You did not enter a basal rate in a time segment in Personal Profiles. A basal rate must be entered in each time segment (rate can be 0 u/hr).

How will the System notify me?

Prompt screen only.


Will the System re-notify me?

No. A basal rate must be entered to save the time segment.

How should I respond?

Tap . Enter a basal rate in the time segment.

15.11 Max Hourly Bolus Alert

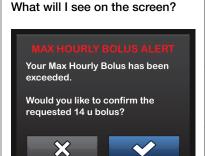

What does it mean?

In the previous 60 minutes, you requested total bolus delivery that is more than 1.5 times your Max Bolus setting.

How will the System notify me?

Prompt screen only.

Will the System re-notify me?

No. You must tap x or to deliver the bolus.

How should I respond?

- Tap to return to the Bolus screen and adjust the bolus delivery amount.
- Tap to confirm the bolus.

15.12 Max Bolus Alert 1

What will I see on the screen?



What does it mean?

You requested a bolus larger than the Max Bolus setting in your active Personal Profile.

How will the System notify me?

Prompt screen only.

Will the System re-notify me?


No. You must tap X or V to deliver the bolus.

How should I respond?

- Tap x to return to the Bolus screen and adjust the bolus delivery amount.
- Tap to confirm the amount of your Max Bolus setting.

15.13 Max Bolus Alert 2

The following applies only if you have Carbs turned on in your active Personal Profile and your Max Bolus amount is set to 25 units.

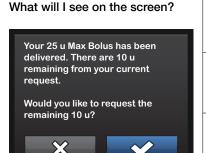

What does it mean?

Your Max Bolus is set to 25 units and you requested a bolus larger than 25 units.

How will the System notify me?

Prompt screen only.

Will the System re-notify me?

No. You must tap x or v to deliver the remaining amount of the bolus request.


How should I respond?

Before responding to this Alert, always consider whether your bolus insulin needs have changed since you requested the original bolus.

- Tap to confirm the remaining amount of the bolus request. A confirmation screen will appear.
- Tap x if you do not want to deliver the remaining amount of the bolus request.

15.14 Max Basal Alert 1

What will I see on the screen?

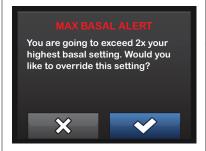

What does it mean?

When entering a basal rate or requesting a temp rate, you requested a basal rate at least two times higher than the Max Basal set for the pump.

How will the System notify me?

Prompt screen only.

Will the System re-notify me?

No.

How should I respond?

- Tap and review the temp rate in the Options menu.
- Tap to dismiss the alert and continue with the request.

15.15 Max Basal Alert 2

What will I see on the screen?

MAX BASAL ALERT

You have exceeded 2x your highest basal setting. Please review your current temp rate in the Options menu.

OK

What does it mean?

The basal rate in the current time segment of the active Personal Profile is higher than the Max Basal set for the pump, which applies to all Personal Profiles. The pump has limited the basal rate to the Max Basal set for the pump.

How will the System notify me?

Prompt screen only.

Will the System re-notify me?


No.

How should I respond?

Tap ox and review the basal rate in the active Personal Profile.

15.16 Min Basal Alert 1

What will I see on the screen?

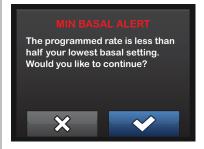

What does it mean?

When entering a basal rate or requesting a temp rate, you requested a basal rate less than half of the lowest basal rate defined in your Personal Profile.

How will the System notify me?

Prompt screen only.

Will the System re-notify me?

No. You must tap X or V to move forward.

How should I respond?

- Tap x to return to the previous screen to adjust the amount.
- Tap to dismiss the alert and continue with the request.

15.17 Min Basal Alert 2

What will I see on the screen?

MIN BASAL ALERT

You have dropped below half your lowest basal setting. Please review your current temp rate in the Options menu.

OK

What does it mean?

An active temp rate dropped below half of your lowest basal setting defined in your Personal Profile.

How will the System notify me?

2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes, every 5 minutes until acknowledged.

How should I respond?

Tap and review your current temp rate in the Options menu.

15.18 Connection Error Alert

What will I see on the screen?

CONNECTION ERROR ALERT Pump cannot connect with the computer. Close this prompt and reconnect the USB cable to try again. OK

What does it mean?

You connected your t:slim X2 Pump to a computer with the USB cable to charge it and a connection could not be made.

How will the System notify me?

2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes, every 5 minutes until acknowledged.

How should I respond?

Tap ox . Disconnect and reconnect the USB cable to try again.

15.19 Pairing Code Timeout



PAIRING CODE TIMEOUT

Close this prompt and press the Pair Device button to generate a new Pairing Code.

CLOSE

What does it mean?

You attempted to connect a mobile device to the pump, but the pairing process took too long (more than 5 minutes) and was unsuccessful.

How will the System notify me?

2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

No.

How should I respond?

Tap ok. Try to pair the mobile device again.

15.20 Power Source Alert

What will I see on the screen?

POWER SOURCE ALERT The pump cannot charge using the current power source. Please try a different power source. OK

What does it mean?

You connected your t:slim X2 Pump to a power source that does not have enough power to charge the pump.

How will the System notify me?

2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes, every 5 minutes until acknowledged.

How should I respond?

Tap ok. Connect the pump to a different power source to charge.

15.21 Data Error Alert

What will I see on the screen?

DATA ERROR ALERT Please verify that your active profile and pump settings are accurate. OK

What does it mean?

Your t:slim X2 Pump encountered a condition that could potentially result in a loss of data.

How will the System notify me?

2 sequences of 3 notes or 2 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes, every 5 minutes until acknowledged.

How should I respond?

Tap ox. Check your Personal Profiles and pump settings to verify that they are accurate. See Section 6.4 Editing or Reviewing an Existing Profile.

THIS PAGE IS INTENTIONALLY LEFT BLANK

Chapter 16

t:slim X2 Pump Alarms

Chapter 16 – t:slim X2 Pump Alarms

A PRECAUTION

CHECK your System regularly for potential alarm conditions that may display. It is important to be aware of conditions that may affect insulin delivery and require your attention so you can respond as soon as possible.

Your t:slim X2 Pump lets you know important information about the System with Reminders, Alerts, and Alarms. Reminders are displayed to notify you of an option that you have set (for example, a reminder to check you BG after a bolus). Alerts display automatically to notify you about safety conditions that you need to know (for example, an alert that your insulin level is low). Alarms display automatically to let you know of an actual or potential stopping of insulin delivery (for example, an alarm that the insulin cartridge is empty). Pay special attention to Alarms.

If multiple Reminders, Alerts, and Alarms happen at the same time, Alarms will be displayed first, Alerts will be displayed second, and Reminders will be displayed third. Each must be confirmed separately until all have been confirmed.

Information in this section will help you learn how to respond to Alarms.

Alarms notify you with 3 sequences of 3 notes or 3 vibrations depending on the volume/vibrate setting selected in Sound Volume. If not acknowledged, alarms escalate to highest volume and vibe. Alarms repeat regularly until the condition that caused the alarm is corrected.

■ NOTE: CGM Alerts and Errors

There is an additional list of Alerts and Errors related to CGM use in Chapter 26 CGM Alerts and Errors.

16.1 Resume Pump Alarm

What will I see on the screen?

RESUME PUMP ALARM

The pump has been stopped for an extended period of time.

Select RESUME INSULIN in the Options menu to continue therapy.

OK

What does it mean?

You tapped STOP INSULIN in the Options menu and insulin delivery has been stopped for more than 15 minutes.

How will the System notify me?

3 sequences of 3 notes or 3 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes.

- If not acknowledged by tapping ox, the system will re-notify you every 3 minutes at highest volume and vibrate.
- If acknowledged by tapping ____, the system will re-notify you in 15 minutes if you still have not resumed insulin.

How should I respond?

To resume insulin, from the Options menu, tap **RESUME INSULIN** and tap **RESUME** to confirm.

16.2 Low Power Alarm 1

What will I see on the screen?



What does it mean?

Your t:slim X2 Pump detected a power level of 1% or less remaining and all deliveries have stopped.

How will the System notify me?

3 sequences of 3 notes or 3 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes, the System will re-notify you every 3 minutes until no power remains and the pump shuts down.

How should I respond?

Tap OK. Charge your pump immediately to resume insulin delivery.

16.3 Low Power Alarm 2

What will I see on the screen?

LOW POWER ALARM

ALL DELIVERIES STOPPED!

Your pump is about to shut down. Please charge your pump immediately.

SILENCE ALARM

What does it mean?

Your t:slim X2 Pump detected a voltage level too low to ensure normal performance and all deliveries have stopped.

How will the System notify me?

A series of rapid beeps at maximum volume for at least 20 seconds before the System shuts down.

Will the System re-notify me?

The pump will power back on once it has been plugged into a charging source and has reached an adequate level of charge.

How should I respond?

Tap SILENCE ALARM. Charge your pump immediately to resume insulin delivery.

16.4 Empty Cartridge Alarm

What will I see on the screen?



What does it mean?

Your t:slim X2 Pump detected that the cartridge is empty and all deliveries have stopped.

How will the System notify me?

3 sequences of 3 notes or 3 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes. The system will re-notify you every 3 minutes until you change the cartridge.

How should I respond?

Tap ok. Change your cartridge immediately by tapping OPTIONS from the Home Screen, then Load and follow the instructions in Section 5.3 Filling and Loading a t:slim Cartridge.

16.5 Cartridge Error Alarm

What will I see on the screen?

CARTRIDGE ALARM ALL DELIVERIES STOPPED! This cartridge cannot be used. Remove and replace with a new cartridge. OK

What does it mean?

Your t:slim X2 Pump detected that the cartridge could not be used and all deliveries have stopped. This can be caused by cartridge defect, not following the proper procedure to load the cartridge, or over filling the cartridge (with more than 300 units of insulin).

How will the System notify me?

3 sequences of 3 notes or 3 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes. The system will re-notify you every 3 minutes until you change the cartridge.

How should I respond?

Tap or . Change your cartridge immediately by tapping **OPTIONS** from the Home Screen, then **Load** and follow the instructions in section 5.3 Filling and Loading a t:slim Cartridge.

16.6 Cartridge Removal Alarm

What will I see on the screen?

ALL DELIVERIES STOPPED! The cartridge cannot be detected.

Press INSTALL to install a new cartridge or press CONNECT to reconnect the current cartridge.

CONNECT

INSTALL

What does it mean?

Your t:slim X2 Pump detected that the cartridge has been removed and all deliveries have stopped.

How will the System notify me?

3 sequences of 3 notes or 3 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes. The system will re-notify you every 3 minutes until you reconnect the current cartridge or change the cartridge.

How should I respond?

Tap CONNECT to reattach the current cartridge. Tap INSTALL to load a new cartridge.

16.7 Temperature Alarm

What will I see on the screen?

TEMPERATURE ALARM ALL DELIVERIES STOPPED! Remove pump from extreme temperatures and then resume insulin delivery. OK

What does it mean?

Your t:slim X2 Pump detected an internal temperature below 2°C (35°F) or above 45°C (113°F) or a battery temperature below 2°C (35°F) or above 52°C (125°F) and all deliveries have stopped.

How will the System notify me?

3 sequences of 3 notes or 3 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes. The system will re-notify you every 3 minutes until a temperature in the operating range is detected.

How should I respond?

Tap ox. Remove the pump from the extreme temperature and then resume insulin delivery.

16.8 Occlusion Alarm 1

What will I see on the screen?



What does it mean?

Your t:slim X2 Pump detected that insulin delivery is blocked and all deliveries have stopped. See Section 30.4 t:slim X2 Pump Performance Characteristics for more information on how long it can take the system to detect an occlusion.

How will the System notify me?

3 sequences of 3 notes or 3 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes. The System will re-notify you every 3 minutes until you resume insulin delivery.

How should I respond?

Tap ox. Check the cartridge, tubing, and infusion site for any sign of damage or blockage and correct the condition. To resume insulin, from the Options menu, tap RESUME INSULIN and tap RESUME to confirm.

■ NOTE: Occlusion During Bolus

If the occlusion alarm occurs during bolus delivery, after tapping a screen will appear letting you know how much of the requested bolus was delivered before the occlusion alarm. When the occlusion is cleared, some or all of the previously requested insulin volume may be delivered. Test your BG at the time of alarm and follow your healthcare provider's instructions for managing potential or confirmed occlusions.

16.9 Occlusion Alarm 2

What will I see on the screen?

OCCLUSION ALARM ALL DELIVERIES STOPPED! Insulin delivery may be blocked. Change your site and check your BG in 1-2 hours.

OK

What does it mean?

Your t:slim X2 Pump detected a second occlusion alarm shortly after the first occlusion alarm and all deliveries have stopped.

How will the System notify me?

3 sequences of 3 notes or 3 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes. The System will re-notify you every 3 minutes until you resume insulin delivery.

How should I respond?


Tap os . Change the cartridge, tubing, and infusion site to ensure proper delivery of insulin. Resume insulin after changing the cartridge, tubing, and infusion site.

■ NOTE: Occlusion During Bolus

If the second occlusion alarm occurs during bolus delivery, after tapping www, a screen will appear letting you know that the amount of bolus delivery could not be determined and was not added to your IOB.

16.10 Screen On/Quick Bolus Button Alarm

What will I see on the screen?

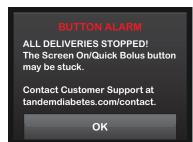

What does it mean?

The Screen On/Quick Bolus Button (on the top of your t:slim X2 Pump) is stuck or not functioning properly and all deliveries have stopped.

How will the System notify me?

3 sequences of 3 notes or 3 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes. The System will re-notify you every 3 minutes until the condition is corrected.

How should I respond?

Tap ok. Contact your local Tandem Diabetes Care® representative.

16.11 Altitude Alarm

What will I see on the screen?

ALTITUDE ALARM ALL DELIVERIES STOPPED! Remove cartridge from pump, reconnect cartridge and then resume insulin. OK

What does it mean?

Your t:slim X2 Pump detected a difference between inside the cartridge and the surrounding air within the validated operating range of -400 meters to 3,000 meters and all deliveries have stopped.

How will the System notify me?

3 sequences of 3 notes or 3 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes. The System will re-notify you every 3 minutes until the condition is corrected.

How should I respond?

Tap ox. Remove the cartridge from the pump (this will allow the cartridge to fully vent) and then reconnect the cartridge.

16.12 Reset Alarm

What will I see on the screen?

PUMP HAS BEEN RESET All active deliveries have been stopped and your IOB and Max Hourly Bolus have been reset. Contact Customer Support at tandemdiabetes.com/contact.

What does it mean?

Your t:slim X2 Pump detected that one if its micro-processors experienced a reset and all deliveries have been stopped.

How will the System notify me?

3 sequences of 3 notes or 3 vibrations depending on the volume/vibrate setting selected in Sound Volume.

Will the System re-notify me?

Yes. The System will re-notify you every 3 minutes until you tap ox.

How should I respond?

Tap ox . Contact your local Tandem Diabetes Care representative.

Chapter 17

t:slim X2 Pump Malfunction

Chapter 17 – t:slim X2 Pump Malfunction

17.1 Malfunction

If your t:slim X2 Pump detects a system error, the MALFUNCTION screen appears and all deliveries are stopped. Contact your local Tandem Diabetes Care® representative.

Malfunctions notify you with 3 sequences of 3 notes at highest volume and 3 vibrations. They repeat at regular intervals until acknowledged by tapping SILENCE ALARM.

A PRECAUTION

ALWAYS check with your healthcare provider for specific guidelines if you want or need to disconnect from the pump for any reason. Depending on the length of time and reason you are disconnecting, you may need to replace missed basal and/or bolus insulin. Check your blood glucose before disconnecting from the pump and again when you reconnect, and treat high blood glucose (BG) levels as recommended by your healthcare provider.

What will I see on the screen?



What does it mean?

Your t:slim X2 Pump detected a system error and all deliveries have been stopped.

How will the System notify me?

3 sequences of 3 notes at highest volume and 3 vibrations.

Will the System re-notify me?

Yes, the System will notify you every 3 minutes until you acknowledge the malfunction by tapping SILENCE ALARM.

How should I respond?

- Write down the Malfunction Code number that appears on the screen.
- Tap SILENCE ALARM. The MALFUNCTION screen will remain on the pump even though the alarm is silenced.
- Contact your local Tandem Diabetes Care representative and provide the Malfunction Code number that you wrote down.

Section 5

Getting Started with Your CGM System

Chapter 18

Dexcom G5 Mobile CGM Overview

18.1 CGM System Overview

This section of the User Guide covers instructions for using Dexcom G5 Mobile Continuous Glucose Monitoring (CGM) with your t:slim X2 Pump. Use of Dexcom G5 Mobile CGM is optional and when used, allows readings from your sensor to be displayed on your t:slim X2 Pump screen. Dexcom G5 Mobile CGM uses two parts, a sensor and a transmitter. You will also need a commercially available blood glucose meter to use with your System.

The Dexcom G5 Mobile Sensor is a disposable device that is inserted under the skin to continuously monitor glucose levels for up to 7 days. The Dexcom G5 Mobile Transmitter connects to the sensor pod and wirelessly sends readings to the pump display every 5 minutes. The display shows sensor glucose readings, trend graph, direction and rate of change arrows. You can also program your pump to alert you when your CGM readings are above or below a given level, or are rising or falling quickly.

Unlike the readings from a standard blood glucose meter, CGM readings allow you to view trends in real time, as well as capture information when you would otherwise be unable to check your blood sugar, such as while you are asleep. This information can be useful for you and your healthcare provider when considering changes to your therapy. In addition, the programmable alerts can help you to spot potential low or high glucose sooner than you would using a only a blood glucose meter.

■ NOTE: Device Connections

The Dexcom G5 Mobile CGM only allows pairing with one medical device at a time (either the t:slim X2 Pump or the Dexcom Receiver), but you can still use the Dexcom Mobile App and your t:slim X2 Pump simultaneously using the same transmitter ID.


18.2 Sensor Overview

For your safety, the sensor is packaged in a sterile sealed pack. When you first open the pack, your sensor looks like one item; however, it's actually three: sensor applicator, sensor pod, and sensor wire.

The applicator helps you insert the sensor wire inside the sensor pod under your skin. After inserting the sensor wire, remove the applicator. The sensor wire stays in the sensor pod with the pod attached to your skin by adhesive; the sensor wire is made of silver and platinum with polymer membranes. Once inserted, the thin and flexible wire measures your glucose levels in the fluid between your cells (interstitial fluid) for up to seven days.

Dexcom G5 Mobile Sensor Diagram

- 1. Plunger
- 2. Applicator Barrel
- 3. Transmitter Latch
- 4. Release Tab
- 5. Safety Lock
- 6. Collar
- 7. Sensor Pod
- 8. Adhesive Patch

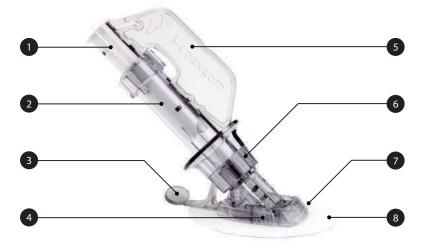

18.3 Receiver (t:slim X2 Pump) Overview

- 1. Time and Date Display: Displays the current time and date.
- Antenna: Indicates communication status between pump and transmitter.
- Battery Level: Displays the level of battery power remaining. When connected for charging, the charging icon (lightning bolt) will display.
- 4. High Glucose Alert Setting
- Glucose Target Range
- Low Glucose Alert Setting
- Plot of Most Recent Sensor Glucose Readings
- 8. Options: Stop/Resume insulin delivery, manage Pump and CGM Settings, program a Temp Rate, Load cartridge, and view History.

- 9. Bolus: Program and deliver a bolus.
- 10. Status: Displays current system settings and insulin delivery status.
- 11. **Insulin Level:** Displays the current amount of insulin in the cartridge.
- 12. Most Recent 5-Minute Glucose Reading
- 13. **Trend Arrow:** Indicates direction and rate of change.
- 14. Trend Graph Time (HRS): 1, 3, 6, 12 and 24 hour views available.
- Insulin On Board (IOB): Amount and time remaining of any active insulin on board.


To view CGM information on the full screen:

From the Home Screen tap anywhere on the CGM trend graph.



Tap the "minimize" icon to return to the home screen.



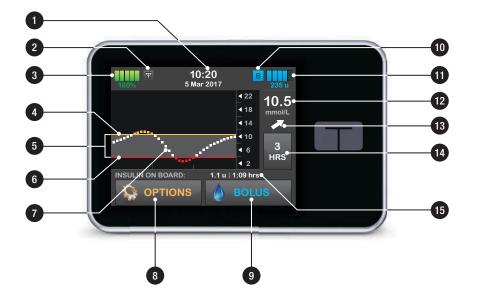

18.4 Transmitter Overview

Snapping into the sensor pod, the transmitter wirelessly sends your glucose information to your display devices—receiver and/or smart device. If you have a new transmitter, only open the package when you are ready to use it.

Transmitter features:

- Reusable
 - Do not discard after sensor session.
 - Only for you, don't share transmitter.
- Water resistant
- Can transmit data to your display devices for up to 6 meters
 - Range is less if you are in or under water.
- Battery lasts approximately 90 days

- Receiver or smart devices prompts you when battery is running low.
- Serial number is on the back

▲ PRECAUTION

AVOID separating the transmitter and pump by more than 6 meters. The transmission range from the transmitter to the pump is up to 6 meters without obstruction. Wireless communication does not work well through water so the range is much less if you are in a pool, bathtub, or on a water bed, etc. Types of obstruction differ and have not been tested. If your transmitter and pump are farther than 6 meters apart or are separated by an obstruction, they might not communicate or the communication distance may be shorter and result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.



The transmitter battery will last at least 90 days. Once you see the Low Transmitter Battery Alert, replace the transmitter as soon as possible. Your transmitter battery may drain as quickly as 7 days after this alert occurs.



Chapter 19

CGM Settings

19.1 About Bluetooth

Bluetooth® Low Energy (BLE) is a type of wireless communication used in cell phones and many other devices. Your t:slim X2 Pump wirelessly pairs with other devices, including the Dexcom G5 Mobile Transmitter, using BLE communication. This allows the pump and paired devices to communicate securely and only with each other.

19.2 Disconnecting from the Dexcom Receiver

The Dexcom G5 Mobile CGM only allows pairing with one medical device at a time. Ensure your CGM transmitter is not connected to the Dexcom Receiver before pairing with the pump by doing the following:

Before entering your CGM transmitter ID into the t:slim X2 Insulin Pump, turn off the Dexcom Receiver and wait 15 minutes. This allows the Dexcom G5

Transmitter to forget the connection currently in place with the Dexcom G5 Receiver.

■ NOTE: Turn Receiver Off

It is not enough to Stop the Sensor Session on your Dexcom receiver prior to pairing to the t:slim X2 Pump. The receiver power must be completely off in order to avoid connection problems.

You may still use a smartphone with the Dexcom Mobile App and your t:slim X2 Insulin Pump simultaneously with the same transmitter ID.

19.3 Entering Your Transmitter ID

To activate the BLE communication, you need to enter the unique transmitter ID into your pump. Once the transmitter ID has been entered into your pump, the two devices can be paired, allowing your sensor glucose readings to be displayed on your t:slim X2 Pump.

If you need to replace your transmitter, you will need to enter the new transmitter ID into your pump. If you need to

replace your pump, you will need to reenter the transmitter ID into your pump.

- 1. Remove the transmitter from its packaging.
- 2. From the Home Screen, tap OPTIONS.
- 3. Tap My CGM.
- 4. Tap CGM Settings.
- 5. Tap Transmitter ID.
- 6. Using the onscreen keypad, enter the unique transmitter ID.

The transmitter ID can be found on the bottom (flat surface) of your transmitter.

The letters I, O, V, and Z are not used in transmitter IDs and should not be entered. If one of these letters is entered, you will be notified that an invalid ID was entered and prompted to enter a valid ID.

7. Tap ____

- 8. To make sure that the correct transmitter ID is entered, you will be prompted to enter it a second time.
- 9. Repeat step 6 above, then tap ...



If the transmitter IDs you entered do not match you will be prompted to start the process again.

✓ Once matching values have been entered, you will be returned to the CGM Settings screen and the transmitter ID vou entered will be highlighted in yellow.

10. Tap _____.

19.4 Setting CGM Volume

You can set the sound pattern and volume for CGM alerts and prompts to meet your individual needs. Reminders, alerts, and alarms for insulin delivery functions are separate from alerts and errors for CGM functions and do not follow the same pattern and volume.

To set your Sound Volume, see Section 11.3 t:slim X2 Sound Volume.

CGM Volume options:

Vibrate

When you want to silence CGM volume and be alerted by vibration. The only exception to this is the fixed low alert at 3.1 mmol/L, which alerts you as a vibration first, followed by beeps 5 minutes later if not confirmed.

Soft

When you want your alert to be less noticeable. This sets all alerts and alarms to lower volume beeps.

Normal

The default profile when you receive vour System. This sets all alerts and alarms to higher volume beeps.

HypoRepeat

Very similar to normal profile, but it continuously repeats the fixed low alert every 5 seconds until your sensor glucose reading rises above 3.1 mmol/L or is confirmed. This can be helpful if you want extra alerts for severe low sensor glucose readings.

The CGM Volume setting that you choose applies to all CGM alerts. errors, and prompts which have their own unique sound pattern, tone and volume. This allows you to identify each alert and error and its meaning.

The fixed low alert at 3.1 mmol/L cannot be turned off or changed.

The Soft, Normal, and HypoRepeat options have the following sequence:

- The first alert is vibrate only.
- If the alert is not confirmed in 5 minutes, the system vibrates and beeps.
- If the alert is not confirmed in 5 more minutes, the system vibrates and beeps louder. This continues at the same volume every 5 minutes until confirmed.

For the HypoRepeat Option Only:

If the alert is confirmed and your sensor glucose readings continue to be at or

Chapter 19 – CGM Settings

below 3.1 mmol/L your system repeats the alert sequence in 30 minutes.

CGM Volume Details

CGM Volume	Vibrate	Soft	Normal	HypoRepeat
High Alert	2 long vibrates	2 long vibrates + 2 low beeps	2 long vibrates + 2 medium beeps	2 long vibrates + 2 medium beeps
Low Alert	3 short vibrates	3 short vibrates + 3 low beeps	3 short vibrates + 3 medium beeps	3 short vibrates + 3 medium beeps
Rise Alert	2 long vibrates	2 long vibrates + 2 low beeps	2 long vibrates + 2 medium beeps	2 long vibrates + 2 medium beeps
Fall Alert	3 short vibrates	3 short vibrates + 3 low beeps	3 short vibrates + 3 medium beeps	3 short vibrates + 3 medium beeps
Out of Range Alert	1 long vibrate	1 long vibrate + 1 low beep	1 long vibrate + 1 medium beep	1 long vibrate + 1 medium beep
Fixed Low Alert	4 short vibrates + 4 medium tone beeps + pause + repeat sequence			
All Other Alerts	1 long vibrate	1 long vibrate + 1 low beep	1 long vibrate + 1 medium beep	1 long vibrate + 1 medium beep

Chapter 19 – CGM Settings

To Select Your CGM Volume:

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My CGM.
- 3. Tap CGM Settings.
- Tap CGM Volume to change your setting. The default setting is Normal.
- 5. Tap Vibrate, Soft, Normal or HypoRepeat to select.
- ✓ Once a value is selected, the pump will return to the previous screen.
- 6. Tap <

19.5 CGM Info

CGM info contains important information about your device. The following can be found in CGM info:

- Firmware Revision
- Hardware Revision

You can view this information at any time.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My CGM.
- 3. Tap CGM Settings.
- 4. Tap CGM Info.

Chapter 20

Setting CGM Alerts

Setting Your CGM Alerts

You can create personal settings for how and when you want the System to tell you what is happening.

The High and Low Alerts tell you when your sensor glucose readings are outside your target glucose range.

Rise and Fall (rate of change) Alerts let you know when your glucose levels are changing fast.

The System also has a 3.1 mmol/L Fixed Low Alert that cannot be changed or turned off. This safety feature tells you your glucose level may be dangerously low.

The Out of Range Alert notifies you when the transmitter and pump are not communicating. Keep the transmitter and the pump within 6 meters of each other without obstruction. When the transmitter and the pump are too far apart, you will not get sensor glucose readings or alerts.

High and Low Glucose Alerts

You can personalize the High and Low Alerts which tell you when your sensor glucose readings are outside of your target glucose range. When you have both your High and Low Alerts turned on, a grey zone on your trend graph shows your target range. The default for the High Alert is on, 11.1 mmol/L. The default for the Low Alert is on, 4.4 mmol/L. Consult with your healthcare provider before setting the High and Low Glucose Alert setting.

20.1 Setting Your High Glucose Alert and Repeat Feature

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My CGM.
- Tap CGM Alerts.
- 4. Tap High and Low.
- 5. To Set the High Alert, tap **High** Alert.

6. Tap Alert Me Above.

The default setting for the High Alert is 11.1 mmol/L.

NOTE: Turning the Alert Off To turn off the High Alert, tap OFF/ON. The

To turn off the High Alert, tap OFF/ON. The screen will indicate that OFF is selected.

- Using the onscreen keypad, enter the value above which you want to be notified. It can be set between 6.7 and 22.2 mmol/L in 0.1 mmol/L increments.
- 3. Tap 🕶.

The repeat feature allows you to set a time for the High Alert to sound again and display on your pump as long as your sensor glucose reading remains above the High Alert value. The default value is: Never (the alert will not sound again). You can set the repeat feature to sound again every 15 minutes, 30 minutes, 1 hour, 2 hours, 3 hours, 4 hours, or 5 hours when your sensor glucose reading remains above the High Alert value.

To Set Up the Repeat Feature:

- 9. Tap Repeat.
- 10. To select the repeat time, tap the time you want the alert to sound again. For instance, if you select 1 hr, the alert will sound every hour as long as your sensor glucose reading remains above the High Alert value.

Use the up and down arrows to view all Repeat options.

- ✓ Once a value is selected, the pump will return to the previous screen.
- 11. Tap 💟.

20.2 Setting Your Low Glucose Alert and Repeat Feature

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My CGM.
- 3. Tap CGM Alerts.
- 4. Tap High and Low.

- 5. To Set the Low Alert, tap Low Alert.
- 6. Tap Alert Me Below.

The default setting for the Low Alert is 4.4 mmol/L.

NOTE: Turning the Alert Off
To turn off the Low Alert, tap OFF/ON. The
screen will indicate that OFF is selected.

- Using the onscreen keypad, enter the value below which you want to be notified. It can be set between 3.3 and 5.6 mmol/L in 0.1 mmol/L increments.
- Tap _____.

The repeat feature allows you to set a time for the Low Alert to sound again and display on your pump as long as your sensor glucose reading remains below the Low Alert value. The default value is: Never (the alert will not sound again). You can set the repeat feature to sound again every 15 minutes, 30 minutes, 1 hour, 2 hours, 3 hours, 4 hours, or 5 hours when your sensor glucose

reading remains below the Low Alert value.

To Set Up the Repeat Feature:

- 9. Tap Repeat.
- 10. To select the repeat time, tap the time you want the alert to sound again. For instance, if you select 1 hr, the alert will sound every hour as long as your sensor glucose reading remains below the Low Alert Value.

Use the up and down arrows to view all repeat options.

- ✓ Once a value is selected, the pump will return to the previous screen.
- 11. Tap 🛶.

20.3 Rate Alerts

Rate alerts tell you when your glucose levels are rising (Rise Alert) or falling (Fall Alert) and by how much. You can choose to be alerted when your sensor glucose reading is rising or falling 0.11 mmol/L or more per minute, or 0.17

Chapter 20 - Setting CGM Alerts

mmol/L or more per minute. The default value for both the Fall Alert and the Rise Alert is off. When turned on, the default is 0.17 mmol/L. Consult with your healthcare provider before setting the Rise and Fall Alerts.

Examples

If you set your Fall Alert to 0.11 mmol/L per minute and your sensor glucose readings fall at this rate or faster, the CGM FALL ALERT with one arrow pointing down shows. The pump vibrates or beeps according to your CGM Volume selection.



If you set your Rise Alert to 0.17 mmol/ L per minute and your sensor glucose readings rise at this rate or faster, the CGM RISE ALERT with two arrows pointing up shows. The pump vibrates or beeps according to your CGM Volume selection.



20.4 Setting Your Rise Alert

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My CGM.
- 3. Tap CGM Alerts.
- 4. Tap Rise and Fall.
- 5. Tap Rise Alert.
- 6. To select the default of 0.17 mmol/L/min, tap

To change your selection, tap Rate.

■ NOTE: Turning the Alert Off

To turn off the Rise Alert, tap OFF/ON.

- 7. Tap 0.11 mmol/L/min to select.
- Once a value is selected, the pump will return to the previous screen.
- 8. Tap 💙.

20.5 Setting Your Fall Alert

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My CGM.
- 3. Tap CGM Alerts.
- 4. Tap Rise and Fall.
- 5. Tap Fall Alert.
- 6. To select the default of 0.17 mmol/ L/min, tap ...

To change your selection, tap Rate.

NOTE: Turning the Alert Off To turn off the Fall Alert, tap OFF/ON.

- 7. Tap 0.11 mmol/L/min to select.
- ✓ Once a value is selected, the pump will return to the previous screen.
- 8. Tap <u></u>

20.6 Setting Your Out of Range Alert

The range from the transmitter to the pump is up to 6 meters without obstruction.

The Out of Range Alert lets you know when your transmitter and pump are not communicating with each other. Keep the transmitter and the pump within 6 meters of each other without obstruction. When the transmitter and pump are not communicating, you will not get sensor glucose readings or alerts. The default value is off; 30 minutes if on.

The Out Of Range symbol appears on the pump Home Screen and on the Out

of Range Alert screen (if turned on) when the transmitter and pump are not communicating. The amount of time out of range also shows on the alert screen. It will continue to re-alert until the transmitter and pump are back in range.

To Set Your Out of Range Alert:

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My CGM.
- 3. Tap CGM Alerts.
- 4. Tap Out of Range.

The default value is off; 30 minutes if on.

5. To select the default of 30 min, tap

To change the time, tap Alert After.

6. Using the onscreen keypad, enter the time after which you want to be alerted (between 20 minutes and 3

hours and 20 minutes) then tap

7. Tap 🛶.

THIS PAGE IS INTENTIONALLY LEFT BLANK

Chapter 21

Inserting the Sensor and Placing the Transmitter

Chapter 21 – Inserting the Sensor and Placing the Transmitter

To use continuous glucose monitoring (CGM) with your t:slim X2 Pump, you need a Dexcom G5 Mobile Sensor and Transmitter, a blood glucose meter, and test strips for calibration. After insertion and calibration, the sensor continuously measures and displays your sensor glucose readings for up to 7 days. The following sections will show you how to insert the sensor, attach the transmitter, and start a new CGM session.

A WARNING

DO NOT ignore sensor fractures. Sensors may fracture on rare occasions. If a sensor breaks and no portion of it is visible above the skin, do not attempt to remove it. Seek professional medical help if you have symptoms of infection or inflammation—redness, swelling or pain—at the insertion site. If you experience a broken sensor, please report this to your local Tandem Diabetes Care representative. For patients undergoing an MRI with a retained wire broken off from a Dexcom G5 Mobile Sensor, in-vitro MRI testing did not detect any safety hazards. There was no significant migration or heating of the wire and imaging artifacts were limited to the area around the wire.

21.1 Before You Start

- Wash your hands thoroughly.
- Make sure that the date and time on your t:slim X2 Pump are correct.
- Make sure that the correct transmitter ID has been entered into your pump (see Section 19.1 About Bluetooth). You only need to enter a new transmitter ID if the pump or transmitter have been replaced.
- Check the expiration date on the sensor package label. The format is YYYY-MM-DD. Insert sensors on or before the end of the expiration date.
- Follow your blood glucose meter manufacturer's instructions to make sure you are getting accurate blood glucose values for calibration and treatment decisions (check code and use control solution).
- Wipe the bottom of the transmitter with a damp cloth or alcohol wipe.
 Place the transmitter on a clean, dry

- cloth to dry before you start a new sensor session.
- Make sure that the time on your pump and your CGM are the same.
- Closely inspect pouch, check it has not been damaged. If yes, do not use
- Remove sensor applicator from pouch.
- Closely inspect sensor, check it has not been damaged. If yes, do not use.
- Keep sensor packaging until sensor session is complete.

A WARNING


STORE the sensor at temperatures between 2°C (36°F) and 25°C (77°F) for the length of the sensor's life. You may store the sensor in the refrigerator if it is within this temperature range. The sensor should not be stored in a freezer. Storing the sensor improperly might cause the sensor glucose readings to be inaccurate, and you might miss severe hypoglycemia (low blood

glucose) or hyperglycemia (high blood glucose) events.

Chapter 21 – Inserting the Sensor and Placing the Transmitter

Review the CGM Sensor Applicator Picture Before Using a New Sensor

- 1. Plunger
- 2. Applicator Barrel
- 3. Transmitter Latch
- 4. Release Tab
- 5. Safety Lock
- 6. Collar
- 7. Sensor Pod
- 8. Adhesive Patch

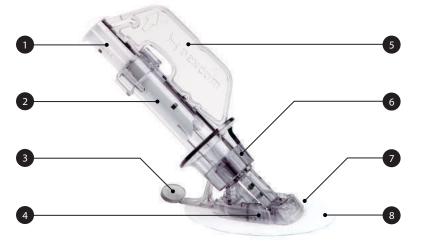

21.2 Removing the Sensor From its Packaging

A PRECAUTION

DO NOT use the sensor if its sterile package has been damaged or opened. Using an unsterile sensor might cause infection.

- Wash your hands thoroughly, and dry them.
- Carefully remove the sensor from its packaging. Look closely at the sensor to make sure it is not damaged.
- The applicator is for single use and is disposable.
- The safety lock prevents you from releasing the needle accidentally before you are ready.

21.3 Choosing an Insertion Site

A PRECAUTION

AVOID inserting the sensor in areas that are likely to be bumped, pushed or compressed or areas of skin with scarring, tattoos, or irritation

as these are not ideal sites to measure glucose. Insertion in those areas might affect accuracy and could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

A PRECAUTION

AVOID injecting insulin or placing an infusion set within 8 centimeters of the sensor. The insulin might affect sensor accuracy and could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

For adults (ages 18+) only sites on the belly are approved for sensor insertion. For children (ages 6–17) both sites on the belly and upper buttocks are approved for sensor insertion.

No other sensor insertion sites have been tested and it is not known how well the sensor will work in other sites. The best areas to insert your sensor are usually flat, "pinchable," and free from where rubbing can occur, such as along the waist band or seat belt strap.

• Choose an area at least 8 centimeters from your insulin pump infusion set or injection site.

- You may need to shave the area where you plan to put the sensor so the adhesive patch sticks securely.
- Make sure there are no traces of lotions, perfumes, or medications on the area.

A PRECAUTION


AVOID using the same spot repeatedly for sensor insertion. Rotate your sensor placement sites, and do not use the same site for 2 sensor sessions in a row. Using the same site might cause scarring or skin irritation.

21.4 Placing the CGM Sensor

- Clean your skin at the sensor placement site with an alcohol wipe. Make sure the area is clean and completely dry before you insert the sensor.
- 2. Remove the adhesive backing from the sensor pod one half at a time, using the white tabs on the backing. Hold the sensor by the applicator

Chapter 21 – Inserting the Sensor and Placing the Transmitter


barrel, and try not to touch the sticky adhesive patch.

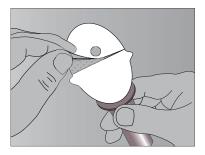


 Place the sensor horizontally (side to side) as shown in the picture. Do not place the sensor vertically (up and down) on your skin.



- 4. Move your fingers around the adhesive patch to remove any wrinkles and secure the tape firmly to your skin.
- Hold the applicator and pull the safety lock straight out away from the applicator, in the direction of the arrows in the picture.

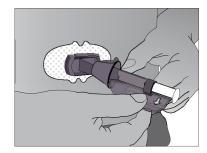

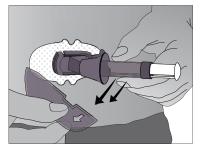


 Save the safety lock to help you remove the transmitter at the end of your sensor session. The safety lock can be used for transmitter removal but is not required.


21.5 Inserting the CGM Sensor

After you place the applicator on your skin and remove the safety lock, you are ready to insert the sensor. Follow the steps below to insert your sensor.


- Place the fingers of one hand at the edge of the white adhesive (at the opposite side of the sensor from the transmitter latch). You may pinch up on the skin using this hand. Do not pinch up in the middle section of the plastic base.
- While still pinching, use your other hand to place two fingers above the collar on the applicator barrel so they are resting above the collar.

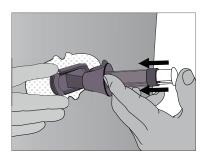


 Place your thumb on the white plunger. Push the plunger down completely, making sure it is flush against the applicator barrel. You should hear 2 clicks. This inserts the needle and sensor under your skin. When you are pushing down on the plunger, do not pull back on the collar.

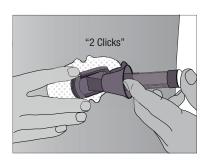

4. Keep pinching up on your skin with one hand. Place two fingers under the collar. Keep your thumb lightly on top of the white plunger, and pull the collar back towards your thumb until you hear 2 clicks or cannot pull back any more. This leaves the sensor under your skin and removes the needle from your body.

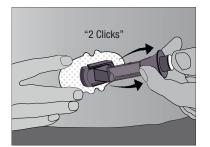

5. Squeeze the center of the ribbed release tabs on the sides of the sensor pod to remove the applicator barrel. Only the sensor pod will be left on your body.

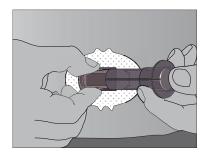

- Make sure the transmitter latch is down (against your body) before squeezing the tabs to remove the applicator barrel.
- Squeeze the center of the ribbed part of the release tabs.
- While squeezing the tabs, rock the applicator barrel forward and out away from your body.

If you have any problems with insertion, save the sensor and applicator and contact your local Tandem Diabetes Care representative.

Chapter 21 – Inserting the Sensor and Placing the Transmitter


21.6 Attaching the CGM Transmitter

The CGM transmitter wirelessly sends your sensor glucose readings to your t:slim X2 Pump. You must snap the transmitter into the sensor pod after you insert your sensor. Follow the steps below to attach your transmitter.

- Wipe and dry the bottom of the transmitter with a damp cloth or an alcohol wipe before every use.
 - Do not touch the metal circles on the bottom of the transmitter with your skin.
 - Do not scratch the bottom of the transmitter as scratches may compromise the waterproof seal.
- Place the transmitter in the sensor pod with the flat side down, and the narrower side away from the transmitter latch.



3. Snap in the transmitter:

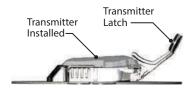

- You may want to pinch up on your skin at the front edge of the white adhesive with one hand.
- Place one finger on the transmitter to keep it in place.

- With your other hand, pull the latch up and forward, over the transmitter to snap it into place.
 The transmitter should lie flat in the sensor pod.
- Make sure you hear 2 clicks when you snap the transmitter in place. If it is not fully snapped in, this may lead to a poor connection, allowing fluid to get under the transmitter. This can lead to inaccurate sensor glucose readings.
- Release your pinch on the adhesive edge at this time.
- Make sure the transmitter is secure by sliding your fingers under each long side of the sensor pod and pressing down on the transmitter with your thumb of the same hand, like you are pinching it.
- Remove the transmitter latch by holding the sides of your sensor pod with one hand and quickly twisting off the latch away from your body with your other hand.

5. Do not remove the transmitter from the sensor pod while the pod is attached to your skin.

The sensor pod should stay on your skin using its own adhesive. But, if the patch is peeling up, you can use medical tape for extra support. If you use tape, only tape over the white adhesive patch on all sides for even support. Do not tape over the transmitter or any of the plastic parts of the sensor pod. Do not tape under the sensor pod or leave any substance on the skin where you insert the sensor.



Chapter 21 – Inserting the Sensor and Placing the Transmitter

THIS PAGE IS INTENTIONALLY LEFT BLANK

Chapter 22

Starting a CGM Sensor Session

22.1 Starting Sensor

Once you have entered your transmitter ID, inserted your sensor, and attached your transmitter, you are ready to start a CGM session. Follow the steps below.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My CGM.
- 3. Tap START SENSOR.
- Once you start a sensor session, the START SENSOR option is replaced with STOP SENSOR.
- 4. Tap to confirm.
- The "SENSOR STARTED" screen will appear to let you know your 2hour sensor startup has begun.
- ✓ Your t:slim X2 will return to the CGM Home Screen with the 3 hour trend graph displayed.

- Check your t:slim X2 Pump CGM
 Home Screen 10 minutes after
 starting your sensor session to
 make sure your pump and
 transmitter are communicating. The
 antenna symbol should be to the
 right of the battery indicator and
 should be white.
- If you see the out of range symbol below the insulin level indicator, and the antenna symbol is grayed out, follow these troubleshooting tips:
 - a. Make sure your t:slim X2 Pump and transmitter are within 6 meters of each other without obstruction. Re-check in 10 minutes to see if the out of range symbol is still active.
 - If the pump and transmitter are still not communicating, check the CGM Info screen to make sure the correct transmitter ID is entered.
 - c. If the correct transmitter ID is entered and the pump and transmitter are still not communicating, contact your local Tan-

dem Diabetes Care representative.

22.2 Sensor Startup Period

The sensor needs a 2-hour startup period to adjust to being under your skin. You will not get sensor glucose readings or alerts until the 2-hour startup period ends and you complete your first calibrations.

During the startup period, the CGM Home Screen on your t:slim X2 Pump shows a 2-hour countdown symbol in the upper right portion of the screen. The countdown symbol fills in over time to show that you are getting closer to the first calibration time.



■ NOTE:

Once your sensor session has completed, you must end the current session before you can start a new session.

Examples

For example, if you started your sensor session 20 minutes ago, you would see this countdown symbol on the CGM Home Screen.



If you started your sensor session 90 minutes ago, you would see this count-down symbol on the CGM Home Screen.



At the end of the 2-hour startup period, you will be prompted to enter 2 calibration values and 2 blood drops will appear in the place where the countdown symbol was.



Follow the instructions in the next section to calibrate your sensor.

■ NOTE:

Once your sensor session has completed, you must end the current session before you can start a new session.

Chapter 23

Calibrating Your CGM System

23.1 Calibration Overview

Your CGM System needs calibrations using blood glucose (BG) values obtained from a commercially available blood glucose meter to display continuous sensor glucose readings and trend information. There are important times when you must calibrate:

- 2-hour startup: 2 hours after you start your sensor session
- 12-hour update: every 12 hours after the 2 hour start up calibration
- More information is needed for other reasons

On the first day of your sensor session, you must enter 2 blood glucose values into your pump to calibrate. You must enter 1 blood glucose value to calibrate every 12 hours after your first startup calibration. The pump will remind you when the System needs these calibrations. In addition, you may be prompted to enter additional blood glucose values to calibrate as needed.

▲ CONTRAINDICATION

Taking medications with acetaminophen/paracetamol while wearing the sensor may falsely raise your sensor glucose readings. The level of inaccuracy depends on the amount of acetaminophen/paracetamol active in your body and may be different for each person.

A WARNING

CALIBRATE your CGM at least once every 12 hours. Calibrating less often than every 12 hours might cause sensor glucose readings to be inaccurate and glucose alerts to become unreliable. This could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

A PRECAUTION

DO NOT calibrate if your blood glucose is changing at a significant rate, typically more than 0.11 mmol/L per minute. Do not calibrate when your receiver screen is showing the rising or falling single arrow or double arrow, which indicates that your blood glucose is rapidly rising or falling. Calibrating during significant rise or fall of blood glucose may affect sensor accuracy and could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

When calibrating, you must enter your blood glucose values into the pump by hand. You can use any commercially available blood glucose meter. You must calibrate with accurate blood glucose meter values to get accurate sensor glucose readings.

Follow these important instructions when obtaining blood glucose values for calibration:

- Blood glucose values used for calibration must be between 2.2 to 22.2 mmol/L and must have been taken within the past 5 minutes.
- Your sensor cannot be calibrated if the glucose value from your meter is less than 2.2 mmol/L. For safety reasons, if your blood glucose is low, first treat your low blood glucose.
- Make sure a sensor glucose reading shows in the upper right portion of the CGM Home Screen before calibrating.
- Make sure the antenna symbol is visible to the right of the battery

indicator on the CGM Home Screen and is active (white, not grayed out) before calibrating.

- Always use the same meter to calibrate that you routinely use to measure your blood glucose. Do not switch your meter in the middle of a sensor session. Blood glucose meter and strip accuracy vary between blood glucose meter brands.
- The accuracy of the blood glucose meter used for calibration may affect the accuracy of sensor glucose readings. Follow your blood glucose meter manufacturer's instructions for blood glucose testing.

23.2 Startup Calibration

2 hours after you start the sensor session, the CALIBRATE CGM screen will appear, letting you know that 2 separate blood glucose values from your meter must be entered. You will not see sensor glucose readings until the pump accepts the blood glucose values.

- 1. From the Calibrate CGM Screen, Tap ok.
- The CGM Home Screen will appear with two blood drops in the upper right portion of the screen. The two blood drops will stay on the screen until you enter 2 separate blood glucose values to calibrate.
- Wash and dry your hands, make sure your glucose test strips are not expired and have been stored properly, and make sure your meter is properly coded (if required).
- Take a blood glucose measurement using your meter. Carefully apply the blood sample to the test strip following your meter manufacturer's instructions.

APRECAUTION

DO NOT use alternative blood glucose site testing (blood from your palm or forearm, etc.) for calibration. Alternative site blood glucose values may be different than those taken from a fingerstick blood glucose value and may not represent the timeliest blood glucose value. Use a blood glucose value taken only from a fingerstick for cali-

bration. Alternative site blood glucose values might affect sensor accuracy and result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

- 4. Tap OPTIONS.
- Tap My CGM.
- Tap Calibrate CGM.
- Using the onscreen keypad, enter the blood glucose value from your meter.

APRECAUTION

To calibrate the System, **D0** enter the exact blood glucose value that your blood glucose meter displays within 5 minutes of a carefully performed blood glucose measurement. Do not enter sensor glucose readings for calibration. Entering incorrect blood glucose values, blood glucose values obtained more than 5 minutes before entry, or sensor glucose readings might affect sensor accuracy and could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

8. Tap 🗸

Chapter 23 – Calibrating Your CGM System

- 9. Tap to confirm the calibration.
 - Tap x if the BG value does not exactly match the reading from your meter. The onscreen keypad will reappear. Enter the exact reading from your meter.
- √ The "CALIBRATION ACCEPTED" screen will appear.
- ✓ The My CGM screen will appear.
- 10. Tap Calibrate CGM to enter your second blood glucose value.
- ✓ The onscreen keypad will appear.
- 11. Wash and dry your hands, make sure your glucose test strips are not expired and have been stored properly, and make sure your meter is properly coded (if required).
- 12. Take a blood glucose measurement using your meter. Carefully apply the blood sample to the test strip following your meter manufacturer's instructions.

- 13. Follow steps 8 –10 to enter your second blood glucose value.
- Carefully apply the blood sample to the test strip following your meter manufacturer's instructions.
- 15. Follow steps 9 –10 to enter your second blood glucose value.

23.3 Calibration Blood Glucose Value and Correction Bolus

Your t:slim X2® Pump uses the blood glucose value entered for calibration to determine if a correction bolus is needed, or to provide other important information about your insulin on board and blood glucose.

 If you enter a calibration value that is above your Target BG in Personal Profiles, a message screen will indicate "YOUR BG IS ABOVE TARGET". To add a correction bolus, tap ... Follow the instructions in Section 7.2 BG Entry and Correction Bolus Calculation to deliver a correction bolus.

- If you enter a calibration value that is below your Target BG in Personal Profiles, a message screen will indicate "YOUR BG IS BELOW TARGET", and other important information will appear on the screen.
- If you enter your Target BG as a calibration value, the pump will return to the CCM Home Screen.

23.4 12 Hour Calibration Update

Calibrate your CGM System at least every 12 hours after your first calibration (2-hour startup calibration) to make sure your sensor glucose readings remain accurate and close to your blood glucose values. You can enter blood glucose values sooner than 12 hours if you want. If you have not entered any blood glucose values in the past 12 hours, the pump will ask you to enter a blood glucose value to update its calibration.

A WARNING

CALIBRATE your CGM at least once every 12 hours. Calibrating less often than every 12

hours might cause sensor glucose readings to be inaccurate and glucose alerts to become unreliable. This could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

2 hours after you start the sensor session, and every 12 hours after that, the "CALIBRATE CGM" screen will appear, letting you know that a blood glucose value from your meter must be entered to calibrate. In addition, a blood drop will appear to the right of the antenna symbol and will remain there until a blood glucose value is entered to calibrate.

1. From the Calibrate CGM Screen, Tap ok .

A PRECAUTION

DO NOT use alternative blood glucose site testing (blood from your palm or forearm, etc.) for calibration. Alternative site blood glucose values may be different than those taken from a fingerstick blood glucose value and may not represent the timeliest blood glucose value. Use a blood glucose value taken only from a fingerstick for calibration. Alternative site blood glucose values might affect sensor accuracy and result

in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

- 2. Tap OPTIONS.
- 3. Tap My CGM.
- Tap Calibrate CGM.
- Using the onscreen keypad, enter the blood glucose value from your meter.

A PRECAUTION

To calibrate the System, **D0** enter the exact blood glucose value that your blood glucose meter displays within 5 minutes of a carefully performed blood glucose measurement. Do not enter sensor glucose readings for calibration. Entering incorrect blood glucose values, blood glucose values obtained more than 5 minutes before entry, or sensor glucose readings might affect sensor accuracy and could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

- 6. Tap 🕶
- 7. Tap to confirm the calibration.

- Tap x if the BG value does not exactly match the reading from your meter. The onscreen keypad will appear. Enter the exact reading from your meter.
- The "CALIBRATION ACCEPTED" screen will appear, followed by the CGM Home Screen or bolus prompt.

23.5 Other Reasons You May Need to Calibrate

You may need to calibrate when your System did not accept the last calibration, or when the blood glucose value you entered for calibration is very different from the sensor glucose reading.

When you see the CALIBRATE CGM screen, calibrate using the instructions in the previous chapters.

If you see the CALIBRATION ERROR screen, you will be prompted to enter a blood glucose value to calibrate in

Chapter 23 – Calibrating Your CGM System

either 15 minutes or 1 hour, depending on the error.

Chapter 24

Viewing CGM Data on Your t:slim X2 Pump

24.1 Overview

During an active sensor session, CGM readings are sent to your t:slim X2 Pump every 5 minutes. This section teaches you how to view your sensor glucose readings and trend information. The trend graph provides additional information that your blood glucose meter does not. It shows your current glucose value, the direction it is changing and how fast it is changing. The trend graph can also show you where your glucose has been over time.

Your blood glucose meter measures glucose in your blood. Your sensor measures glucose from interstitial fluid (the fluid under your skin). Because glucose from different fluids is measured, readings from your blood glucose meter and sensor may not match.

The greatest benefit you get from using continuous glucose monitoring will come from trending information. It is important that you focus on the trends and rate of change on your receiver rather than the exact glucose reading.

▲ CONTRAINDICATION

Taking medications with acetaminophen/paracetamol while wearing the sensor may falsely raise your sensor glucose readings. The level of inaccuracy depends on the amount of acetaminophen/paracetamol active in your body and may be different for each person.

Press the Screen On button to turn the screen on. If a CGM session is active, you will see the CGM Home Screen with the 3 hour trend graph displayed.



- The current time and date are shown at the top of the screen in the middle.
- Each "dot" on the trend graph is a sensor glucose reading reported every 5 minutes.

- Your High Alert setting shows as an orange line across the trend graph.
- Your Low Alert setting shows as a red line across the trend graph.
- The gray zone highlights your target glucose range, based on your High and Low Alert settings.
- Sensor glucose readings are shown in millimoles per liter (mmol/L).
- If your sensor glucose reading is in between your High and Low Alert settings, it is shown in white.
- If your sensor glucose reading is above your High Alert setting, it is shown in orange.
- If your sensor glucose reading is below your Low Alert setting, it is shown in red.
- If the Low Alert is not set and your glucose reading is 3.1 mmol/L or lower, it is shown in red.
- The dots on the trend graph also change colors based on your High

and Low Alert settings: white if between High and Low Alert settings, orange if above High Alert setting, red if below Low Alert setting.

24.2 CGM Trend Graphs


You can view your past sensor glucose trend information on your CGM Home Screen.

1, 3, 6, 12, and 24 hour trend views can be seen. The 3 hour Trend Graph is the default view and will be shown on the Home Screen even if a different trend graph was shown when the screen turned off.

Sensor glucose information is only reported for values between 2.2 to 22.2 mmol/L. Your trend graph shows a flat line or dots at 2.2 to 22.2 mmol/L when your glucose is outside this range.

To view different Trend Graph times, tap on the Trend Graph Time (HRS) to cycle through the options.

3 Hour Trend Graph (default view) shows you your current glucose reading along with the last 3 hours of sensor glucose readings.

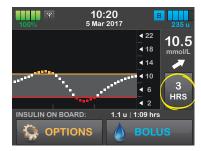

6 Hour Trend Graph shows you your current glucose reading along with the last 6 hours of sensor glucose readings.



12 Hour Trend Graph shows you your current glucose reading along with the last 12 hours of sensor glucose readings.



24 Hour Trend Graph shows you your current glucose reading along with the last 24 hours of sensor glucose readings.



Chapter 24 – Viewing CGM Data on Your t:slim X2 Pump

1 Hour Trend Graph shows you your current glucose reading along with the last 1 hour of sensor glucose readings.



LOW shows when your most recent sensor glucose reading is less than 2.2 mmol/L.



HIGH shows when your most recent sensor glucose reading is greater than 22.2 mmol/L.



24.3 Rate of Change Arrows

Your rate of change arrows add detail about the direction and speed of glucose change over the last 15–20 minutes.

The trend arrows show below your current sensor glucose reading.



Do not overreact to the rate of change arrows. Consider recent insulin dosing, activity, food intake, your overall trend graph and your blood glucose value before taking action.

If there are missed communications between the sensor and your t:slim X2 Pump during the last 15–20 minutes due to being out of range or due to an error condition, an arrow may not display. If the trend arrow is missing, and you are concerned that your blood glucose level may be rising or falling, take a blood glucose measurement using your blood glucose meter.


A WARNING

DO NOT ignore symptoms of high and low glucose. If your sensor glucose alerts and readings do not match your symptoms, measure your blood glucose with a blood glucose meter even if your sensor is not reading in the high or low range. Solely relying on the sensor glucose alerts and readings for treatment decisions could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

Chapter 24 – Viewing CGM Data on Your t:slim X2 Pump

The table below shows the different trend arrows your receiver displays:

→	Constant: Your glucose is steady (not increasing/decreasing more than 0.06 mmol/L each minute). Your glucose could increase or decrease by up to 0.9 mmol/L in 15 minutes.
*	Slowly rising: Your glucose is rising 0.06 – 0.11 mmol/L each minute. If it continued rising at this rate, your glucose could increase up to 1.7 mmol/L in 15 minutes.
1	Rising: Your glucose is rising 0.11 – 0.17 mmol/L each minute. If it continued rising at this rate, your glucose could increase up to 2.5 mmol/L in 15 minutes.
11	Rapidly rising: Your glucose is rising more than 0.17 mmol/L each minute. If it continued rising at this rate, your glucose could increase more than 2.5 mmol/L in 15 minutes.

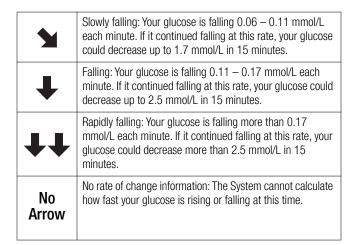


24.4 Additional CGM Status Symbols

The following symbols may also appear on your CGM Home Screen in the place where sensor glucose readings normally show. They let you know important information about the System.

You will not get sensor glucose readings during the time a status symbol shows except during the regular 12-hour calibration prompt.

Chapter 24 – Viewing CGM Data on Your t:slim X2 Pump

	Startup Calibration Needed This symbol in the place where sensor glucose readings normally show means that you need to enter startup calibrations. It remains on the screen until you calibrate with BG values.
	Additional Calibration Required This symbol in the place where sensor glucose readings normally show means that you need to enter another startup calibration. It remains on the screen until you calibrate with a BG value.
۵	Calibration Needed This symbol to the right of the antenna means that you need to enter a calibration. It shows when it is time for your 12-hour calibration update or any other time an additional calibration is needed.
Z •	Wait 15 Minute Calibration Error This symbol in the place where sensor glucose readings normally show means that the sensor cannot calibrate right now. If you see this, enter at least one more calibration blood glucose value after 15 minutes. If the sensor still cannot calibrate after that, the sensor needs to be removed and a new sensor needs to be inserted.
 mmol/L	Glucose Reading Error This symbol in the place where sensor glucose readings normally show means that the pump does not understand the sensor signal but it likely to recover. This symbol relates to the sensor only. You should wait for more prompts and do not enter any calibrations when you see this symbol.

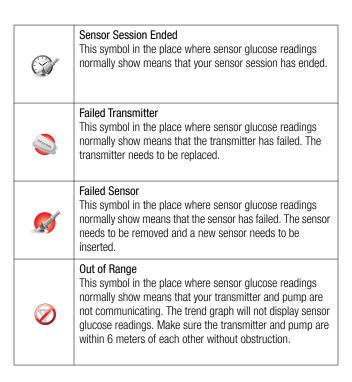

24.5 CGM History

CGM History displays the historical log of CGM events. At least 90 days of data can be viewed in History. When the maximum number of events is reached, the oldest events are removed from the history log and replaced with the most recent events. The following history sections can be viewed:

- Sessions and Calibrations
- Alerts and Errors
- Complete

Each section above is organized by date. If there are no events associated with a date, the day will not be shown in the list.

The Sessions and Calibrations section includes the start time and date for each Sensor Session, the stop time and date for each Sensor Session, and all calibration blood glucose values entered.

The Alerts and Errors section includes the date and time for all Alerts and Errors that occurred. The letter "D" (D: Alert) before an Alert or Alarm indicates the time it was declared. The letter "C" (C: Alert) indicates the time it was cleared.

The Complete history section includes all information from the Sessions and Calibrations and Alerts and Errors sections as well as any changes to Settings.

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap the Down Arrow.
- 3. Tap CGM History.
- Tap section you want to view. Each section is organized by date. Tap the date to view events from that day. Use the **Down Arrow** to scroll to more dates.

THIS PAGE IS INTENTIONALLY LEFT BLANK

Chapter 25

Ending Your CGM Sensor Session

Chapter 25 – Ending Your CGM Sensor Session

Ending Your Sensor Session

Your sensor gives you sensor glucose readings for up to 7 days. The performance of the sensor has not been tested beyond 7 days. When the sensor session ends, you will need to replace the sensor and start a new sensor session. In some cases your sensor session may end before you have finished a full 7 day period. You may also choose to end the sensor session early.

Glucose alerts and alarms do not work after the sensor session ends.

25.1 Automatic Sensor Shut-Off

Your t:slim X2 Pump tells you how much time you have left until your sensor session is complete. The "SENSOR EXPIRING SOON" screen shows at 6 hours remaining, 2 hours remaining, and 30 minutes remaining before your 7 day session ends. You will continue to receive sensor glucose readings after each reminder.

When you see the SENSOR EXPIRING SOON Screen:

- 1. Tap ok to return to the previous screen.
- The "SENSOR EXPIRING SOON" screen will show again when there are 2 hours remaining, and when there are 30 minutes remaining.
- After the final 30 minutes, the "REPLACE SENSOR" screen is displayed.
- 2. Tap ок .
- The Home Screen will appear with the Replace Sensor icon in the place where sensor glucose readings normally show.

New sensor glucose readings do not show on your t:slim X2 Pump after your sensor session ends. You must remove your sensor and insert a new sensor.

25.2 Ending a Sensor Session Before Automatic Shut-Off

You can end your sensor session at any time before the automatic sensor shutoff. To end your sensor session early:

- 1. From the Home Screen, tap OPTIONS.
- 2. Tap My CGM.
- 3. Tap STOP SENSOR.
- 4. Tap to confirm.
- The "SENSOR STOPPED" screen is temporarily displayed.
- ✓ The Home Screen will appear with the Replace Sensor icon in the place where sensor glucose readings normally show.

New sensor glucose readings do not show on your t:slim X2 Pump after your sensor session ends. You must remove your sensor and insert a new sensor.

25.3 Removing the Sensor and Transmitter

WARNING

DO NOT ignore sensor fractures. Sensors may fracture on rare occasions. If a sensor breaks and no portion of it is visible above the skin, do not attempt to remove it. Seek professional

medical help if you have symptoms of infection or inflammation—redness, swelling or pain—at the insertion site. If you experience a broken sensor, please report this to your local Tandem Diabetes Care representative.

Do not remove the transmitter from the sensor pod while the pod is attached to your skin. When you remove the sensor, make sure to pull out the sensor pod while the transmitter is still attached.

Gently peel up the sensor pod adhesive patch from your skin. This will pull out your sensor.

A PRECAUTION

DO NOT discard your transmitter. It is reusable. The same transmitter is used for each session until you have reached the end of the transmitter battery life. After the sensor pod is off your body, you must remove the transmitter to reuse it. Use either of the two transmitter removal methods below.


Method 1

The safety lock that you removed from the applicator (see Section 21.4 Placing the CGM Sensor) can be used as a tool to remove the transmitter.

- Grasp the end of the adhesive patch. Peel adhesive patch up and away from your body to remove sensor pod and transmitter.
- 2. Put sensor pod on flat surface.
- 3. Place the sensor pod on a flat surface.
- 4. Hold the rounded edge of the safety lock.
- 5. Make sure the jagged edge of the safety lock is facing down, with the arrow pointing up.

Push the safety lock down until it will go no further. The transmitter will "pop" out of the sensor pod.

Safety Lock

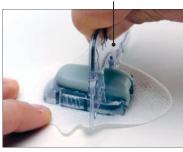

Transmitter with Safety Lock snapped in



Chapter 25 – Ending Your CGM Sensor Session

Method 2

If you did not save the safety lock when you inserted the sensor, you can use your fingers to spread out the tabs at the back of the sensor pod.

- Grasp the end of the adhesive patch. Peel adhesive patch up and away from your body to remove sensor pod and transmitter.
- 2. Put sensor pod on flat surface.
- Grasp sensor pods wide end with two fingers and place fingers in sides' open slots.
- 4. Pull tabs away from transmitter.

THIS PAGE IS INTENTIONALLY LEFT BLANK

Section 6

CGM Alerts, Errors and Troubleshooting

Chapter 26

CGM Alerts and Errors

Chapter 26 – CGM Alerts and Errors

This section describes CGM alerts and errors that appear on your t:slim X2 Home Screen. It applies only to the CGM portion of your System. CGM alerts and errors do not follow the same pattern of vibration and beeps as insulin delivery reminders, alerts, and alarms.

For information on insulin delivery reminders, alerts, and alarms, see Chapters 14 User Settable Alerts & Reminders, 15 t:slim X2 Pump Alerts, and 16 t:slim X2 Pump Alarms.

26.1 Startup Calibration Alert

What will I see on the screen?



What does it mean?

2-hour CGM startup period is complete.

How will the System notify me?

1 vibration, then vibration/beep every 5 minutes until confirmed.

Will the System re-notify me?

Yes, every 15 minutes until you calibrate.

How should I respond?

Tap and enter 2 separate blood glucose values to calibrate the system and start your CGM session.

26.2 Second Startup Calibration Alert

What will I see on the screen?



What does it mean?

The System needs an additional blood glucose value to complete startup calibration.

How will the System notify me?

1 vibration, then vibration/beep every 5 minutes until confirmed.

Will the System re-notify me?

Yes, every 15 minutes until second calibration is entered.

How should I respond?

Tap ox and enter a blood glucose value to calibrate the System and start your CGM session.

26.3 12 Hour Calibration Alert

What will I see on the screen?



What does it mean?

The System needs a blood glucose value to calibrate.

How will the System notify me?

On screen only with no vibration or beep.

Will the System re-notify me?

Yes, every 15 minutes.

How should I respond?

Tap and enter a blood glucose value to calibrate the system.

A WARNING

CALIBRATE your CGM at least once every 12 hours. Calibrating less often than every 12 hours might cause sensor glucose readings to be inaccurate and glucose alerts to become unreliable. This could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

26.4 Incomplete Calibration

What will I see on the screen?



What does it mean?

If you start to enter a calibration value using the keypad and do not complete the entry within 90 seconds, this screen appears.

How will the System notify me?

2 beeps or vibrations depending on Sound Volume selected.

Will the System re-notify me?

Yes, every 5 minutes until confirmed.

How should I respond?

Tap ox and complete your calibration by entering the value using the onscreen keypad.

26.5 Calibration Timeout

What will I see on the screen?

CALIBRATION TIMEOUT

You have exceeded the maximum time to calibrate your CGM.

Please use a new BG reading for CGM calibration.

OK

What does it mean?

If you start to enter a calibration value using the keypad and do not complete the entry within 5 minutes, this screen appears.

How will the System notify me?

2 beeps or vibrations depending on Pump Volume selected.

Will the System re-notify me?

Yes, every 5 minutes until confirmed.

How should I respond?

Tap ox and obtain a new blood glucose value using your meter. Enter the value using the onscreen keypad to calibrate the System.

26.6 Wait 15 Minute Calibration Error Alert

What will I see on the screen?



What does it mean?

The sensor cannot calibrate.

How will the System notify me?

1 vibration, then vibration/beep every 5 minutes until confirmed.

Will the System re-notify me?

No.

How should I respond?

Tap ok to confirm. Wait 15 minutes then enter 1 more blood glucose value. Wait 15 more minutes. If error screen still appears, enter 1 more blood glucose value. Wait 15 minutes. If no sensor glucose readings appear, the sensor needs to be replaced.

26.7 Calibration Required Alert

What will I see on the screen?



What does it mean?

The System needs a blood glucose value to calibrate. Sensor glucose readings will not be displayed at this time.

How will the System notify me?

1 vibration, then vibration/beep every 5 minutes until confirmed.

Will the System re-notify me?

Yes, every 15 minutes.

How should I respond?

Tap and enter a blood glucose value to calibrate the System.

26.8 CGM High Alert

What will I see on the screen?



What does it mean?

Your most recent sensor glucose reading is at or above the High Alert setting.

How will the System notify me?

2 vibrations, then 2 vibrations/beeps every 5 minutes until confirmed or your glucose value drops below the Alert level.

Will the System re-notify me?

No, unless you have turned on the Repeat feature.

How should I respond?

26.9 CGM Low Alert

What will I see on the screen?



What does it mean?

Your most recent sensor glucose reading is at or below the Low Alert setting.

How will the System notify me?

3 vibrations, then 3 vibrations/beeps every 5 minutes until confirmed or your glucose value goes above the Alert level.

Will the System re-notify me?

Only if you have turned on the Repeat feature.

How should I respond?

26.10 CGM Fixed Low Alert

What will I see on the screen?



What does it mean?

Your most recent sensor glucose reading is at or below 3.1 mmol/L.

How will the System notify me?

4 Vibrations, then 4 vibrations/beeps every 5 minutes until confirmed or your glucose value goes above 3.1 mmol/L.

Will the System re-notify me?

Yes, 30 minutes after each confirmation until your glucose value goes above 3.1 mmol/ $\rm L.$

How should I respond?

26.11 CGM Rise Alert

What will I see on the screen?



What does it mean?

Your glucose levels are rising at 0.11 mmol/L per minute or faster (at least 1.7 mmol/L in 15 minutes).

How will the System notify me?

2 vibrations, then 2 vibrations/beeps every 5 minutes or until confirmed.

Will the System re-notify me?

No.

How should I respond?

26.12 CGM Rapid Rise Alert

What will I see on the screen?



What does it mean?

Your glucose levels are rising at 0.17 mmol/L per minute or faster (at least 2.5 mmol/L in 15 minutes).

How will the System notify me?

2 vibrations, then 2 vibrations/beeps every 5 minutes or until confirmed.

Will the System re-notify me?

No.

How should I respond?

26.13 CGM Fall Alert

What will I see on the screen?



What does it mean?

Your glucose levels are falling at 0.11 mmol/L per minute or faster (at least 1.7 mmol/L in 15 minutes).

How will the System notify me?

3 vibrations, then 3 vibrations/beeps every 5 minutes or until confirmed.

Will the System re-notify me?

No.

How should I respond?

26.14 CGM Rapid Fall Alert

What will I see on the screen?



What does it mean?

Your glucose levels are falling at 0.17 mmol/L per minute or faster (at least 2.5 mmol/L in 15 minutes).

How will the System notify me?

3 vibrations, then 3 vibrations/beeps every 5 minutes or until confirmed.

Will the System re-notify me?

No.

How should I respond?

26.15 Unknown Sensor Reading

What will I see on the screen?



What does it mean?

The sensor is sending sensor glucose readings that the System does not understand. You will not receive sensor glucose readings.

How will the System notify me?

On screen only with no vibration or beep.

Will the System re-notify me?

The 3 dashes will remain on the screen until a new glucose reading is received and displayed in their place.

How should I respond?

Wait 30 minutes for more information from the system. Do not enter blood glucose values for calibration. The system will not use blood glucose values for calibration when "- - -" appears on the screen.

26.16 Out of Range Alert

What will I see on the screen?



What does it mean?

The transmitter and pump are not communicating and you will not receive sensor glucose readings.

How will the System notify me?

1 vibrate, then vibration/beep every 5 minutes until the transmitter and pump are back in range.

Will the System re-notify me?

Yes, if the transmitter and pump remain out of range.

How should I respond?

Tap ok to confirm and move the transmitter and pump closer together, or remove the obstruction between them.

26.17 Low Transmitter Battery Alert

What will I see on the screen?



What does it mean?

Transmitter battery is low.

How will the System notify me?

1 vibration, then vibration/beep every 5 minutes until confirmed.

Will the System re-notify me?

Yes, the alarm will notify you when there are 21, 14, and 7 days of transmitter battery life remaining.

How should I respond?

Tap to confirm. Replace the transmitter as soon as possible.

26.18 Transmitter Error

What will I see on the screen?



What does it mean?

The transmitter has failed and the CGM session has stopped.

How will the System notify me?

1 vibration, then vibration/beep every 5 minutes.

Will the System re-notify me?

No.

How should I respond?

Tap MORE INFO. A screen notifying you that your CGM session has stopped but insulin delivery continues will appear.

Replace the transmitter immediately.

26.19 Failed Sensor Error

What will I see on the screen?



What does it mean?

The sensor is not working properly and the CGM session has stopped.

How will the System notify me?

1 vibration, then vibration/beep every 5 minutes.

Will the System re-notify me?

No.

How should I respond?

Tap MORE INFO. A screen notifying you that your CGM session has stopped but insulin delivery continues will appear.

Replace the sensor and begin a new CGM session.

26.20 CGM System Error

What will I see on the screen?



What does it mean?

Your CGM System is not working properly; the CGM session has stopped and the system can no longer be used.

How will the System notify me?

1 vibration, then vibration/beep every 5 minutes.

Will the System re-notify me?

No.

How should I respond?

Tap MORE INFO. A screen notifying you that your CGM System cannot operate but insulin delivery continues will appear. Contact your local Tandem Diabetes Care representative.

THIS PAGE IS INTENTIONALLY LEFT BLANK THIS PAGE IS INTENTIONALLY LEFT BLANK

Chapter 27

CGM Troubleshooting

Chapter 27 - CGM Troubleshooting

This chapter provides helpful tips and instructions to help you fix issues you may have while using the CGM portion of your System.

If the troubleshooting steps in this chapter do not fix your issue, contact your local Tandem Diabetes Care representative.

27.1 CGM Pairing Troubleshooting

Possible issue:

Difficulty pairing your Dexcom G5 Mobile CGM with your t:slim X2 Insulin Pump.

Troubleshooting tip:

The Dexcom G5 Mobile CGM only allows pairing with one medical device at a time. Ensure your CGM is not connected to the Dexcom Receiver before pairing with the pump. You may still use a smartphone with the Dexcom Mobile App and your t:slim X2 Insulin Pump simultaneously with the same transmitter ID. Reference Section 19.2 Disconnecting from the Dexcom Receiver.

27.2 Calibration Troubleshooting

To ensure proper calibration of your CGM, follow these important tips.

Before you take a blood glucose value for calibration, wash your hands, make sure your glucose test strips have been stored properly and are not expired and make sure that your meter is properly coded (if required). Carefully apply the blood sample to the test strip following the instructions that came with your meter or test strips.

Do not calibrate if you see the Out of Range symbol in the place where your sensor glucose readings are normally shown on the screen.

Do not calibrate if you see "- - -" in the place where you sensor glucose readings are normally shown on the screen.

Do not calibrate if your blood glucose value is below 2.2 mmol/L or above 22.2 mmol/L.

Make sure you have not taken any medications containing acetaminophen/paracetamol.

27.3 Sensor Insertion Troubleshooting

Possible issue:

You are having trouble taking the safety lock off the sensor applicator barrel.

Troubleshooting tip:

 Make sure to pull the safety lock straight out away from your body. Use the arrows on the safety lock as a quide.

Possible issue:

You are not able to pull the collar up the applicator barrel after inserting the sensor.

Troubleshooting tips:

 Make sure the white plunger is completely pressed down before pulling the collar up. Try using more force when pulling the collar up.

Possible Issue:

You are not able to remove the applicator barrel from the sensor pod.

Troubleshooting tips:

- Make sure the collar is pulled all the way up. When pulling the collar up you should hear 2 "clicks". The collar should be as close to the top of the applicator as possible.
- Make sure the transmitter latch is flat against the adhesive on your body before squeezing the release tabs.
- Squeeze the center of the ribbed release tabs on the sides of the sensor pod and lift the applicator away from your body.

Possible Issue:

You are not able to remove the transmitter latch from the sensor pod.

Troubleshooting tips:

- Hold the sensor pod with one hand and twist the transmitter latch with the other hand to remove it.
- Do not try to snap it straight off.

Possible Issue:

The sensor pod is not sticking long enough on your body.

Troubleshooting tips:

- Make sure your skin is clean, clear of any cream or lotion, and completely dry before you insert the sensor.
- If hair is preventing the sensor pod from sticking, shave your skin before you insert the sensor.
- You may use medical tape (such as Blenderm[™], Tegaderm[™], IV 3000, 3M tape) over the white adhesive patch of the sensor pod, but do not place the tape over the transmitter

or the plastic parts of the sensor pod.



27.4 Unknown Sensor Reading Troubleshooting

When your CGM cannot provide a sensor glucose reading "- - -" shows in the place where your sensor glucose is normally shown on the screen. This means that the System does not understand the sensor signal temporarily.

Often the System can correct the problem and continue providing sensor glucose readings. If it has been at least 3 hours since your last sensor glucose reading, contact your local Tandem Diabetes Care representative.

Chapter 27 - CGM Troubleshooting

Do not enter any blood glucose values for calibration when you see "- - -" on your screen. The System will not use a blood glucose value for calibration when this symbol is on your screen.

If you see "- - -" often during a sensor session, follow the troubleshooting tips below before inserting another sensor.

- Make sure your sensor is not expired.
- Make sure your sensor pod is not dislodged or peeling up.
- Make sure your transmitter is snapped in completely.
- Make sure nothing is rubbing the sensor pod (i.e. clothing, seat belts, etc.).
- Make sure to select a good insertion site.
- Make sure your insertion site is clean and dry before sensor insertion.

 Wipe the bottom of the transmitter with a damp cloth or isopropyl alcohol wipe. Place the transmitter on a clean, dry cloth and air dry for 2–3 minutes.

27.5 Out of Range/No Antenna Troubleshooting

A PRECAUTION

AVOID separating the transmitter and receiver by more than 6 meters. The transmission range from the transmitter to the receiver is up to 6 meters without obstruction. Wireless communication does not work well through water so the range is much less if you are in a pool, bathtub, or on a water bed, etc. Types of obstruction differ and have not been tested. If your transmitter and receiver are farther than 6 meters apart or are separated by an obstruction, they might not communicate or the communication distance may be shorter and result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

If you see the Out of Range icon on your screen in the place where your sensor glucose reading normally shows, then your t:slim X2 Pump is not communicating with your transmitter

and sensor glucose readings will not show on your screen. Each time you start a new sensor session, wait 10 minutes for your t:slim X2 Pump to start communicating with your transmitter. When a sensor session is active, you may sometimes experience loss of communication for 10 minutes at a time. This is normal.

If you see the Out of Range icon for more than 10 minutes, move your t:slim X2 Pump and CGM transmitter closer together and remove any obstructions. Wait 10 minutes and communication should be restored.

You must enter your transmitter ID correctly into your pump to receive sensor glucose readings (see Section 19.1 About Bluetooth). Make sure you have removed your sensor and stopped your sensor session before checking or changing your transmitter ID. You cannot change your transmitter ID during a sensor session.

If you are still having trouble getting sensor glucose readings, contact your local Tandem Diabetes Care representative.

27.6 Failed Sensor Troubleshooting

The System may detect issues with your sensor where it cannot determine your glucose reading. The sensor session ends and the "FAILED SENSOR" screen shows on your t:slim X2 Pump. If you see this screen, it means your CGM session has ended.

- Remove your sensor and insert a new sensor.
- To help improve future sensor performance, follow the troubleshooting tips below.
- Make sure your sensor is not expired.
- Make sure your sensor pod is not dislodged or peeling up.
- Make sure your transmitter is snapped in completely.
- Make sure nothing is rubbing the sensor pod (i.e. clothing, seat belts, etc.).

 Make sure you have selected a good insertion site.

27.7 Sensor Inaccuracies

Inaccuracies are usually related to your sensor only and not to your transmitter or pump. Your sensor glucose readings are meant to be used for trending purposes only. The sensor measures glucose in the fluid under the skin—not in blood, and sensor glucose readings are not identical to readings from your blood glucose meter.

A WARNING

CALIBRATE your CGM at least once every 12 hours. Calibrating less often than every 12 hours might cause sensor glucose readings to be inaccurate and glucose alerts to become unreliable. This could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events.

A PRECAUTION

To calibrate the System, **DO** enter the exact blood glucose value that your blood glucose meter displays within 5 minutes of a carefully performed blood glucose measurement. Do not enter sensor glucose values for calibration.

Entering incorrect blood glucose values, blood glucose values obtained more than 5 minutes before entry, or sensor glucose readings might affect sensor accuracy and could result in you missing severe hypoglycemia (low blood glucose) or hyperglycemia (high blood glucose) events

If the difference between your sensor glucose reading and blood glucose value is greater than 20% of the blood glucose value for sensor readings > 4.4 mmol/L or greater than 20 points for sensor readings < 4.4 mmol/L, wash vour hands and take another blood glucose measurement. If the difference between this second blood alucose measurement and the sensor is still greater than 20% for sensor readings > 4.4 mmol/L or greater than 20 points for sensor readings < 4.4 mmol/L, recalibrate vour sensor using the second blood glucose value. The sensor glucose reading will correct over the next 15 minutes. If you see differences between your sensor glucose readings and blood glucose values outside of this acceptable range, follow the troubleshooting tips below before inserting another sensor:

Chapter 27 - CGM Troubleshooting

- Make sure your sensor is not expired.
- Make sure you do not calibrate when "- - -" or the Out of Range icon are on the screen.
- Do not use alternative blood glucose site testing (blood from your palm or forearm, etc.) for calibration as alternative site readings may be different than those from a blood glucose value. Use a blood glucose value only from your fingers for calibration.
- Use only blood glucose values between 2.2–22.2 mmol/L for calibration. If one or more of your values is outside of this range, the receiver will not calibrate.
- Use the same meter you routinely use to measure your blood glucose to calibrate. Do not switch your meter in the middle of a sensor session. Blood glucose meter and strip accuracy vary between blood glucose meter brands.

- Before taking a blood glucose measurement for calibration, wash your hands, make sure your glucose test strips have been stored properly and are not expired and make sure that your meter is properly coded (if required).
 Carefully apply the blood sample to the test strip following the instructions provided with your meter or test strips.
- Make sure you are using your blood glucose meter following the manufacturer's instructions to get accurate blood glucose values for calibration.
- Make sure you have not taken any medications containing acetaminophen/paracetamol to ensure you are getting accurate blood glucose values for calibration.

THIS PAGE IS INTENTIONALLY LEFT BLANK

Section 7

Living With and Caring for Your System

Chapter 28

Lifestyle Issues and Travel

28.1 Overview

While the convenience and flexibility of the t:slim X2 System allow most users to participate in a variety of activities, some lifestyle changes may be required. Additionally, your insulin needs may change in response to lifestyle changes.

A PRECAUTION

CONSULT your healthcare provider about lifestyle changes such as weight gain or loss, and starting or stopping exercise. Your insulin needs may change in response to lifestyle changes. Your basal rate(s) and other settings may need adjustment.

Physical Activity

The t:slim X2 System can be worn during most forms of exercise, such as running, cycling, hiking, and resistance training. During exercise, the pump can be worn in the provided case, your pocket, or other third-party "sport cases."

For activities where contact is a concern, such as hockey, martial arts, or

basketball, you can disconnect from your pump for short periods of time. If planning to disconnect from your pump, discuss a plan with your health-care provider to compensate for any basal insulin delivery you miss while disconnected, and be sure to continue to check your blood glucose levels. Even if you disconnect your tubing from your infusion site, the t:slim X2 Pump should continue to receive data from the transmitter as long as it is within the 6 meter range without obstruction.

Aquatic Activities

A PRECAUTION

AVOID submersing your pump in fluid beyond a depth of 1 meter or for more than 30 minutes (IPX7 rating). If your pump has been exposed to fluid beyond these limits, check for any signs of fluid ingress. If there are signs of fluid entry, discontinue use of the pump and contact your local Tandem Diabetes Care representative.

Your t:slim X2 Pump is watertight to a depth of 1 meter for up to 30 minutes (IPX7 rating), but it is not waterproof. Your pump should not be worn while swimming, scuba diving, surfing, or during any other activities that could

submerge the pump for an extended period of time. Your pump should not be worn in hot tubs or Jacuzzis.

The Dexcom G5 Mobile Sensor is water resistant when showering, bathing or swimming if the transmitter is fully snapped in. The sensor has been tested to be water resistant when submerged for up to 2.5 meters and up to 24 hours. Underwater use will affect the ability to communicate with the t:slim X2 pump, so the range will be much less than during normal use. Extended contact with water could weaken the adhesive used by your infusion sets and Dexcom CGM Sensors and cause them to fall off prematurely.

Extreme Altitudes

Some activities, such as hiking, skiing or snowboarding, could expose your pump to extreme altitudes. The t:slim X2 Pump has been tested at altitudes up to 3,000 meters at standard operating temperatures.

Extreme Temperatures

You should avoid activities which could expose your System to temperatures below 5°C (40°F) or above 37°C (99°F), as insulin can freeze at low temperatures or degrade at high temperatures.

Other Activities Which Require Removing Your Pump

There are other activities, such as bathing and intimacy, when it may be more convenient for you to remove your pump. It is safe to do so for short periods of time. If planning to disconnect from your pump, discuss a plan with your healthcare provider for compensating for any basal delivery you miss while disconnected, and be sure to check your blood glucose levels frequently. Missing basal delivery could cause your blood sugar to rise.

Travel

The flexibility afforded by an insulin pump can simplify some aspects of travel, but it still requires planning. Be sure to order your pump supplies before your trip so that you have

enough supplies with you while you're away from home. In addition to pump supplies, you should also always bring the following items:

- The items listed in the Emergency Kit described in 2.8 Emergency Kit.
- A prescription for both rapid-acting and long-acting insulin of the type recommended by your healthcare provider in case you need to take insulin by injection.
- A letter from your HCP explaining the medical need for your insulin pump and other supplies.

Traveling by Air

A PRECAUTION

DO NOT expose your System to X-ray screening used for carry-on and checked luggage. Newer full body scanners used in airport security screening are also a form of X-ray and your System should not be exposed to them. Notify an airport security agent that your System cannot be exposed to X-ray machines and request an alternate means of screening.

Your System has been designed to withstand common electromagnetic interference including airport metal detectors.

The System is safe for use on U.S. commercial airlines. The Dexcom G5 Mobile Transmitter is an M-PED with emission levels that meet FAA standards and may be used on board the aircraft without any further testing by the operator.

Pack your pump supplies in your carryon luggage. DO NOT pack your supplies in checked luggage as it could get delayed or lost.

If traveling internationally, contact your local Tandem Diabetes Care representative prior to your trip to obtain a travel loaner pump in case your pump malfunctions outside of your replacement area.

THIS PAGE IS INTENTIONALLY LEFT BLANK

Chapter 29

Taking Care of Your t:slim X2 System

29.1 Overview

This section provides information on caring for and maintaining your System.

Cleaning Your System

When cleaning your t:slim X2 Pump, use a damp lint-free cloth. Do not use household or industrial cleaners, solvents, bleach, scouring pads, chemicals, or sharp instruments. Never submerge the pump in water or use any other liquid to clean it. Do not place the pump in the dishwasher or use hot water to clean it. If needed, use only a very mild detergent, such as a bit of liquid soap with warm water. When drying your pump, use a soft towel; never place your pump in a microwave oven or baking oven to dry it.

Wipe the outside of the transmitter with a damp lint-free cloth or isopropyl alcohol wipe between uses.

Inspecting Your System for Damage

A PRECAUTION

DO NOT use your pump if you think it might be damaged due to dropping it or hitting it against a hard surface. Check that the System is working properly by plugging a power source into the USB port and confirming that the display turns on, you hear audible beeps, feel the pump vibrate, and see the green LED light blinking around the edge of the Screen On Button. If you are unsure about potential damage, discontinue use of the System and contact your local Tandem Diabetes Care representative.

If you drop your t:slim X2 Pump or it has been hit against something hard, ensure that it is still working properly. Check that the touch screen is working and clear, and that the cartridge and infusion set are properly in place. Check for leaks around the cartridge and at the tubing connector to the infusion set. Immediately contact your local Tandem Diabetes Care representative if you notice any cracks, chips, or other damage.

If your transmitter is damaged or cracked, do not use it. Immediately contact your local Tandem Diabetes Care representative if you notice any cracks or other damage.

Do not use the sensor if its sterile package has been damaged or opened.

Storing Your System

If you need to stop using your t:slim X2 Pump for a long period of time, you can place the pump in storage mode. To place the pump in storage mode, connect the pump to a power source and then press and hold down the Screen On/Quick Bolus Button for 25 seconds. The pump will beep 3 times before going into storage mode. Disconnect the pump from the power source.

Keep the pump protected when not in use. Store at temperatures between -20°C (-4°F) and 60°C (140°F) and at relative humidity levels between 20% and 90%.

To bring the pump out of storage mode, simply connect the pump to a power source.

Keep the sensor in its sterile packaging until you are ready to use it. Store at

temperatures between 2°C (36°F) and 25°C (77°F) and at relative humidity levels between 15% and 85%. Storing outside of this temperature range may result in reduced sensor response to glucose and may cause inaccurate CGM readings. You may store your sensors in the refrigerator if it is within the temperature range above. Sensors should not be stored in a freezer.

Do not insert sensors past the expiration date. The expiration date format is YYYY-MM-DD. Insert sensors on or before the end of the calendar day printed on the sensor package label.

Keep the transmitter protected when not in use. Store at temperatures between 0°C (32°F) and 45°C (113°F) and humidity levels between 10% and 95% relative humidity.

Disposing of System Components

Consult your healthcare provider for instructions for disposal of devices containing electronic waste such as your pump and transmitter and for instructions for disposal of potentially biohazardous materials such as used

cartridges, needles, syringes, infusion sets. and sensors.

Section 8

Technical Specifications

Chapter 30

Technical Specifications

30.1 Overview

This section provides tables of technical specifications, performance characteristics, options, settings, and electromagnetic compliance information for the t:slim X2 System. The specifications in this section meet the international standards set forth in IEC 60601-1 and IEC 60601-2-24.

30.2 t:slim X2 System Specifications

t:slim X2 System Specifications (including the t:slim X2 Pump, and Dexcom G5 Mobile Sensor and Transmitter)

Specification Type	Specification Details
Operating Conditions	Temperature: 5°C (41°F) to 37°C (98.6° F) Humidity: 20% to 90% RH non-condensing
Storage Conditions	Temperature: 2°C (36°F) to 25°C (77°F) Humidity: 20% to 85% RH non-condensing
Operating Altitude	-396 meters to 3,048 meters
Moisture Protection	IPX7: Watertight to a depth of 1 meter for up to 30 minutes
Protection Against Electrical Shock	Type BF applied part

t:slim X2 Pump Specifications

Specification Type	Specification Details
Classification	External PSU: Class II, Infusion Pump. Internally-powered equipment, Type BF applied part. The risk of ignition of flammable anesthetics and explosive gases by the pump is remote. While this risk is remote, it is not recommended to operate the t:slim X2 Pump in the presence of flammable anesthetics or explosive gases.
Size	7.95 cm x 5.08 cm x 1.52 cm (3.13" x 2.0" x 0.6") (L x W x H)
Weight (with full disposable)	3.95 ounces (112 grams)
Operating Conditions	Temperature: 5°C (41°F) to 37°C (98.6°F) Humidity: 20% to 90% RH non-condensing
Storage Conditions	Temperature: -20°C (-4°F) to 60°C (140°F) Humidity: 20% to 90% RH non-condensing
Atmospheric Pressure	-396 meters to 3048 meters
Moisture Protection	IPX7: Watertight to a depth of 1 meter for up to 30 minutes
Reservoir Volume	3.0 mL or 300 units
Cannula Fill Amount	0.1 to 1.0 units of insulin
Insulin Concentration	U-100
Alarm Type	Visual, audible, and vibratory

t:slim X2 Pump Specifications (Part 1 of 3)

Specification Type	Specification Details
Basal Delivery Accuracy at all Flow Rates (tested per IEC 60601-2-24)	$\pm 5\%$ The System is designed to vent automatically when there is a pressure difference between inside the cartridge and the surrounding air. In certain conditions, such as a gradual elevation change of 300 meters, the System may not vent immediately and delivery accuracy can vary up to 15% until 3 units have been delivered or elevation changes by more than 305 meters.
Bolus Delivery Accuracy at all Volumes (tested per IEC 60601-2-24)	±5%
Patient Protection from Air Infusion	The pump provides subcutaneous delivery into interstitial tissue and does not deliver intravenous injections. Clear tubing aids in detecting air.
Maximum Infusion Pressure Generated and Occlusion Alarm Threshold	2 bar
Frequency of Basal Delivery	5 minutes for all Basal Rates
Retention Time of Electronic Memory when Internal System Battery is Fully Discharged (including Alarm Settings and Alarm History)	Greater than 30 days
Infusion Set used for Testing	Unomedical Comfort™ Infusion Set
Internal System Battery Life	4 years minimum under normal use conditions
Typical Operating Time when System is Operating at Intermediate Rate	During normal use, the intermediate rate is 2 units/hr; battery charge can be reasonably expected to last up to 7 days (5 days if using CGM) from a fully charged state to a totally discharged state

t:slim X2 Pump Specifications (Part 2 of 3)

Specification Type	Specification Details
Handling of Over-Infusion or Under- Infusion	The method of delivery isolates the insulin chamber from the patient and the software performs frequent monitoring of system status. Multiple software monitors provide redundant protection against unsafe conditions.
	Over-infusion is mitigated by continuous self-tests, layering of redundancies and confirmations, and numerous other safeguard alarms. Users are required to review and confirm the details of all bolus deliveries, basal rates, and temp rates to ensure certainty before initiating a delivery. In addition, once bolus deliveries are confirmed, the user is given 5 seconds to cancel the delivery before it is started. An optional Auto-Off alarm triggers when the user has not interacted with the pump's user interface for a predefined period of time. Under-infusion is mitigated by occlusion detection and blood glucose monitoring as blood glucose entries are recorded. Users are prompted to treat high blood glucose conditions with a correction bolus.
Bolus Volume at Release of Occlusion (2 units per hour Basal)	Less than 3 units with VariSoft™ (110 cm) Infusion Set
Residual Insulin Remaining in the Cartridge (unusable)	Approximately 15 units
Minimum Audible Alarm Volume	45 dBA at 1 meter

t:slim X2 Pump Specifications (Part 3 of 3)

Dexcom G5 Mobile CGM Product Specifications

User is the single use operator in the home environment. Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation. Do not touch the metal connectors on the bottom of the transmitter and other open connectors on the receiver, charging cable and charger.

Dexcom G5 Mobile Sensor Specifications

Specification Type	Specification Detail
Glucose Range	2.2 – 22.2 mmol/L
Sensor Life	Up to 7 days
Calibration	Commercially available blood glucose meter
Calibration Range	2.2 – 22.2 mmol/L
Storage Condition	Temperature: 2°C – 25°C (36°F to 77°F) Humidity: 15% – 85% RH
Sterilization	Sterile by radiation

Dexcom G5 Mobile Transmitter Specifications

Specification Type	Specification Detail
Part Number	P/N 9438-06
Dimensions (including sensor pod)	Length: 38.1 millimeters Width: 22.8 millimeters Thickness: 12.7 millimeters
Weight (including sensor pod)	0.4 ounces
Power Supply	Silver oxide batteries (not replaceable)
Operational Conditions	Temperature: 10°C to 42°C (50°F to 107.6°F) Humidity: 10% to 95% RH
Storage Conditions	Temperature: 0°C to 45°C (32°F to 113°F) Humidity: 10% to 95% RH
Operating Altitude	-396 meters to 4,206 meters
Limited Warranty	3 months
Moisture Protection	IP28: temporary submersion
Protection Against Electrical Shock	Type BF applied part

USB Charging/Download Cable Specifications

Specification Type	Specification Detail
Tandem P/N	004113
Length	2 meters
Туре	USB A to USB Micro B

Power Supply/Charger, AC, Wall Mount, USB Specifications

Specification Type	Specification Detail
Tandem P/N	007866
Input	100 to 240 Volts AC, 50/60 Hz
Output Voltage	5 Volts DC
Max Output Power	5 Watts
Output Connector	USB type A

Car Adapter (not included), Specifications

Specification Type	Specification Detail
Tandem P/N	003934
Input	12 Volts DC
Output Voltage	5 Volts DC
Max Output Power	5 Watts minimum
Output Connector	USB type A

PC, USB Connector, Specifications

Specification Type	Specification Detail
Output Voltage	5 Volts DC
Output Connector	USB type A
Safety Standard Compliance	60950-1 or 60601-1 or equivalent

PC, Requirements for Charging

The t:slim X2 System is designed to be connected to a host PC for battery charging and data transfer. The following minimum characteristics are required of the host PC:

- USB 1.1 port (or later)
- PC compliant with 60950-1 or equivalent safety standard

Connecting the t:slim X2 System to a host PC that is attached to other equipment could result in previously unidentified risks to the patient, operator, or a third party. The user should identify, analyze, evaluate, and control these risks.

Subsequent changes to the host PC could introduce new risks and require additional analysis. These changes can include but are not limited to changing the configuration of the PC, connecting additional items to the PC, disconnecting items from the PC, and updating or

upgrading equipment connected to the PC.

30.3 t:slim X2 Pump Options and Settings

Option/Setting Type	Option/Setting Detail
Time	12 or 24 hour clock
Maximum Basal Rate	15 units/hr
Insulin Delivery Profiles (Basal and Bolus)	6
Basal Rate Segments	16 per delivery profile
Basal Rate Increment	0.001 at programed rates equal to or greater than 0.1 units/hr
Temp Basal Rate	15 minutes to 72 hours with 1 minute resolution with a range of 0% to 250%
Bolus Setup	Can deliver based on carb input (grams) or insulin input (units). Default is units. (The range for carbs is 1 to 999 grams; the range for insulin is 0.05 to 25 units)
Insulin-to-Carb (IC) Ratio	16 time segments per 24-hour period; Ratio: 1 unit of insulin per x grams of carbs; 1:1 to 1:300 (can be set by 0.1 below 10)
BG Correction Target Value	16 time segments. 3.9–13.9 mmol/L in 0.1 mmol/L increments
Insulin Sensitivity Factor (ISF)	16 time segments; Ratio: 1 unit of insulin reduces glucose x mmol/L; 1:0.1 to 1:33.3 (0.1 mmol/L increments)
Duration of Insulin Action	1 time segment; 2 to 8 hours in 1-minute increments (default is 5 hrs)
Bolus Increment	0.01 at volumes greater than 0.05 units
Option/Setting Type	Option/Setting Details

t:slim X2 Pump Options and Settings (Part 1 of 3)

Option/Setting Type	Option/Setting Detail
Quick Bolus Increments	When set to units: 0.5, 1, 2, 5 units (default is 0.5 units); or when set to grams/carbs: 2, 5, 10, 15 grams (default is 2 g)
Maximum Extended Bolus Time	8 hours
Maximum Bolus Size	25 units
Low Reservoir Volume Indicator	Status indicator visible on Home Screen; Low Insulin Alert is user adjustable from 10 to 40 units (default is 20 units).
Auto-Off Alarm	On or Off (default is On); user-adjustable (5 to 24 hours; default is 12 hours, which you can change when option is set to On).
History Storage	At least 90 days of data
Language	English
Feature Lock	Blocks access to insulin delivery screens and pump setting screens (default is Off).
Screen Lock	Protects from unintentional taps.
Site Reminder	Prompts user to change infusion set. Can be set for 1 to 3 days at a time selected by user (default is Off).
Missed Meal Bolus Reminder	Prompts user if a bolus has not occurred during the period of time the reminder is set for. 4 reminders available (default is Off).
After Bolus Reminder	Prompts user to test BG at a selected time period after a Bolus has been delivered. Can be set between 1 to 3 hours (default is Off).
High BG Reminder	Prompts user to retest BG after a High BG has been entered. User selects High BG value and time for reminder. (default is Off).

t:slim X2 Pump Options and Settings (Part 2 of 3)

Option/Setting Type	Option/Setting Detail
Low BG Reminder	Prompts user to retest BG after a Low BG has been entered. User selects Low BG value and time for reminder. (default is Off).

t:slim X2 Pump Options and Settings (Part 3 of 3)

30.4 t:slim X2 Pump Performance Characteristics

Rate of Delivery

Characteristic	Value
25 Unit Bolus Delivery Speed	2.97 Units/min Typical
2.5 Unit Bolus Delivery Speed	1.43 Units/min Typical
20 Unit Prime	9.88 Units/min Typical

Bolus Duration

Characteristic	Value
25 Unit Bolus Duration	8 minutes 26 seconds Typical
2.5 Unit Bolus Duration	1 minute 45 seconds Typical

Time to Occlusion Alarm*

Operating Rate	Typical	Maximum
Bolus (3 units or Greater)	1 Minute 2 Seconds	3 Minutes
Basal (2 units/hr)	1 Hour 4 Minutes	2 Hours
Basal (0.1 units/hr)	19 Hours 43 Minutes	36 Hours

^{*}The time to occlusion alarm is based on insulin volume not delivered. During an occlusion event, boluses of less than 3 units may not trigger an occlusion alarm if no basal insulin is being delivered. The bolus amount will reduce the time to occlusion depending on the basal rate.

■ NOTE: Low Basal Rates and Occlusions

At very low basal rates checking your blood glucose frequently may give you an early indication of an occlusion.

30.5 Dexcom G5 Mobile CGM Transmitter Performance Characteristics

Parameter	Performance Characteristics
TX/RX Frequencies	2.402-2.480 GHz
Bandwidth	1.02 MHz
Maximum Output Power	1.0 mW EIRP
Modulation	Gaussian Frequency-Shift Keying
Data Rate	1 Mbps
Data Communication Range	6 meters

The Dexcom G5 Mobile Continuous Glucose Monitoring System is safe for use on U.S. commercial airlines. The Dexcom G5 Mobile Transmitter is an M-PED with emission levels that meet RTCA/DO160, Section 21, Category M. Per FAA Advisory, Circular #91-21, 1B, dated 8/25/06, any M-PED that meets this standard in all modes may be used onboard the aircraft without any further testing by the operator. This device can withstand exposure to common electrostatic (ESD) and electromagnetic interference (EMI).

30.6 Dexcom G5 Mobile CGM Sensor Performance Characteristics

We recommend that you review the information in this chapter with your healthcare provider to understand how well Dexcom G5 Mobile CGM performs.

Dexcom G5 Mobile CGM uses a glucose sensor to continuously measure and monitor your glucose levels. The sensor is "calibrated" using a commercially available blood glucose meter; Once calibrated, the CGM reports glucose readings up to every 5 minutes. The CGM was evaluated in a clinical study in which CGM sensor readings were compared to blood glucose values to assess its performance and how well the CGM readings compare to a laboratory test method that measures blood glucose values. Additionally, patients performed self-monitoring blood glucose meter tests at home to assess the CGM performance in real use environment.

Although the performance characteristics of the CGM are presented in the following, there is no commonly accepted statistical approach for capturing performance of continuous glucose monitors (CGMs), such as the Dexcom G5 Mobile CGM System.

Clinical Study Overview

The CGM performance was evaluated in four separate prospective clinical studies. In all four studies, subjects were required to confirm glucose readings with their SMBG meter before making any treatment decisions. Two studies included adults, and two studies included pediatrics. In the following sections and tables, the studies will be identified as follows:

Adult Studies (18 years and older) Original Adult Study: the Receiver included software version SW10050

Software 505 Adult Study: the Receiver included software version SW10505 Pediatric Studies (2 to 17 years) Original Pediatric Study: the Receiver included software version SW10050

Software 505 Pediatric Study: the Receiver included software version SW10505

The Dexcom G5 Mobile CGM System incorporates the algorithm from software version SW10505 and has a new software number.

Overview of Adult Studies

The CGM performance for adults was evaluated in two separate prospective clinical studies:

Original Adult Study (software SW10050) and the Software 505 Adult Study (software SW10505).

Differences between the studies include the number of subjects enrolled, the number of sensors worn by each participant, the SMBG meter used, and the number of clinic days each subject participated in during the study. An overview of each study is provided here.

The Original Adult Study enrolled 72 subjects, and the Software 505 Adult Study enrolled 51 subjects. All subjects had Type 1 or Type 2 diabetes mellitus, and required insulin or oral medication to manage their diabetes. In the Original Adult Study, 83% of subjects had Type 1 diabetes, and 17% of subjects had Type 2 diabetes. In the Software 505 Adult Study, 86% of subjects had Type 1 diabetes, and 14% of subjects had Type 2 diabetes. Both studies included subjects greater than 18 years of age.

Subjects in both studies used the CGM for seven days. In the Original Adult Study, thirty-six subjects each wore 2 sensors; in the Software 505 Adult Study, all subjects wore 1 sensor only. Throughout the 7-day wear period, the sensor was calibrated with an average of 2 fingersticks per day (approximately once every 12 hours). In the Original Adult Study, subjects used the LifeScan® OneTouch® Ultra® 2 meter and in the Software 505 Adult Study, subjects used Bayer's CONTOUR® NEXT USB meter.

In the Original Adult Study, all subjects were evaluated in a controlled clinic environment on all three clinic days:

Day 1, Day 4, and Day 7 of the 7-day wear period. In the Software 505 Adult Study, subjects were evaluated in one of the three clinic days so there are fewer data samples than in the Original Adult Study. While using the CGM in the clinic, subjects had their blood glucose measured every 15 minutes with a reliable laboratory method, the Yellow Springs Instrument 2300 STAT Plus™ Glucose Analyzer.

This instrument is referred to as the "YSI." Readings from the CGM were reported every 5 minutes and paired with YSI values in order to characterize how well the CGM readings agreed with laboratory standard blood glucose results. The remainder of the study took place at home, and the CGM performance was also paired with the comparative meter results, referred to as the "SMBG."

Overview of Pediatric Studies

The CGM performance for children and adolescents was evaluated in two separate prospective clinical studies: the Original Pediatric Study (SW10050) and the Software 505 Pediatric Study (SW10505). Differences between the studies include the number of subjects enrolled, the number of CGMs worn by each participant, the SMBG meter used, the length of time subjects were evaluated in a controlled clinic environment and whether or not subjects ages 13-17 had their glucose levels intentionally manipulated during the study. An overview of each study is provided here.

The Original Pediatric Study enrolled 176 subjects, with 16% of subjects younger than 6-years old, and the Software 505 Pediatric Study enrolled 79 subjects, with 20% of subjects younger than 6-years old. All subjects had Type 1 or Type 2 diabetes mellitus and required insulin or oral medication to manage their diabetes. In the Original Pediatric Study, about 99% of subjects had Type 1 diabetes and 1% had Type 2 diabetes. In the Software 505 Pediat-

ric Study, all subjects had Type 1 diabetes. Sensors were inserted in either the abdomen or upper buttocks.

Subjects in all studies used the CGM for seven days. In the Original Pediatric Study, all subjects wore 2 sensors; in the Software 505 Pediatric Study, all subjects wore 1 sensor only. Throughout the 7-day wear period, the sensors were calibrated with an average of 2 fingersticks per day (approximately once every 12 hours), using self-monitoring blood glucose (SMBG) meter values. The Original Pediatric Study used the LifeScan® OneTouch® Verio® IQ meter; the Software 505 Pediatric Study used Bayer's CONTOUR® NEXT USB meter.

All subjects were evaluated in a controlled clinic environment on Day 1, Day 4 or Day 7 of the 7-day wear period. While using the CGM in the clinic, subjects provided at least two fingerstick measurements per hour, and subjects ages 6-17 also provided venous blood for comparison to a laboratory method, the Yellow Springs Instrument 2300 STAT Plus™ Glucose Analyzer. This instrument is referred to as the "YSI." In the Original Pediatric Study, subjects'

alucose levels were not intentionally manipulated during this study; in the Software 505 Pediatric Study, subjects ages 13-17 had their glucose levels intentionally manipulated during the clinic session. Readings from the CGM were reported every 5 minutes and paired with YSI values collected every 15 minutes in order to characterize how well the CGM readings agreed with laboratory standard blood glucose results. The remainder of the study took place at home, and the CGM performance was also paired with the comparative meter results, referred to as the "SMBG."

Table 1-A. Agreement to YSI within CGM Glucose Ranges (Adult)

CGM Glucose Range ¹ (mmol/L)	Study ²	Number of paired CGM-YSI	Percent within 15/15% YSI	Percent within 20/20% YSI	Percent within 30/30% YSI	Percent Greater than 40/40% YSI
Overall	Original	9152	71%	82%	92%	3%
Overall	Software 505	2263	86%	93%	98%	1%
2.2–3.3	Original	512	67%	78%	88%	6%
2.2–3.3	Software 505	120	89%	94%	98%	0%
3.4–4.4	Original	781	73%	85%	94%	2%
3.4–4.4	Software 505	226	91%	96%	99%	0%
4.5–10	Original	3853	67%	78%	91%	3%
4.5-10	Software 505	738	84%	92%	98%	1%
10.1–16.6	Original	2784	72%	84%	93%	4%
10.1–10.0	Software 505	798	86%	93%	98%	1%
16.7.10.4	Original	775	82%	91%	97%	2%
16.7–19.4	Software 505	229	86%	94%	98%	1%
19.5–22.2	Original	447	74%	84%	91%	5%
19.0-22.2	Software 505	152	80%	92%	97%	0%

^{1.} CGM readings are within 2.2–22.2 mmol/L, inclusive.

^{2.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Table 1-B. Agreement to YSI within CGM Glucose Ranges (Pediatric)

CGM Glucose Range ¹ (mmol/L)	Study ²	Number of paired CGM-YSI	Percent within 15/15% YSI	Percent within 20/20% YSI	Percent within 30/30% YSI	Percent Greater than 40/40% YSI
Overall	Original	2922	55%	68%	85%	7%
Overall	Software 505	2262	81%	91%	96%	2%
2.2–3.3	Original	19	63%	74%	79%	21%
2.2-3.3	Software 505	86	54%	74%	91%	3%
3.4–4.4	Original	76	61%	82%	92%	4%
3.4-4.4	Software 505	142	77%	82%	90%	3%
4.5–10	Original	1155	56%	69%	84%	6%
4.5-10	Software 505	805	78%	88%	97%	1%
10.1.16.6	Original	1380	55%	68%	85%	7%
10.1–16.6	Software 505	957	89%	96%	99%	1%
16.7–19.4	Original	206	48%	62%	80%	11%
10.7-19.4	Software 505	209	81%	91%	94%	5%
19.5–22.2	Original	86	48%	61%	79%	12%
19.5–22.2	Software 505	63	64%	81%	83%	8%

^{1.} CGM readings are within 2.2–22.2 mmol/L, inclusive.

^{2.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Agreement Relative to YSI

Agreement between the CGM and blood glucose values is characterized using paired CGM and YSI values. The CGM and YSI results were compared by pairing the YSI blood glucose value to a CGM glucose reading that occurred immediately after the YSI was collected.

The agreement of the CGM to blood glucose values was assessed by calculating the percentage of CGM readings that were within 15%, 20%, 30% and greater than 40% of the YSI values. For readings less than or equal to 4.4 mmol/L the absolute difference in mmol/L between the two alucose results was calculated. For values greater than 4.4 mmol/L the absolute percent difference (%) from the YSI values was calculated. The percentages of total readings within 0.8 mmol/L or 15%, 1.1 mmol/L or 20%, 1.7 mmol/L or 30% or greater than 2.2 mmol/L or 40% are provided in Tables 1-A and 1-B. The tables are categorized within CGM glucose ranges. When you see a CGM reading on your receiver, this table shows you how likely that reading

matches your blood glucose level (measured by YSI in the study).

For example, in the SW10505 Adult Study (Table 1-A), the total number of data pairs considered in the analysis was 2263. Of these, 93% of the CGM readings fall within \pm 1.1 mmol/L of the YSI blood glucose values 4.4 mmol/L and within \pm 20% of YSI blood glucose values > 4.4 mmol/L.

Table 2-A. Number and Percentage of YSI Values When CGM Readings are "Low" or "High" (Adult)

CGM Readings	Study ¹	CGM-YSI Pairs						Total
		n	66	84	123	142	13	155
"LOW"	Original	Cumulative Percent	42%	54%	79%	92%	8%	
LOW		n	11	16	17	18	0	18
	Software 505	Cumulative Percent	61%	89%	94%	100%	0%	

CGM Readings	Study ¹	CGM-YSI Pairs	> 18.9	> 17.8	> 15.5	> 13.3	≤ 13.3	Total
		n	189	220	238	246	2	248
"HIGH"	Original	Cumulative Percent	76%	89%	96%	99%	1%	
Illuit		n	40	43	45	45	0	45
	Software 505		89%	96%	100%	100%	0%	

^{1.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Table 2-B. Number and Percentage of YSI Values When CGM Readings are "Low" or "High" (Pediatric)

CGM Readings	Study ¹ CGM-YSI Pairs		< 3.1	< 3.3	< 3.9	< 4.4	> 4.4	Total
		n	0	0	0	0	13	13
"LOW"	Original	Cumulative Percent	0%	0%	0%	0%	100%	
LOW	Software 505	n	3	5	10	15	1	16
5		Cumulative Percent	19%	31%	63%	94%	6%	

CGM Readings	Study ¹	CGM-YSI Pairs	> 18.9	> 17.8	> 15.5	> 13.3	≤ 13.3	Total
			38	51	68	69	1	70
"HIGH"	Original	Cumulative Percent	54%	73%	97%	99%	1%	
HIGH	пічп	n	14	19	22	23	1	24
	Software 505	Cumulative Percent	58%	79%	92%	96%	4%	

^{1.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Agreement When CGM Reads "LOW" or "HIGH"

The CGM reports glucose readings between 2.2 and 22.2 mmol/L. When the CGM determines the glucose reading is below 2.2 mmol/L, it displays "LOW" on the pump home screen. When the CGM determines that the glucose level is above 22.2 mmol/L, it displays "HIGH" on the pump home screen. Because the CGM does not display glucose values below 2.2 mmol/ L or above 22.2 mmol/L, the comparisons to the actual blood alucose levels (as determined by the YSI analyzer) when CGM is classified as "LOW" or "HIGH" are included separately in Tables 2-A and 2-B. The tables include the numbers and the cumulative percentages when YSI values were less than certain glucose levels (for "LOW"), and when YSI values were greater than certain glucose levels (for "HIGH").

For example, in the Software 505 Adult Study (Table 2-A), when the CGM displayed "LOW" (18 occasions), 100% (18 out of 18) of the YSI values were less than 4.4 mmol/L, and 94% (17 out of 18) of the YSI values were less than

3.9 mmol/L. When the CGM displayed "HIGH" (45 occasions), 100% (45 out of 45) of the YSI values were greater than 13.3 mmol/L, and 100% (45 out of 45) of the YSI values were greater than 15.5 mmol/L.

Table 3-A. Concurrence of CGM Readings and YSI Values (Original Adult Study)

CGM		YSI (mmol/L) Row Percentage of Matched Pairs in Each CGM Glucose Range											
mmol/L	< 2.2	2.2- 3.3	3.4- 4.4	4.5- 6.7	6.7-8.9	8.9- 11.1	11.1- 13.9	13.9- 16.7	16.7- 19.5	19.5- 22.2	> 22.2	Paired CGM- YSI	
< 2.2	6%	48%	37%	7%	1%	0%	0%	0%	0%	0%	0%	155	
2.2–3.3	4%	49%	36%	11%	1%	0%	0%	0%	0%	0%	0%	512	
3.4-4.4	0%	22%	51%	24%	1%	0%	0%	0%	0%	0%	0%	781	
4.5–6.7	0%	2%	17%	66%	13%	1%	0%	0%	0%	0%	0%	1706	
6.7–8.9	0%	0%	1%	25%	60%	13%	2%	0%	0%	0%	0%	1492	
8.9–11.1	0%	0%	0%	2%	28%	53%	16%	2%	0%	0%	0%	1240	
11.1– 13.9	0%	0%	0%	0%	3%	21%	51%	21%	3%	1%	0%	1181	
13.9– 16.7	0%	0%	0%	0%	0%	4%	19%	49%	24%	3%	0%	1018	
16.7– 19.4	0%	0%	0%	0%	0%	0%	3%	28%	51%	16%	1%	775	
19.5– 22.2	0%	0%	0%	0%	0%	0%	3%	10%	43%	38%	7%	447	
> 22.2	0%	0%	0%	0%	0%	0%	1%	6%	21%	57%	15%	248	

Table 3-B. Concurrence of CGM Readings and YSI Values (Software 505 Adult Study)

CGM	YSI (mmol/L) Row Percentage of Matched Pairs in Each CGM Glucose Range											
mmol/L	< 2.2	2.2- 3.3	3.4- 4.4	4.5- 6.7	6.7-8.9	8.9- 11.1	11.1- 13.9	13.9- 16.7	16.7- 19.4	19.5- 22.2	> 22.2	Paired CGM- YSI
< 2.2	6%	83%	11%	0%	0%	0%	0%	0%	0%	0%	0%	18
2.2–3.3	2%	74%	22%	3%	0%	0%	0%	0%	0%	0%	0%	120
3.4-4.4	0%	19%	68%	13%	0%	0%	0%	0%	0%	0%	0%	226
4.5–6.7	0%	0%	19%	72%	8%	1%	0%	0%	0%	0%	0%	347
6.7–8.9	0%	0%	0%	17%	72%	11%	0%	0%	0%	0%	0%	246
8.9–11.1	0%	0%	0%	0%	25%	59%	16%	0%	0%	0%	0%	286
11.1– 13.9	0%	0%	0%	0%	0%	16%	70%	13%	1%	0%	0%	376
13.9– 16.7	0%	0%	0%	0%	0%	2%	16%	61%	14%	7%	0%	281
16.7– 19.4	0%	0%	0%	0%	0%	0%	2%	28%	59%	10%	1%	229
19.5– 22.2	0%	0%	0%	0%	0%	0%	0%	4%	47%	45%	5%	152
> 22.2	0%	0%	0%	0%	0%	0%	0%	0%	20%	38%	42%	45

Table 3-C. Concurrence of CGM Readings and YSI Values (Original Pediatric Study)

CGM		YSI (mmol/L) Row Percentage of Matched Pairs in Each CGM Glucose Range											
mmol/L	< 2.2	2.2- 3.3	3.4- 4.4	4.5- 6.7	6.7-8.9	8.9- 11.1	11.1- 13.9	13.9- 16.7	16.7- 19.4	19.5- 22.2	> 22.2	Paired CGM- YSI	
< 2.2	0%	0%	0%	54%	31%	15%	0%	0%	0%	0%	0%	13	
2.2–3.3	0%	21%	58%	16%	5%	0%	0%	0%	0%	0%	0%	19	
3.4-4.4	0%	21%	45%	30%	4%	0%	0%	0%	0%	0%	0%	76	
4.5–6.7	0%	1%	20%	66%	12%	1%	0%	0%	0%	0%	0%	338	
6.7–8.9	0%	0%	1%	36%	54%	7%	1%	0%	0%	0%	0%	511	
8.9–11.1	0%	0%	0%	4%	40%	48%	6%	1%	0%	0%	0%	596	
11.1– 13.9	0%	0%	0%	1%	9%	44%	41%	5%	0%	0%	0%	658	
13.9– 16.7	0%	0%	0%	0%	2%	7%	50%	36%	3%	0%	2%	432	
16.7– 19.4	0%	0%	0%	0%	0%	2%	18%	59%	21%	0%	0%	206	
19.5– 22.2	0%	0%	0%	0%	0%	0%	3%	28%	50%	16%	2%	86	
> 22.2	0%	0%	0%	0%	0%	0%	1%	14%	41%	36%	7%	70	

Table 3-D. Concurrence of CGM Readings and YSI Values (Software 505 Pediatric Study)

CGM		YSI (mmol/L) Row Percentage of Matched Pairs in Each CGM Glucose Range											
mmol/L	< 2.2	2.2- 3.3	3.4- 4.4	4.5- 6.7	6.7-8.9	8.9- 11.1	11.1- 13.9	13.9- 16.7	16.7- 19.4	19.5- 22.2	> 22.2	Paired CGM- YSI	
< 2.2	6%	25%	63%	6%	0%	0%	0%	0%	0%	0%	0%	16	
2.2–3.3	0%	33%	60%	6%	1%	0%	0%	0%	0%	0%	0%	86	
3.4-4.4	0%	8%	64%	26%	2%	0%	0%	0%	0%	0%	0%	142	
4.5–6.7	0%	1%	15%	69%	13%	1%	1%	0%	0%	0%	0%	314	
6.7-8.9	0%	0%	0%	15%	66%	18%	1%	0%	0%	0%	0%	313	
8.9–11.1	0%	0%	0%	1%	18%	66%	15%	0%	0%	0%	0%	355	
11.1– 13.9	0%	0%	0%	0%	1%	17%	68%	14%	0%	0%	0%	444	
13.9– 16.7	0%	0%	0%	0%	0%	0%	26%	58%	16%	0%	0%	336	
16.7– 19.4	0%	0%	0%	0%	0%	0%	4%	40%	46%	9%	0%	209	
19.5– 22.2	0%	0%	0%	0%	0%	0%	3%	14%	62%	21%	0%	63	
> 22.2	0%	0%	0%	0%	0%	0%	4%	13%	29%	38%	17%	24	

Concurrence of CGM and Laboratory Reference

Table 3-A (Original Adult Study), 3-B (Software 505 Adult Study), 3-C (Original Pediatric Study) and 3-D (Software 505 Pediatric Study) are categorized by ranges of CGM glucose readings. These tables describe, for each range of CGM glucose readings, what percentage of paired YSI values were in the same glucose range (shaded) or in glucose ranges above and below the paired CGM readings.

For example, based on the Software 505 Adult Study, when CGM readings are within 4.5 to 6.7 mmol/L, you can expect your blood glucose levels are within 4.5 to 6.7 mmol/L 72% of time.

Table 4-A. Difference to YSI within CGM Glucose Ranges (Adult)

CGM Glucose Range ¹ (mmol/L)	Study ²	Number of paired CGM-YSI	Mean Percent Difference	Median Percent Difference	Mean Absolute Percent Difference	Median Absolute Percent Difference
Overall	Original	9152	2.9%	1.7%	13.3%	9.8%
Overall	Software 505	2263	2.5%	2.4%	9.0%	7.0%
*2.2–3.3	Original	512	-10.0	-8.2	13.5	9.7
2.2-3.3	Software 505	120	-3.3	-2.1	6.9	4.8
*3.4–4.4	Original	781	-2.4	-0.4	11.4	8.6
3.4-4.4	Software 505	226	0.8	1.4	6.7	5.4
4.5–10.0	Original	3853	4.8%	3.0%	13.8%	9.8%
4.5-10.0	Software 505	738	3.9%	4.1%	9.6%	8.2%
10.1–16.6	Original	2784	2.1%	0.0%	11.9%	9.2%
10.1-10.0	Software 505	798	0.6%	0.4%	8.0%	6.1%
16.7–19.4	Original	775	3.8%	2.8%	9.8%	7.9%
10.7-19.4	Software 505	229	4.1%	3.4%	8.0%	5.8%
10 5 22 2	Original	447	10.4%	7.7%	12.8%	9.1%
19.5–22.2	Software 505	152	7.2%	6.3%	9.2%	7.2%

^{1.} CGM readings are within 2.2–22.2 mmol/L, inclusive. 2. Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

^{*} For CGM 4.4 mmol/L, the difference and absolute difference in mmol/L are included instead of percent differences (%).

Table 4-B. Difference to YSI within CGM Glucose Ranges (Pediatric)

CGM Glucose Range ¹ (mmol/L)	Study ²	Number of paired CGM-YSI	Mean Percent Difference	Median Percent Difference	Mean Absolute Percent Difference	Median Absolute Percent Difference
Overall	Original	2922	13.5%	11.6%	17.4%	13.5%
Overall	Software 505	2262	1.8%	1.2%	10.4%	7.9%
*2.2–3.3	Original	19	-18.1	-9.1	19.2	9.1
2.2-3.3	Software 505	86	-15.3	-13.2	16.1	13.2
*3.4–4.4	Original	76	-3.7	-2.3	13.4	10.6
3.4-4.4	Software 505	142	-4.8	-1.0	11.8	7.7
4.5–10.0	Original	1155	11.9%	9.7%	17.0%	13.0%
4.5-10.0	Software 505	805	1.9%	0.7%	10.6%	8.1%
10.1.10.0	Original	1380	14.8%	12.4%	17.4%	13.3%
10.1–16.6	Software 505	957	2.2%	1.0%	8.1%	6.5%
10.7.10.4	Original	206	19.2%	15.9%	19.4%	15.9%
16.7–19.4	Software 505	209	7.8%	6.5%	11.0%	7.9%
10 5 22 2	Original	86	18.5%	15.5%	19.1%	15.5%
19.5–22.2	Software 505	63	14.9%	11.6%	15.2%	11.6%

^{1.} CGM readings are within 2.2–22.2 mmol/L, inclusive. 2. Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

^{*} For CGM ≤ 4.4 mmol/L, the difference and absolute difference in mmol/L are included instead of percent differences (%).

Accuracy Relative to YSI

Accuracy between matched pairs was also estimated by calculating the percent difference between the CGM reading and the YSI value. For example, if the YSI value is 5.5 mmol/L and the CGM reading is 5.0 mmol/L, a 10% difference between the CGM and the YSI is reported. The CGM and YSI values were compared by pairing the CGM reading that fell immediately after the YSI value was collected.

In the example above, the CGM reading is less than the YSI value, so the percent difference reading is negative. The mean percent difference is the average of all positive and negative percent differences between the two devices; it tells you if the CGM reads higher or lower on average than the YSI within each glucose range.

Another estimate used to show the accuracy of the CGM is the absolute percent difference. The absolute percent difference tells you the percent difference or "distance" between the CGM and YSI values, but does not tell you whether the CGM is reading, on aver-

age, higher or lower than the YSI laboratory standard. The mean absolute percent difference is the average "distance" (regardless if positive or negative) between CGM readings and YSI values.

Accuracy measures in differences for both the Original Adult and Software 505 Adult Studies are summarized in Table 4-A. Accuracy measures in differences for both the Original Pediatric and Software 505 Pediatric Studies are summarized in Table 4-B. Tables 4-A and 4-B are categorized within CGM glucose ranges.

For example, in the Software 505 Adult Study (Table 4-A), overall, on average, the CGM reads 2.5% different (Mean Percent Difference) than the reference and 9.0% absolute different (Mean Absolute Difference) than the reference values. The Median Percent Difference shows that half of the time the CGM reads 2.4% or less than the YSI blood glucose values and the Median Absolute Percent Difference shows that half of the time the CGM reads about 7.0% or less than the YSI blood glucose values.

Table 5-A. Hypoglycemic Alert and Detection Rate Evaluation in Reference to YSI 15 Minutes Before and After (Adult)

Hypoglycemic Alert Level (mmol/L)	Study ¹	True Alert Rate	False Alert Rate	Hypoglycemia Detection Rate	Hypoglycemia Missed Detection Rate
3.0	Original	50%	50%	71%	29%
3.0	Software 505	71%	29%	68%	32%
3.3	Original	64%	36%	75%	25%
3.3	Software 505	85%	15%	83%	17%
3.9	Original	79%	21%	83%	17%
3.9	Software 505	92%	8%	91%	9%
4.4	Original	87%	13%	86%	14%
4.4	Software 505	95%	5%	90%	10%
5.0	Original	90%	10%	89%	11%
5.0	Software 505	96%	4%	94%	6%

^{1.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Table 5-B. Hypoglycemic Alert and Detection Rate Evaluation in Reference to YSI 15 Minutes Before and After (Pediatric, Ages 6-17 Years)

Hypoglycemic Alert Level (mmol/L)	Study ¹	True Alert Rate	False Alert Rate	Hypoglycemia Detection Rate	Hypoglycemia Missed Detection Rate
3.0	Original	0%	100%	0%	100%
3.0	Software 505	22%	78%	75%	25%
3.3	Original	11%	89%	25%	75%
3.3	Software 505	42%	58%	78%	23%
3.9	Original	47%	53%	50%	50%
3.9	Software 505	68%	32%	75%	25%
4.4	Original	55%	45%	55%	45%
4.4	Software 505	86%	14%	91%	9%
5.0	Original	69%	31%	62%	38%
3.0	Software 505	90%	10%	93%	7%
5.5	Original	75%	25%	62%	38%
0.5	Software 505	91%	9%	93%	7%

^{1.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Table 5-C. Hypoglycemic Alert and Detection Rate Evaluation in Reference to SMBG 30 Minutes Before and After (Pediatric, Ages 2-5 Years)

Hypoglycemic Alert Level (mmol/L)	Study ¹	True Alert Rate	False Alert Rate	Hypoglycemia Detection Rate	Hypoglycemia Missed Detection Rate
3.0	Original	3%	97%	57%	43%
3.0	Software 505	25%	75%	100%	0%
3.3	Original	11%	89%	62%	38%
3.3	Software 505	20%	80%	100%	0%
3.9	Original	29%	71%	77%	23%
3.9	Software 505	20%	80%	100%	0%
4.4	Original	35%	65%	85%	15%
4.4	Software 505	61%	39%	100%	0%
F.0	Original	51%	49%	89%	11%
5.0	Software 505	78%	22%	100%	0%
5.5	Original	64%	36%	91%	9%
0.5	Software 505	82%	18%	100%	0%

^{1.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Low and High Glucose Alerts

The ability of the CGM to detect high and low glucose levels is assessed by comparing CGM results to YSI results at low and high blood glucose levels and determining if the alert may have sounded. The CGM and YSI values were compared by pairing the CGM reading that occurred immediately after the YSI value was collected. We suggest that you ask your doctor what alert settings would be best for you.

The Low Glucose Alert

Estimates of how well the adjustable Low Glucose Alert performs are presented in Tables 5-A, 5-B and 5-C. Table 5-A represents the hypoglycemic alert evaluation within 15 minutes of the YSI value in the adult studies. Table 5-B represents the alert evaluation within 15 minutes of the YSI value for a subset of the pediatric population—subjects age 6 to 17 years who had YSI measurements every 15 minutes. Table 5-C represents the alert evaluation within 30 minutes of an SMBG reading for 2- to 5-year old subjects in the pediatric studies.

Hypoglycemia Alert Rate

The Alert Rate shows how often the alert is right or wrong. The True Alert Rate is the % of time the device alarmed when the blood glucose level was at or below the alert setting within 15 or 30 minutes before or after the device alarmed. The False Alert Rate is the % of time the device alarmed when the blood glucose level was above the alert setting within 15 or 30 minutes before or after the device alarmed.

For example, if you set the Low Glucose Alert to 3.9 mmol/L and your alarm sounds, how often can you expect your blood sugar to actually be low? In the Software 505 Adult Study (Table 5-A), when your alarm sounds, you can expect your blood sugar to be below 3.9 mmol/L approximately 92% of the time and above 3.9 mmol/L approximately 8% of the time within the 15 minute period before or after your alarm sounds.

Hypoglycemia Detection Rate

The Detection Rate shows how often the device recognizes and alerts you to an episode of hypoglycemia or how often it misses such an event. The Hypoglycemia Detection Rate is the % of time the blood glucose level was at or below the alert setting and device alarmed within 15 or 30 minutes before or after the blood glucose was at or below the alert settings. The Hypoglycemia Missed Detection Rate is the % of time the blood glucose was at or below the alert setting, but the device did not alarm within 15 or 30 minutes before or after the blood glucose was at or below the alert setting.

For example, if you set the Low Glucose alert to 3.9 mmol/L, how often will your alarm alert you if your blood glucose goes below 3.9 mmol/L? In the Software 505 Adult Study (Table 5-A), when your blood sugar goes below 3.9 mmol/L, you can expect your alarm to sound 91% of the time and not to sound approximately 9% of time within the 15 minute period before or after your blood sugar goes below 3.9 mmol/L.

Table 6-A. Hyperglycemic Alert and Detection Rate Evaluation in Reference to YSI 15 Minutes Before and After (Adult)

Hyperglycemic Alert Level (mmol/L)	Study ¹	True Alert Rate	False Alert Rate	Hyperglycemia Detection Rate	Hyperglycemia Missed Detection Rate
6.7	Original	95%	5%	98%	2%
0.7	Software 505	98%	2%	100%	0%
7.8	Original	94%	6%	97%	3%
7.0	Software 505	97%	3%	99%	1%
10.0	Original	92%	8%	97%	3%
10.0	Software 505	97%	3%	99%	1%
11.1	Original	92%	8%	97%	3%
11.1	Software 505	96%	4%	98%	2%
12.2	Original	91%	9%	95%	5%
12.2	Software 505	94%	6%	98%	2%
10.0	Original	91%	9%	94%	6%
13.3	Software 505	93%	7%	95%	5%
10.7	Original	82%	18%	86%	14%
16.7	Software 505	86%	14%	90%	10%

^{1.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Table 6-B. Hyperglycemic Alert and Detection Rate Evaluation in Reference to YSI 15 Minutes Before and After (Pediatric, Ages 6-17 Years)

Hyperglycemic Alert Level (mmol/L)	Study ¹	True Alert Rate	False Alert Rate	Hyperglycemia Detection Rate	Hyperglycemia Missed Detection Rate
6.7	Original	91%	9%	98%	2%
0.7	Software 505	98%	2%	99%	1%
7.8	Original	87%	13%	99%	1%
7.0	Software 505	97%	3%	98%	2%
10.0	Original	75%	25%	99%	1%
10.0	Software 505	94%	6%	98%	2%
11.1	Original	71%	29%	98%	2%
11.1	Software 505	94%	6%	97%	3%
12.2	Original	67%	33%	97%	3%
12.2	Software 505	93%	7%	96%	4%
13.3	Original	62%	38%	96%	4%
13.3	Software 505	88%	12%	94%	6%
16.7	Original	43%	57%	93%	7%
16.7	Software 505	69%	31%	84%	16%

^{1.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Table 6-C. Hyperglycemic Alert and Detection Rate Evaluation in Reference to SMBG 30 Minutes Before and After (Pediatric, Ages 2-5 Years)

Hyperglycemic Alert Level (mmol/L)	Study ¹	True Alert Rate	False Alert Rate	Hyperglycemia Detection Rate	Hyperglycemia Missed Detection Rate
6.7	Original	92%	8%	98%	2%
0.7	Software 505	97%	3%	99%	1%
7.8	Original	90%	10%	98%	2%
7.0	Software 505	98%	2%	100%	0%
10.0	Original	87%	13%	96%	4%
10.0	Software 505	99%	1%	93%	7%
11.1	Original	85%	15%	96%	4%
11.1	Software 505	98%	2%	93%	7%
12.2	Original	81%	19%	95%	5%
12.2	Software 505	100%	0%	97%	3%
13.3	Original	80%	20%	95%	5%
13.3	Software 505	99%	1%	98%	2%
16.7	Original	71%	29%	90%	10%
16.7	Software 505	95%	5%	96%	4%

^{1.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

The High Glucose Alert

Estimates of how well the adjustable High Glucose Alert performs are presented in Tables 6-A, 6-B and 6-C. Table 6-A represents the hyperglycemic alert evaluation within 15 minutes of the YSI value in the adult studies. Table 6-B represents the alert evaluation within 15 minutes of the YSI value for a subset of the pediatric population—subjects age 6 to 17 years who had YSI measurements every 15 minutes. Table 6-C represents the alert evaluation within 30 minutes of an SMBG reading for 2- to 5-year old subjects in the pediatric studies.

Hyperglycemia Alert Rate

The Alert Rate shows how often the alert is right or wrong. The True Alert Rate is the % of time the device alarmed when the blood glucose level was at or above the alert setting within 15 or 30 minutes before or after the device alarmed. The False Alert Rate is the % of time the device alarmed when the blood glucose level was below the alert setting within 15 or 30 minutes before or after the device alarmed.

For example, if you set the High Glucose alert to 11.1 mmol/L and your alarm sounds, how often can you expect your blood sugar to actually be high? In the Software 505 Adult Study (Table 6-A), when your alarm sounds, you can expect your blood sugar to be at or above 11.1 mmol/L approximately 96% of the time and not be above 11.1 mmol/L approximately 4% of the time within the 15 minute period before or after your alarm sounds.

Hyperglycemia Detection Rate

The Detection Rate shows how often the device recognizes and alerts you to an episode of hyperglycemia or how often it misses such an event. The Hyperglycemia Detection Rate is the % of time the blood glucose level was at or above the alert setting and the device alarmed within 15 or 30 minutes before or after the blood glucose was at or above the alert settings. The Hyperglycemia Missed Detection Rate is the % of time the blood glucose was at or above the alert setting, but the device did not alarm within 15 or 30 minutes before or after the blood glucose was at or above the alert setting.

For example, if you set your High Glucose alert to 11.1 mmol/L, how often will your alarm alert you if your blood glucose goes at or above 11.1 mmol/L? In the Software 505 Adult Study (Table 6-A), when your blood sugar goes above 11.1 mmol/L, you can expect your alarm to sound 98% of the time and not to sound approximately 2% of time within the 15 minute period before or after your blood sugar goes above 11.1 mmol/L.

Table 7-A. Percentage of CGM Readings¹ Within YSI Values With Data Stratified in 2-Hour Increments After Calibration (Adult)

Time From Calibration	Study ²	Number of paired CGM-YSI	Percent within 15/15% YSI	Percent within 20/20% YSI	Percent within 30/30% YSI	Percent Greater than 40/40% YSI
0-2 hours	Original	1929	78%	88%	96%	2%
0-2 110015	Software 505	469	93%	97%	99%	0%
2-4 hours	Original	1516	69%	81%	91%	4%
2-4 110015	Software 505	389	90%	97%	99%	0%
4-6 hours	Original	1547	69%	79%	91%	5%
4-0 110015	Software 505	383	85%	91%	97%	2%
6-8 hours	Original	1520	68%	79%	92%	3%
0-0 110018	Software 505	380	79%	90%	97%	2%
8-10 hours	Original	1555	71%	82%	92%	4%
0-10 110015	Software 505	347	83%	92%	98%	0%
10.10 hours	Original	1068	65%	77%	91%	4%
10-12 hours	Software 505	295	80%	90%	98%	0%
12-14 hours	Original	17	65%	76%	82%	12%
12-14 110018	Software 505	0				

^{1.} CGM readings are within 2.2–22.2 mmol/L, inclusive.

^{2.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Table 7-B. Percentage of CGM Readings¹ Within YSI Values With Data Stratified in 2-Hour Increments After Calibration (Pediatric)

Time From Calibration	Study ²	Number of paired CGM-YSI	Percent within 15/15% YSI	Percent within 20/20% YSI	Percent within 30/30% YSI	Percent Greater than 40/40% YSI
0-2 hours	Original	648	65%	75%	87%	7%
0-2 110015	Software 505	545	83%	91%	97%	1%
2-4 hours	Original	649	51%	67%	86%	7%
Z-4 110u15	Software 505	460	72%	89%	96%	2%
4-6 hours	Original	630	51%	61%	80%	10%
4-6 Hours	Software 505	428	77%	88%	95%	2%
6 0 hours	Original	409	52%	68%	85%	5%
6-8 hours	Software 505	325	88%	92%	94%	3%
8-10 hours	Original	296	53%	69%	84%	7%
0-10 110018	Software 505	305	86%	93%	97%	1%
10.10 hours	Original	253	58%	74%	89%	5%
10-12 hours	Software 505	198	89%	94%	98%	0%
10.14 hours	Original	37	32%	38%	65%	22%
12-14 hours	Software 505	1	100%	100%	100%	0%

^{1.} CGM readings are within 2.2–22.2 mmol/L, inclusive.

^{2.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Calibration Stability

The CGM must be calibrated every 12 hours. To demonstrate performance of the CGM over a 12-hour calibration period, sensors were evaluated to verify that performance remains consistent over the 12-hour calibration period. CGMs were evaluated in 2-hour increments after calibration. Performance was estimated at each 2-hour interval and stratified by glucose values by calculating the percentage of CGM readings within 0.8 mmol/L or 15%, 1.1 mmol/L or 20%, 1.7 mmol/L or 30%, 2.2 mmol/L or 40% and greater than 2.2 mmol/L or 40% of the YSI values in Tables 7-A and 7-B.

Table 8-A. Sensor Stability Relative to YSI (Accuracy Over Time¹) - (Adult)

Day of Wear	Study ²	Number of paired CGM-YSI	Mean Absolute Percent Difference	Median Absolute Percent Difference	Percent within 15/15% YSI	Percent within 20/20% YSI	Percent within 30/30% YSI	Percent Greater than 40/40% YSI
Day 1	Original	3023	16.7%	13.2%	59%	71%	86%	6%
Day 1	Software 505	680	10.7%	7.9%	77%	84%	96%	2%
Day 4	Original	3108	11.4%	8.2%	77%	87%	95%	2%
Day 4	Software 505	777	8.0%	6.4%	89%	96%	99%	0%
Day 7	Original	3021	11.9%	8.9%	76%	87%	95%	2%
Day 7	Software 505	806	8.5%	7.2%	90%	97%	99%	0%

^{1.} CGM readings are within 2.2–22.2 mmol/L, inclusive.

^{2.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Table 8-B. Sensor Stability Relative to YSI (Accuracy Over Time¹) - (Pediatric, Ages 6-17 Years)

Day of Wear	Study ²	Number of paired CGM-YSI	Mean Absolute Percent Difference	Median Absolute Percent Difference	Percent within 15/15% YSI	Percent within 20/20% YSI	Percent within 30/30% YSI	Percent Greater than 40/40% YSI
Day 1	Original	1016	21.2%	15.8%	48%	61%	78%	15%
Day I	Software 505	740	12.7%	8.5%	75%	83%	91%	4%
Day 4	Original	810	16.0%	13.9%	52%	66%	87%	3%
Day 4	Software 505	795	8.1%	6.7%	89%	97%	100%	0%
Day 7	Original	1096	15.1%	11.3%	63%	76%	89%	4%
	Software 505	727	10.4%	8.4%	80%	91%	98%	1%

^{1.} CGM readings are within 2.2–22.2 mmol/L, inclusive.

^{2.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Table 8-C. Sensor Stability Relative to SMBG (Accuracy Over Time¹) - (Pediatric, Ages 2-17 Years)

Day of Wear	Study ²	Number of paired CGM-SMBG	Mean Absolute Percent Difference	Median Absolute Percent Difference	Percent within 15/15% SMBG	Percent within 20/20% SMBG	Percent within 30/30% SMBG	Percent Greater than 40/40% SMBG
Day 1	Original	3216	18.8%	14.2%	53%	65%	81%	10%
Day I	Software 505	893	14.8%	10.7%	64%	79%	91%	5%
Day 2	Original	2148	16.2%	12.4%	60%	74%	87%	6%
Day 2	Software 505	436	13.2%	10.4%	69%	81%	95%	3%
Day 2	Original	1977	15.2%	11.0%	63%	76%	89%	5%
Day 3	Software 505	441	13.8%	11.3%	66%	77%	91%	2%
Day 4	Original	2830	14.0%	10.9%	66%	79%	91%	4%
Day 4	Software 505	850	10.7%	8.5%	79%	91%	97%	1%
Doy 5	Original	1768	15.4%	10.7%	67%	78%	90%	5%
Day 5	Software 505	374	11.4%	8.7%	74%	86%	96%	1%
Doy 6	Original	1704	14.3%	9.8%	68%	79%	90%	4%
Day 6	Software 505	410	12.3%	9.2%	72%	80%	93%	2%
D. 7	Original	2675	12.4%	9.2%	72%	83%	94%	3%
Day 7	Software 505	860	11.3%	8.6%	79%	90%	96%	2%

^{1.} CGM readings are within 2.2–22.2 mmol/L, inclusive.

^{2.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Sensor Stability

Relative to YSI

Sensors can be worn for up to 7 days. Performance was estimated by calculating the percentage of CGM readings within 0.8 mmol/L or 15%, 1.1 mmol/L or 20%, 1.7 mmol/L or 30%, 2.2 mmol/L or 40% and greater than 2.2 mmol/L or 40% of the YSI values at the beginning (Day 1), middle (Day 4) and end (Day 7) of the CGM lifecycle. The average and median of the absolute percent differences are included in Tables 8-A and 8-B showing consistent accuracy and sensor stability over the 7-day life of the sensor.

Relative to SMBG (Pediatric Study)
Performance was also estimated by
calculating the percentage of CGM
readings within various percentages of
the SMBG values at each day of the
sensor wear period (Table 8-C). The
average and median of the absolute
percent differences are included in the
table.

Precision of Sensor Readings

A subset of subjects wore two sensors at the same time. This was to look at how similarly two sensors function on the same subject (sensor precision). Precision was evaluated by comparing the glucose readings from the two sensors worn on the same subject at the same time.

In the Original Adult Study, 36 subjects wore two sensors. Results showed that sensor readings from the two sensors generally agreed with each other within 9% (absolute percent difference) with a 7% coefficient of variation. In the Original Pediatric Study, all subjects wore two sensors. Results showed that sensor readings from the two sensors generally agreed with each other within 10% (absolute percent difference) with a 7% coefficient of variation. Only one sensor was worn in the Software 505 Adult and Software 505 Pediatric Studies so precision data was not collected.

Sensor Life

Sensors may be worn for up to 7 days (168 hours). To estimate how long a sensor will work over 7 days, all sensors worn were evaluated to determine how many days/hours of readings each sensor provided.

In the Original Adult Study, 108 sensors were evaluated. Ninety-four percent (94%) of the sensors lasted until Day 7 (145-168 hours). There were 6 (6%) sensors that ended early, four of which lasted more than 3 days.

In the Software 505 Adult Study, 51 sensors were evaluated. Ninety-eight percent (98%) of the sensors lasted until Day 7 (145-168 hours). There was 1 (2%) sensor that ended early, which lasted until day 5 of the sensor wear.

In the Original Pediatric Study, 351 sensors were evaluated. Eighty-five percent (85%) of the sensors lasted until Day 7 (145-168 hours).

In the Software 505 Pediatric Study, 77 sensors were evaluated. Ninety-four

percent (94%) of the sensors lasted until Day 7 (145-168 hours).

Table 9-A. Number of Readings Provided by Each Sensor Over 7-Days (Adult)

% of Total Possible Readings Provided	Study ¹	Total Readings Provided (Min-Max)	% of Systems Providing That Number of Readings
0–25%	Original	167–491	2%
0-2370	Software 505	0	0%
26–50%	Original	719–914	4%
20-50%	Software 505	856–856	2%
51–75%	Original	1267–1267	1%
31-73%	Software 505	1253–1253	2%
76–100%	Original	1811–1992	94%
70-100%	Software 505	1497–1992	96%

^{1.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Table 9-B. Number of Readings Provided by Each Sensor Over 7-Days (Pediatric)

% of Total Possible Readings Provided	Study ¹	Total Readings Provided (Min-Max)	% of Systems Providing That Number of Readings	
0–25%	Original	103–427	3%	
0-23 /0	Software 505	60–223	4%	
26–50%	Original	569–954	3%	
20-30%	Software 505	877–891	3%	
51–75%	Original	1006–1484	9%	
	Software 505	1131–1342	3%	
76–100%	Original	1518–1992	86%	
	Software 505	1623–1990	91%	

^{1.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Table 10-A. CGM Readings Within Wear Days (Adult)

Statistic	Study ¹	Day 1	Day 2	Day 3	Day 4	Day 5	Day 6	Day 7	All Days ²
Mean	Original	98%	98%	98%	98%	97%	99%	95%	97%
IVIEdIT	Software 505	98%	99%	98%	98%	96%	99%	97%	98%
Median	Original	100%	100%	100%	100%	100%	100%	100%	100%
IVICUIAII	Software 505	99%	100%	100%	100%	100%	100%	100%	100%
Standard	Original	5%	3%	9%	8%	10%	3%	11%	8%
Deviation	Software 505	3%	2%	8%	11%	15%	2%	13%	9%

^{1.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

^{2.} A total of 108 sensors were included with the Original Study and 51 sensors were included with the Software 505 Study.

Table 10-B. CGM Readings Within Wear Days (Pediatric)

Statistic	Study ¹	Day 1	Day 2	Day 3	Day 4	Day 5	Day 6	Day 7	All Days ²
Mean	Original	97%	96%	96%	95%	94%	94%	92%	97%
IVICALI	Software 505	96%	96%	95%	96%	93%	95%	93%	96%
Median	Original	99%	99%	99%	99%	99%	99%	98%	99%
IVIEUIAII	Software 505	99%	98%	99%	99%	97%	97%	98%	99%
Standard	Original	6%	10%	9%	12%	14%	14%	17%	6%
Deviation	Software 505	9%	6%	12%	10%	15%	7%	12%	9%

^{1.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

^{2.} A total of 108 sensors were included with the Original Study and 77 sensors were included with the Software 505 Study.

Number of Readings Provided

The CGM is capable of providing a reading up to every 5 minutes, or up to 288 readings per day. For a variety of reasons, the CGM may not display a glucose reading and readings are "skipped." Tables 9-A and 9-B estimate the number of readings you can expect to receive from the CGM over the entire 7-day period after calibration. Tables 10-A and 10-B show the number of readings you can expect to receive from the sensor within each system wear day.

For the Software 505 Adult Study (SW10505), 96% of CGM provided between 1,497 and 1,992 valid glucose readings (or more than 75% of the expected number of readings) as seen in Table 9-A. Adjusted within each sensor wear-day, the CGM in the Software 505 Adult Study provided an average of 98% of all expected glucose readings (288) as seen in Table 10-A.

Table 11-A. Agreement to SMBG within CGM Glucose Ranges (Adult)

CGM Glucose Range ¹ (mmol/L)	Study ²	Number of paired CGM-SMBG	Percent within 15/15% SMBG	Percent within 20/20% SMBG	Percent within 30/30% SMBG	Percent Greater than 40/40% SMBG
0	Original	7508	69%	81%	94%	2%
Overall	Software 505	2992	77%	87%	96%	1%
2.2–3.3	Original	731	75%	84%	92%	4%
2.2-3.3	Software 505	221	73%	80%	87%	7%
2444	Original	968	78%	86%	95%	1%
3.4–4.4	Software 505	336	77%	85%	95%	1%
45.400	Original	3141	65%	78%	93%	2%
4.5–10.0	Software 505	1362	74%	85%	96%	1%
10.1–16.6	Original	1960	68%	81%	94%	3%
	Software 505	826	80%	90%	97%	1%
16.7–19.4	Original	450	77%	88%	98%	1%
	Software 505	161	83%	93%	99%	0%
19.5–22.2	Original	258	75%	85%	95%	2%
	Software 505	86	90%	93%	98%	1%

^{1.} CGM readings are within 2.2–22.2 mmol/L, inclusive.

^{2.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Table 11-B. Agreement to SMBG within CGM Glucose Ranges (Pediatric)

CGM Glucose Range ¹ (mmol/L)	Study ²	Number of paired CGM-SMBG	Percent within 15/15% SMBG	Percent within 20/20% SMBG	Percent within 30/30% SMBG	Percent Greater than 40/40% SMBG
Overall	Original	16318	64%	76%	89%	5%
Overall	Software 505	4264	73%	84%	94%	2%
2.2–3.3	Original	487	44%	55%	68%	19%
2.2–3.3	Software 505	240	54%	71%	86%	7%
3.4–4.4	Original	1340	59%	70%	85%	7%
3.4-4.4	Software 505	399	64%	76%	92%	2%
4.5–10.0	Original	7084	62%	74%	90%	5%
4.5-10.0	Software 505	1650	72%	84%	95%	2%
10.1–16.6	Original	5627	69%	80%	90%	5%
10.1-10.0	Software 505	1526	79%	89%	97%	2%
16.7.10.4	Original	1176	65%	77%	90%	4%
16.7–19.4	Software 505	319	72%	83%	94%	2%
10.5.22.2	Original	604	58%	72%	86%	6%
19.5–22.2	Software 505	130	69%	79%	86%	8%

^{1.} CGM readings are within 2.2–22.2 mmol/L, inclusive.

^{2.} Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

Table 12-A. CGM Difference to SMBG Within CGM Glucose Ranges (Adult)

CGM Glucose Ranges ¹ (mmol/L)	Study ²	Number of Paired CGM-SMBG	Mean Percent Difference	Median Percent Difference	Mean Absolute Percent Difference	Median Absolute Percent Difference
Overall	Original	7508	-0.4%	-1.4%	14.0%	11.0%
Overall	Software 505	2992	-2.6%	-2.7%	11.3%	8.6%
*2.2–3.3	Original	731	-9.3	-8.0	11.7	8.0
2.2–3.3	Software 505	221	-10.3	-6.0	13.0	8.0
*3.4–4.4	Original	968	-1.0	1.0	10.7	8.0
3.4–4.4	Software 505	336	-4.0	-2.0	10.1	7.0
4.5–10.0	Original	3141	1.4%	0.0%	14.2%	11.0%
4.5-10.0	Software 505	1362	-2.6%	-3.1%	11.4%	8.9%
10.1.16.6	Original	1960	-0.7%	-2.8%	13.0%	10.3%
10.1–16.6	Software 505	826	-1.4%	-2.0%	9.5%	7.4%
16.7–19.4	Original	450	-0.7%	-2.6%	10.5%	8.6%
10.7-19.4	Software 505	161	0.0%	0.0%	8.3%	6.0%
10 5 22 2	Original	258	5.0%	3.0%	11.9%	8.6%
19.5–22.2	Software 505	7508	-0.4%	-1.4%	14.0%	11.0%

^{1.} CGM readings are within 2.2–22.2 mmol/L, inclusive. 2. Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

^{*} For CGM \(\leq 4.4 \) mmol/L, the differences in mmol/L are included instead of percent differences (%).

Table 12-B. CGM Difference to SMBG Within CGM Glucose Ranges (Pediatric)

CGM Glucose Ranges ¹ (mmol/L)	Study ²	Number of Paired CGM-SMBG	Mean Percent Difference	Median Percent Difference	Mean Absolute Percent Difference	Median Absolute Percent Difference
Overall	Original	7508	-0.4%	-1.4%	14.0%	11.0%
Overall	Software 505	2992	-2.6%	-2.7%	11.3%	8.6%
*2.2–3.3	Original	731	-9.3	-8.0	11.7	8.0
2.2–3.3	Software 505	221	-10.3	-6.0	13.0	8.0
*3,4-4,4	Original	968	-1.0	1.0	10.7	8.0
3.4-4.4	Software 505	336	-4.0	-2.0	10.1	7.0
4.5–10.0	Original	3141	1.4%	0.0%	14.2%	11.0%
4.5-10.0	Software 505	1362	-2.6%	-3.1%	11.4%	8.9%
10.1–16.6	Original	1960	-0.7%	-2.8%	13.0%	10.3%
10.1–10.0	Software 505	826	-1.4%	-2.0%	9.5%	7.4%
16.7.10.4	Original	450	-0.7%	-2.6%	10.5%	8.6%
16.7–19.4	Software 505	161	0.0%	0.0%	8.3%	6.0%
19.5–22.2	Original	258	5.0%	3.0%	11.9%	8.6%
19.0-22.2	Software 505	7508	-0.4%	-1.4%	14.0%	11.0%

^{1.} CGM readings are within 2.2–22.2 mmol/L, inclusive. 2. Both sets of study data are presented and are labeled as Original (SW10050) or Software 505 (SW10505).

^{*} For CGM \(\leq 4.4 \) mmol/L, the differences in mmol/L are included instead of percent differences (%).

Agreement and Accuracy Relative to SMBG

Agreement between the CGM and blood glucose values is also characterized using paired CGM and SMBG results (Tables 11-A/B to 12-A/B).

The CGM and SMBG values were compared by pairing the comparative SMBG value to a CGM glucose reading that occurred immediately after the SMBG was collected. These results characterize the performance subjects expect during real-time use of the CGM in their daily diabetes management when comparing the CGM readings to their home blood glucose meter results. For readings less than or equal to 4.4 mmol/L. the absolute difference in mmol/L between the two alucose results was calculated. For values greater than 4.4 mmol/L, the absolute percent difference (%) from the SMBG values was calculated. The percentages of total readings within 0.8 mmol/ L or 15%, 1.1 mmol/L or 20%, 1.7 mmol/L or 30%, 2.2 mmol/L or 40% or greater than 2.2 mmol/L or 40% were then calculated.

For example, if the CGM reads 5.5 mmol/L, it is between 4.5–10.0 mmol/L range and you can expect the CGM readings to be within 20% of the SMBG values 85% of the time for the Software 505 Adult Study, as seen in Table 11-A.

Overall, the CGM in the Software 505 Adult Study reads, on average, 2.6% lower (Mean Percent Difference) than SMBG values and 11.3% absolute different (Mean Absolute Percent Difference) than the SMBG values. The Median Percent Difference shows that half of the time the CGM reads lower in 2.7% or less than the SMBG values and the Median Absolute Percent Difference shows that half of the time the CGM reads about 8.6% or less different than SMBG values, as seen in Table 12-A.

Adverse Events

No serious adverse events or devicerelated serious adverse events occurred during the studies. Mild to moderate skin irritation, such as erythema or edema, occurred at the sensor needle insertion area or around the adhesive area. No infection, bruising, or bleeding occurred at the sensor needle insertion area or the adhesive area.

30.7 Electromagnetic Compatibility

The information contained in this section is specific to the t:slim X2 System. This information provides reasonable assurance of normal operation, but does not guarantee such under all conditions. If the t:slim X2 System must be used in close proximity with other electrical equipment, the System should be observed in this environment to verify normal operation. Special precautions for electromagnetic compatibility must be taken when using medical electrical equipment. The System shall be placed into service with adherence to the EMC information provided here. Using cables and accessories not specified in this User Guide may adversely impact safety, performance, and electromagnetic compatibility, including increased emissions and/or decreased immunity.

For IEC 60601-1 testing, Essential Performance for the t:slim X2 System is defined as follows:

The System will not over deliver a clinically significant amount of insulin.

- The System will not under deliver a clinically significant amount of insulin without notification to the user.
- The System will not deliver a clinically significant amount of insulin after occlusion release.
- The System will not discontinue reporting CGM data without notification to the user.

This section contains the following tables of information:

- Electromagnetic Emissions
- Electromagnetic Immunity
- Distances Between the t:slim X2 System and RF Equipment

30.8 Electromagnetic Emissions

The t:slim X2 System is intended for use in the electromagnetic environment specified below. Always make sure that the System is used in such an environment.

Guidance and Manufacturer's Declaration – Electromagnetic Emissions

Emissions Test	Compliance	Electromagnetic Environment – Guidance
RF Emissions, CISPR 11	Group 1	The System uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF Emissions, CISPR 11	Class B	The System is suitable for use in all establishments,
Harmonic Emissions, IEC 61000-3-2	N/A	including domestic establishments and those directly connected to the public low-voltage power supply
Voltage Fluctuations/Flicker Emissions, IEC 61000-3-3	N/A	network that supplies buildings used for domestic purposes.

30.9 Electromagnetic Immunity

The t:slim X2 System is intended for use in the electromagnetic environment specified below. Always make sure that the System is used in such an environment.

Guidance and Manufacturer's Declaration - Electromagnetic Immunity

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment – Guidance
Electrostatic Discharge (ESD) IEC 61000-4-2	± 8 kV contact ± 15 kV air	± 8 kV contact ± ± 15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical Fast Transient/burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/ output lines (100 kHz repetition frequency)	± 2 kV for power supply lines ± 1 kV for input/output lines (100 kHz repetition frequency)	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1 kV differential mode ± 2 kV common mode	± 1 kV differential mode ± 2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.

Guidance and Manufacturer's Declaration - Electromagnetic Immunity (Part 1 of 4)

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment – Guidance
Voltage Dips, Short Interruptions, and Voltage Variations on Power Supply Input Lines IEC 61000-4-11	70% UR (30% dip in Ur) for 25 cycles 0% Ur (100% dip in Ur) for 1 cycle at 0 degrees 0% Ur (100% dip in Ur) for 0.5 cycles at 0, 45, 90, 135, 180, 225, 270, and 315 degrees 0% Ur (100% dip in Ur) for 250 cycles	70% UR (30% dip in Ur) for 25 cycles 0% Ur (100% dip in Ur) for 1 cycle at 0 degrees 0% Ur (100% dip in Ur) for 0.5 cycles at 0, 45, 90, 135, 180, 225, 270, and 315 degrees 0% Ur (100% dip in Ur) for 250 cycles	Mains power quality should be that of a typical commercial or hospital environment. If the user of the pump requires continued operation during power mains interruptions, it is recommended that the pump be powered from an uninterruptible power supply or a battery. NOTE: Ur is the a.c. mains voltage prior to application of the test level.
Power Frequency (50/60 Hz) Magnetic Field IEC 61000-4-8	30 A/m	400 A/m (IEC 60601-2-24)	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

Guidance and Manufacturer's Declaration - Electromagnetic Immunity (Part 2 of 4)

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	10 V	Portable and mobile RF
Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2,7 GHz	30 V/m	communications equipment should be used no closer to any part of the pump, including cables, than the
Proximity Field from Wireless Transmitters	385 MHz: 27 V/m @ 18 Hz pulse modulation 450 MHz: 28 V/m @ FM modulation 710 MHz, 745 MHz, 780 MHz: 9 V/m @ 217 Hz pulse modulation 810 MHz, 870 MHz, 930 MHz: 28 V/m @ 18 Hz Pulse modulation 1720 MHz, 1845 MHz, 1970 MHz: 28 V/m @ 217 Hz Pulse Modulation 2450 MHz: 28 V/m @ 217 Hz Pulse modulation 5240 MHz, 5500 MHz, 5785 MHz: 9 V/m @ 217 Hz Pulse modulation	385 MHz: 27 V/m @ 18 Hz pulse modulation 450 MHz: 28 V/m @ FM modulation 710 MHz, 745 MHz, 780 MHz: 9 V/m @ 217 Hz pulse modulation 810 MHz, 870 MHz, 930 MHz: 28 V/m @ 18 Hz Pulse modulation 1720 MHz, 1845 MHz, 1970 MHz: 28 V/m @ 217 Hz Pulse Modulation 2450 MHz: 28 V/m @ 217 Hz Pulse modulation 5240 MHz, 5500 MHz, 5785 MHz: 9V/m @ 217 Hz Pulse modulation	recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance: 150 kHz to 80 MHz, d = 1.20 VP 80 MHz to 800 MHz, d = 1.20 VP 800 MHz to 2.5GHz, d = 2.30 VP Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site surveya, should be less than the compliance level in each frequency rangeb. Interference may occur in the vicinity of equipment marked with the following symbol: ((v))

Guidance and Manufacturer's Declaration - Electromagnetic Immunity (Part 3 of 4)

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment – Guidance
---------------	----------------------	------------------	--

NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

^a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast, and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the pump is used exceeds the applicable RF compliance level above, the pump should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the System.

^b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 10 V/m.

Guidance and Manufacturer's Declaration - Electromagnetic Immunity (Part 4 of 4)

30.10 Distances Between the t:slim X2 Pump and RF Equipment

The t:slim X2 System is intended for use in an electromagnetic environment typically found in the home, at work, retail stores, and places of leisure, where daily activities occur. The chart below can be used as a guideline for determining the recommended minimum distance to maintain between a radio frequency (RF) transmitter and the t:slim X2 System. For specific concerns about a particular RF transmitter interfering with your System's operation, please contact the transmitter manufacturer for its rated power and frequency.

Poted Maximum Output Power	Separation Distance According to Frequency of Transmitter in Meters			
Rated Maximum Output Power of Transmitter in Watts	150 kHz to 80 MHz d = 1.2√P	80 MHz to 800 MHz d = 1.2√P	800 MHz to 2.5 GHz d = 2.3√P	
0.01 W	0.12 m	0.12 m	0.23 m	
0.1 W	0.38 m	0.38 m	0.73 m	
1 W	1.2 m	1.2 m	2.3 m	
10 W	3.8 m	3.8 m	7.3 m	
100 W	12 m	12 m	23 m	

For transmitters rated at a maximum output power not listed above, the recommended separation distance (d) in meters (m) can be determined using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

The table below provides a list of typical devices for various levels of transmitter power and frequency, and the recommended separation distances from the transmitter and the System.

Rated Maximum Output Power of Transmitter in Watts	Typical Devices	Recommended Separation Distance in Meters	
0.001 W	Bluetooth Class 3 (standard 1 meter range). Commonly used as Bluetooth® headset.	0.07 m	■ -(((• ⊘
0.01 W	Internet to music adapter. Commonly used for FM wireless music streaming	0.13 m	
0.1 W	Bluetooth Class 1 (100 meter range). Wireless router (WiFi).	0.73 m	
1 W	Typical microwave oven RF leakage.	2.3 m	((()

30.11 Quality of Wireless Service

The manufacturer defines the quality of service of the t:slim X2 System as the percent of readings successfully received by the display, where the transmitter and display attempt to communicate every 5 minutes. One of the t:slim X2 System's essential performance requirements states that the system will not discontinue reporting data and/or information from the Dexcom G5 transmitter to the user without notification.

The system notifies the user of a missed reading, or when the transmitter and pump are out of range of one another in several ways. The first is when a bullet point is missed on the CGM graph which will occur within five minutes of the previous reading. The second indication occurs after 10 minutes when the Out of Range Icon is displayed on the home screen. The third is a user settable alert that will notify the user when the pump and transmitter is out of range of one another. Setting this alert is defined in Section 20.6 Setting Your Out of Range Alert.

The system performance requirements state that 90% of readings will be successfully transferred to the display while the transmitter and display are within 6 meters of each other, and no more than 12 consecutive readings (1 hour) will be missed.

To improve quality of service when other devices operating in the 2.4 GHz band are around, the t:slim X2 System uses the built-in coexistence features provided by Bluetooth® technology.

30.12 FCC Notice Concerning Interference

The transmitter covered by this user's guide has been certified under FCC ID: PH29433.

Although the transmitter has been approved by the Federal Communications Commission, there is no guarantee that it will not receive interference or that any particular transmission from the transmitter will be free from interference.

Compliance Statement (Part 15.19)

This device complies with Part 15 of the FCC Rules.

Operation is subject to the following two conditions:

- 1. This device may not cause harmful interference, and
- This device must accept any interference received, including interference that may cause undesired operation.

Warning (Part 15.21)

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

FCC Interference Statement (Part 15.105 (b))

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protec-

tion against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no quarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This portable transmitter with its antenna complies with FCC/IC RF exposure limits for general population/uncontrolled exposure.

THIS PAGE IS INTENTIONALLY LEFT BLANK

Index

A

AC power adapter 74
Accessories 40, 74
Active Bolus Icon 50, 52
Additional Pump Operations 108
After Bolus BG Reminder 125
Air Bubbles

Check Tubing 87
Removing 83
Removing Before Delivery 87

Airport Security 261

Alarms 127, 130, 155
Altitude Alarm 167, 168
Cartridge Error Alarm 161
Cartridge Removal Alarm 162
Empty Cartridge Alarm 160
Low Power Alarms 158
Occlusion Alarms 164
Resume Pump Alarm 157
Screen On/Quick Bolus Button
Alarm 166

Temperature Alarm 163

Time to Occlusion Alarm,

Specifications 282 Alerts 127, 130, 131, 133 Alert Icon. Where to Find 50 Alert Settings 64 Basal Rate Required Alert 142 CGM Rise and Fall 188 CGM. 12 Hour Calibration Alert 231 CGM, Calibrate CGM Alert 235 CGM, Calibration Error Alert 234 CGM, Calibration Timeout Alert 233 CGM. Error 248 CGM. Failed Sensor 247 CGM, Fall Alert 241, 242 CGM, High Alert 236 CGM, Incomplete Calibration Alert 232 CGM. Low Alert 237, 238 CGM, Low Transmitter Battery Alert 245 CGM, Out of Range Alert 244 CGM, Rise Alert 239, 240 CGM, Second Startup Calibration Alert 230 CGM, Startup Calibration Alert 229

CGM. System Error 249 CGM. Transmitter Error 246 Connection Frror Alert 150 Data Error Alert 153 High Glucose Alert, Set 186 Incomplete Bolus Alert 136 Incomplete Cartridge Change Alert 138 Incomplete Fill Cannula Alert 140 Incomplete Fill Tubing Alert 139 Incomplete Personal Profile Alert 141 Incomplete Setting Alert 141 Incomplete Temp Rate Alert 137 Low Glucose Alert, Set 187 Low Insulin Alert 128 Low Power Alerts 134 Max Basal Alerts 146 Max Bolus Alerts 144 Max Hourly Bolus Alert 143 Min Basal Alerts 148 Out of Range Alert, Set 189 Out of Range, Set 189 Power Source Alert 152 Alternate Site BG Testing 42 Altitude 260

Altitude Alarm 167, 168 Target BG 40, 92 C Target BG in Personal Profiles 94 Aguatic Activities, CGM 260 Calculation 60 **BG Reminder** 125 Aquatic Activities, Pump 260 Calibrate CGM Alert 235 Bluetooth 37, 180 Arrows Calibrate Your CGM 205 Bluetooth, Recommended Distance Menu Arrows 62 Calibrate, Reasons to 209 Between Devices 334 Up/Down Arrows 62 Calibration Alert, 12 Hour 231 **Bolus** 40, 99 Calibration Error Alert 234 В Active Bolus Icon 50. 52 Calibration Timeout Alert 233 After Bolus BG Reminder 125 Basal 40 Cancel a Bolus 104 Bolus Overview 100 Basal Rate Required Alert 142 Cannula 41 Bolus Screen 60 Current Basal Rate 58 Cannula, Fill Cannula 89 Cancel a Bolus 104 Delivery Accuracy 271 Car Adapter, Specifications 276 Correction Bolus 41 Frequency of Delivery 271 Carb Ratio 41 Deliver a Quick Bolus 114 in Personal Profiles 94 in Personal Profiles 94 Delivery Accuracy 271 Set a Temp Rate 110 Timed Settings 93 Extended Bolus 41, 103 Stop a Temp Rate 110 Carbohydrates 41, 58 Food Bolus Using Grams 102 Temp Basal Rate 41 Food Bolus Using Units 102 Carbs, in Personal Profiles 93 Timed Settings 93 Food Bolus Using Grams 102 in Personal Profiles 95 Basal Rate Required Alert 142 Food Bolus, on Bolus Screen 60 Quick Bolus 41, 113 **Battery** turn on in Personal Profiles 96 Rate of Delivery Specifications Battery Level 50, 54 281 Carbs 41 Charging Tips 75 Stop a Bolus 104 Carbs. on Bolus Screen 60 Battery, Charge 74 Timed Settings 93 Caring for Your Pump 263 **BG** 40 Cartridge 81 High BG Reminder 125 Cartridge Tubing 54 Low BG Reminder, Set 124

Change Cartridge 86 CGM Settings 180 Pair Your CGM 180 Fill Cartridge 84 CGM System Error 249 Rate of Change Arrows 214 Load Cartridge 82, 86 Clinical Studies, Sensor 284 Reasons to Calibrate 209 Default Volume 181 Re-Calibrate 209 Load the Cartridge 41 Cartridge Error Alarm 161 Distance to Pump and other Receiver 176 devices 333 Rise and Fall Alerts 188 Cartridge Removal Alarm 162 End a Sensor Session 222 Second Startup Calibration Alert CGM Enter Transmitter ID 180 230 12 Hour Calibration Alert 231 Failed Sensor 247 Sensor Expiration Date 265 12 Hour Calibration Update 208 Failed Sensor, Troubleshooting Sensor Inaccuracies. Alerts and Errors 227 Troubleshooting 255 255 Automatic Sensor Shut-Off 222 Glucose Trend Arrows 214 Sensor Specifications 273 Calibrate Blood Glucose Value Glucose Trend Graphs 213 Sensor Startup Period 203 208 High BG Alert Default 186 Set Correction Bolus 208 Calibrate CGM Alert 235 High BG Alert Repeat 186 Set Volume 181 Calibrate Your CGM 205 High Glucose Alert, Set 186 Start a Sensor Session 201 Calibration Error Alert 234 History, View 219 Start Calibration 207 Calibration Overview 206 Incomplete Calibration Alert 232 Startup Calibration Alert 229 Calibration Prompts 218 Low BG Repeat 187 Status Symbols 218 Calibration Timeout Alert 233 Low Glucose Alert, Set 187 System Overview 174 CGM Alerts 185 Low Transmitter Battery Alert 245 Transmitter Frror 246 CGM Frror 248 My CGM Screen 66 Transmitter ID 180 CGM Fall Alert 241, 242 Out of Range Alert 244 Transmitter Performance CGM High Alert 236 Specifications 283 Out of Range Alert, Set 189 CGM Info 184 Transmitter Specifications 274 Out of Range/No Antenna, CGM Low Alert 237, 238 Troubleshooting 254 Troubleshooting 251 CGM Rise Alert 239, 240

Unknown Sensor Reading 243	Contraindications 24	Site Reminder 126
Unknown Sensor Reading,	Correction Bolus 41	Delete Personal Profile 98
Troubleshooting 253	Correction Factor 41, 93	Delivery Summary 122
View Data on Pump, Overview	in Personal Profiles 94	Disconnect 27, 31
212	Timed Settings 93	Disconnect When Filling 87
CGM Error 248	_	Display Settings 118
CGM Fall Alert 241, 242	D	Disposing of System Components
CGM High Alert 236	Data Error Alert 153	265
CGM ID 180	Data, View CGM Overview 212	_
CGM Low Alert 237, 238	Date	E
CGM Rise Alert 239, 240	Date and Time Display 50	Edit
CGM Settings 179	Edit Date 77	Edit Date 77
CGM System Error 249	Default	Edit Time 77
Charge	Auto-Off Alarm 128	Site Reminder 89
AC Power Outlet 74	Basal Temp Rate 110	Electromagnetic Compatibility 326
Car Adapter 74	CGM Default Volume 181	Electromagnetic Emissions 328
Charging Tips 75	CGM Fall Alert 188	Electromagnetic Immunity 329
Personal Computer 75	CGM Out of Range Alert 189	Emergency Kit 36
Charge the Pump 74	CGM Rise Alert 188	Empty Cartridge Alarm 160
Cleaning Your System 264	Extended Bolus 103	End CGM Sensor Session 222
Colors	High BG Reminder 125	Extended Bolus 41, 103
Explanation of System Colors 48	High Glucose Alert 186	Default 103
Computer Charging Specifications	Low BG Reminder 124	
277	Low Insulin Alert 128	F
Connection Error Alert 150	Quick Bolus 114	Failed Sensor Error 247
Contents, of System 40	Screen Timeout 118	

Fall Alert, Set 188 Feature Lock Pediatric 19 Fill Fill Cannula 89 Fill Cartridge 84 Fill Port 82, 84	Icons Explanation of Icons 46 Incomplete Bolus Alert 136 Incomplete Calibration Alert 232 Incomplete Cartridge Change Alert 138	Insulin on Board (IOB), in Personal Profiles 93 Interference, FCC Notice 335 K Keypad 68 Letter Keypad 70
Fill Tubing 87 G Glucose Rate of Change Arrows 214 Glucose Trend Graphs 213 Grams Food Bolus, on Bolus Screen 60 Food Bolus, using 102	Incomplete Fill Cannula Alert 140 Incomplete Fill Tubing Alert 139 Incomplete Setting Alert 141 Incomplete Temp Rate Alert 137 Indications for Use 24 Infusion Set Risks 35, 80 Infusion Site Care 80 Pediatric 19 Infusion Site Care, Pediatric 19	Number Keypad 68 L LED 48 LED, location on Home Screen 54 Lifestyle Issues 259 Load Cartridge 82, 86 Low BG Reminder 124 Low Insulin Alert 128, 133 Low Power Alarms 158
Healthcare Provider 36 High BG Reminder 125 History CGM History 219 Pump History 122 Home Screen 54 Home Screen, CGM 56	Insulin Insulin Duration 93 Insulin Level Display 88 Insulin On Board (IOB) 41, 50 Insulin On Board (IOB) Display 50 Resume Insulin Delivery 106 Stop Insulin Delivery 106 Insulin Duration, in Personal Profiles 93	Low Power Alerts 134 Low Transmitter Battery Alert 245 M Maintaining Your Pump 263 Malfunction 170 Max Basal Alerts 146 Max Bolus Alerts 144

Max Bolus, in Personal Profiles 93	Create a New Profile 92	R	
Max Hourly Bolus Alert 143	Delete a Profile 98	Receiver, CGM 176	
Min Basal Alerts 148	Edit or View 96	Reminder	
Missed Meal Bolus Reminder 126	Personal Profiles Overview 92	Site Reminder 89	
My Pump Screen 64	Program a Personal Profile 94	Reminders 123, 127	
_	Rename a Profile 98	High BG 125	
0	Power Adapter, AC 74	Low BG 124	
Occlusion Alarms 164	Power Source Alert 152	Missed Meal Bolus 126	
Occlusion Alarms, Specifications	Precautions	Site Reminder 126	
282	Using t:slim X2 Insulin Pump 29	Resume Insulin Delivery 106	
Options Screen 62	Pump History 122	Resume Pump Alarm 157	
Out of Range Alert 244	Pump History, Delivery Summary	Rise Alert, Set 188	
Overview	122	Risks from Using the System 35	
CGM Overview 174	Pump Info 122	Thisks from Osing the System 33	
System Description 18	Pump Info, Serial Number 122	S	
	Pump Performance, Specifications	Safety Information 23	
P	281	Contraindications 24	
Pediatric	Pump Settings 117	Emergency Kit 36	
Feature Lock 19	Pump Settings, Specifications 278	Indications for Use 24	
Important Pediatric User	Pump Specifications 270	Precautions 29	
Information 19	Pump Volume 181	Verification of Proper Functionality	
Infusion Site Care 19, 80		37	
Personal Profiles	Q	Screen Lock 76	
Activate a Profile 97	Quick Bolus 19, 41, 113	Screen On/Quick Bolus Button Alarm	
Add Profiles 96	Deliver a Quick Bolus 114	166	
Copy an Existing 97	Pediatric 19	Screen Protector 40	

Screen Timeout 118	Out of Range/No Antenna,	Distance Between CGM, Pump,	
Screen Timeout, set 118	Troubleshooting 254	and other devices 333	
Screens	Reasons to Calibrate 209	Electromagnetic Compatibility 326	
Bolus Screen 60	Re-Calibrate 209	Electromagnetic Emissions 328	
CGM Home Screen 56	Troubleshooting 251	Electromagnetic Immunity 329	
CGM Screen Lock 52	Troubleshooting Sensor Reading	Pump 270	
Home Screen 54	253	Pump Performance 281	
Letter Keypad Screen 70	Unknown Reading 243	System Specifications 269	
My CGM Screen 66	Sensor Startup Period 203	Time to Occlusion Alarm 282	
My Pump Screen 64	Sensor, Start a Session 202	USB Cable 275	
Number Keypad Screen 68	Sensor, Start Calibration 207	Water Resistance 270	
Options Screen 62	Serial Number 44, 122	Start a CGM Sensor Session 201	
Screen Lock 50	Set CGM Volume 181	Startup Calibration Alert 229	
Status Screen 58	Settings, Pump Setting	Status Screen 58	
Unlock 76	Specifications 278	Stop a Bolus 104	
Second Startup Calibration Alert,	Site Reminder	Stop a Temp Rate 110	
CGM 230	Set Site Reminder 126	Stop CGM Sensor Session 222	
Sensor	Site Reminder, Set 89	Stop Insulin Delivery 106	
Applicator 42	Sound Volume 119, 181	Storing Your System 264	
Automatic Shut-Off 222	Specifications	Supply Reordering 40	
Calibration Update 208	Car Adapter 276	Symbols, Explanation of 44	
CGM Clinical Studies 284	CGM Sensor 273	System Contents 40	
Expiration Date 265	CGM Transmitter 274		
Failed Sensor, Troubleshooting	CGM Transmitter Performance	Т	
255	283	Tandem Logo 54, 76	
Out of Range Alert 244	Computer Charging 277	Target BG 40	

in Personal Profiles 94 Target BG, in Personal Profiles 92 Timed Settings 93 **Technical Specifications** 267 Temp Rate Stop a Temp Rate 110 Temp Rate, Set a Temp Basal Rate 110 Temperature Alarm 163 Temperature, Extreme 261 Time Date and Time Display 50 Edit Time 77 Time Segments 92 Time Segments, in Personal Profiles 94 **Time Segments** add to Personal Profile 95 Time to Occlusion Alarm, Specifications 282 Timed Settings 93 in Personal Profiles 94 Transmitter Safety Lock 42 Transmitter Error 246 Transmitter ID 180

Travel 261 Traveling With Your Pump 259 Traveling, by Air 261 Trend Graphs, Glucose Trends, Arrows 213 Troubleshooting CGM 251 **Tubing** Cartridge Tubing 54 Fill Tubing 87 Tubing Connector 54, 82, 87 U Units 41 Food Bolus, on Bolus Screen 60 Food Bolus, Using Units 102 Units, on Bolus Screen 60 Unknown Sensor Reading 243 Unlock the Screen 76 USB Cable Specifications 275 USB Adapter 74 USB Cable 40, 74 USB Port 54, 74



Verification of Proper Functionality
37
View Calculation 60
Volume
Pump Volume 119



Water Resistance, CGM 260
Water Resistance, Pump 260
Watertight, Pump 260



San Diego, California 92121 USA www.tandemdiabetes.com

© 2018 Tandem Diabetes Care, Inc. All rights reserved. Tandem Diabetes Care, touch simplicity, t:slim are registered trademarks in the US, and t:slim X2 is a trademark, registered trademark, and/or the subject of pending trademark applications in other territories worldwide.

For a list of patents, please see www.tandemdiabetes.com/legal/patents.

